




NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

PRODUCTS WITH A MISSION®



*Increasing
Survivability
on the Battlefield
and Beyond...*

PRODUCTS WITH A MISSION®

ICONS



MADE IN USA: Compliant with the Buy American Act mandating "substantially all" materials, manufacturing & labor are of U.S. content



CE: Meets the requirements of the European Medical Device Directive and appropriate Quality System standards



LATEX FREE: Product is considered latex free and therefore safe for latex sensitive individuals



STERILE: Product is sterile - sterility guaranteed unless packaging is opened or damaged



RED TIP TECHNOLOGY®: Specifically branded with NAR's signature Red Tip Technology® packaging. Product packaging has been engineered for less-than-optimal, often hostile conditions that exist when our equipment is needed the most



NEW PRODUCTS: Marks the newest products throughout the catalog, making them easily recognizable as our newest innovations



Table of Contents

volume 18 • issue 1

M.A.R.C.H. Products

Products designed to address Massive Hemorrhage, Airway Management, Respiratory Support, Circulation/Vascular Access and Hypothermia Prevention

4

Combat Casualty Response Kits • CCRK®

Kits designed to meet the needs/challenges of our customers and to address the leading causes of potentially preventable death

22

DoD Bleeding Control Response Kits • D-BCRK

Individual kits, packs and stations designed to equip bystanders and first responders with the essential medical gear needed to help save lives

48

TAC-EVAC

The most effective and highest quality evacuation platforms and accessories for use in the tactical environment

52

Medical/Surgical Products

Medical and surgical equipment for use in austere environments, during prolonged field care situations and in forward operating field hospitals

60

Rescue Equipment

Includes an array of rescue products ranging from lightweight, versatile Pegasus® Ropes to tactical lights and accessories

74

Treatment Support

From personal protection products to splinting and immobilization products to help from the point of injury throughout the entire continuum of care

86

Training Products & Services

NAR Training Course Overview • Products designed to train users on TCCC protocols to reduce preventable death at the point of wounding

97

Product Index

Product offerings indexed by name, NSN number, and Item number for quick reference and information

111



M•A•R•C•H



pg 6-9



pg 10-12



pg 13-15



pg 16-20



pg 21



MASSIVE HEMORRHAGE | Tourniquets

C-A-T® Combat Application Tourniquet - Gen 7

- Scientifically proven 100% effective in occluding blood flow in both upper and lower extremities by the USAISR and recommended by CoTCCC
- Affordable, lightweight, rugged design for true circumferential pressure featuring patented free-moving internal band, proprietary Red Tip Technology® elliptical tip, and new single routing buckle for simplified training with single protocol for all applications
- Reinforced windlass and optimized windlass securing system with writeable time stamp to document application time

SPECIFICATIONS

Packaged: H 6.5 in. x W 2.4 in. x D 1.5 in.
Device: L 37.5 in.
Weight: 2.7 oz



online video available

< scan to view



CAT TOURNIQUET

COLOR	STYLE	ITEM #	NSN #
Black	Tactical	30-0001	6515-01-521-7976
Orange	Rescue	30-0023	6515-01-612-7129
Blue	Trainer	30-0033	6910-01-560-2972

AAJT-S

Abdominal Aortic and Junctional Tourniquet -STABILIZED



- Easy to use and quick to apply – Lightest, most compact junctional tourniquet on the market
 - Stabilizes the pelvis in addition to stopping bleeding
 - Large area of effective pressure for greater reliability
 - Provides stable and complete occlusion of blood flow to the lower extremities, upper extremities, or pelvis
 - Can be applied at junctional sites (inguinal, axillary) or mid-abdomen
 - Demonstrated ability to remain in place and effective during patient movement
 - Similar hemostatic efficacy to Zone 3 REBOA¹
- ¹Rall, Jason M., et al. "Comparison of Zone 3 Resuscitative Endovascular Balloon Occlusion of the Aorta and the Abdominal Aortic and Junctional Tourniquet in a Model of Junctional Hemorrhage in Swine." *Journal of Surgical Research*, vol. 226, 2018, pp. 31-39.
- Lower training burden compared to REBOA procedure

AAJT-S

ITEM #	NSN #
30-0248	6515-01-695-4457

SPECIFICATIONS

Packaged: L 9 in. x W 8 in. x D 2 in.
Stored: L 7.5 in. x W 6.5 in. x D 2 in.
Weight: 17 oz

JETT® Junctional Emergency Treatment Tool

- Immediately treat high extremity wounds in the thigh/groin area that is not amenable to TQ application
- Unilateral or bilateral treatment solution
- Simultaneous occlusion of blood flow to both lower limbs that will not impede respiration

JETT

ITEM #	NSN #
30-0088	6515-01-616-5841



online video available

< scan to view

SPECIFICATIONS

Packaged: H 5.5 in. x W 7 in. x D 3.75 in.
Overall Device: L 56.75 in.
Buttocks Pad: L 19 in. x W 6.5 in.
Compression Pad: L 4.85 in. x W 2.5 in. x D 11 in.
T-Handle: H 4.62 in. x W 4.04 in. x D 0.63 in.
Weight: 1 lb 9.6 oz



See NARescue.com for more info about the **JETT Pack: Item# 80-0355 & 80-0356** a nylon pouch designed specifically for carrying the JETT



INTRO
M.A.R.C.H
KITS
TAC-EVAC
MED-SURG
RESCUE
TREATMENT SUPPORT
TRAINING
INDEX

RIGID C-A-T® TQ CASE *with Shirt Shield*

- Allows carriage of a Gen 7 C-A-T® Tourniquet on your belt while protecting your shirt from undue wear caused by friction from the windlass clip
- Blade-Tech TekLok Belt Attachment fits up to a 2.25 in. belt
- Can be mounted either vertically or horizontally
- Injection Molded Nylon Case stands up to rugged usage
- Cross medical identifier patch included



RIGID CAT TQ CASE

COLOR	ITEM #
Black	80-0944
Coyote	80-0945
Ranger Green	80-0946

SPECIFICATIONS

H 5.88 in. x W 2.5 in. x D 2 in.
Weight: 5.36 oz

Explore other TQ carrier options available at NARescue.com

Nylon C-A-T® HOLDER



- PALS/MOLLE-style connectors provide versatile attachment capabilities (mounts horizontally or vertically)
- Includes easy-open pull elastic tab and tourniquet "TQ" ID patch
- Compatible with both Gen 6 & Gen 7 C-A-T® Tourniquets (sold separately)

SPECIFICATIONS

L 6.25 in. x W 1.75 in. x D 1.5 in.
Weight: 2.17 oz

NYLON CAT HOLDER

NATO NSN*

COLOR	ITEM #	NSN #
Black	30-0057	----
Coyote	30-0044	----
MultiCam	30-0059	----
Ranger Gr.	30-0208	----
OD Green	30-0053	6530-01-605-9432*

Cummerbund C-A-T® HOLDER

- Allow attachment to tactical armor such as assault and plate carriers
- Utilizes 4 in. elastic coupled with hook and loop attachment to external cummerbund flap
- Rapid, easy access with ambidextrous one handed TQ removal
- Compatible with both Gen 6 & Gen 7 C-A-T® Tourniquets (sold separately)



CUMMERBUND CAT HOLDER

COLOR	ITEM #
Black	30-0212
Coyote	30-0213

SPECIFICATIONS

L 5 in. x W 4 in. x D 0.25 in.
Weight: 0.6 oz

Elastic C-A-T® HOLDER

- PALS/MOLLE-style strap with tuck finish eliminating the need for a button snap
- Versatile attachment to body armor, gear or belt for rapid & easy access
- Easy-open quick release hood & extended elastic body for TQ release outwards under pressure
- Compatible with both Gen 6 & Gen 7 C-A-T® Tourniquets (sold separately)



SPECIFICATIONS

L 7 in. x W 1.75 in. x D 1.5 in.
Weight: 2 oz

ELASTIC TQ HOLDER

COLOR	ITEM #
Black	30-0184
Coyote	30-0185
Multi-Cam	30-0186
Ranger Green	30-0191



MASSIVE HEMORRHAGE | Hemostatics

COMBAT GAUZE® QuikClot Z-Fold Hemostatic

- Sterile, non-woven Z-fold gauze impregnated with kaolin for temporary external control of bleeding from penetrating trauma
- Effective, simple and safe to use on wounds not amenable to tourniquet application
- Flexible & pliable, contours to all wounds
- Does not generate heat
- Individually wrapped in easy-rip, military-grade pouch
- Features an embedded x-ray detectable strip

SPECIFICATIONS

Standard Gauze: L 12 ft x W 3 in.
 Packaged: H 5 in. x W 5.5 in. x D 0.625 in. Weight: 0.8 oz
 XL Gauze: L 4 yd x W 4 in. (2 ply)
 Packaged: H 6.5 in. x W 6.5 in. x D 0.75 in. Weight: 1.75 oz



COMBAT GAUZE

SIZE	ITEM #	NSN #
Standard	30-0039	6510-01-562-3325
XL	30-0070	6510-01-600-4705



See page 105 for information on the Combat Gauze® Trainer Version

CHITOGAUZE®


SPECIFICATIONS

Gauze: L 12 ft x W 3 in.
 Packaged: H 6 in. x W 5 in. x D 0.65 in.
 Weight: 0.8 oz



CHITOGAUZE

ITEM #	NSN #
30-0049	6510-01-591-7740

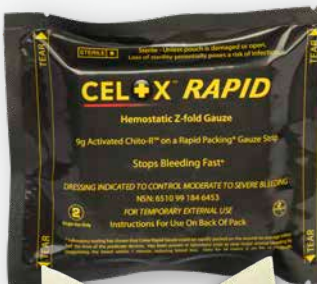

CELOX™ RAPID

SPECIFICATIONS

Gauze: W 3 in. x 5 ft
 Packaged:
 H 5.25 in. x W 4.25 in. x D 0.5 in.
 Weight: 11 oz



CELOX RAPID

ITEM #	NSN #
30-0069	6510-01-623-9903



NUSTAT®

- Patented dual-fiber technology controls bleeding through absorption at the injury site, which results in a physical aggregation of blood cells and clotting factors
- Compact packaging allows for convenient storage and rapid application

SPECIFICATIONS

Pads: 4 in. x 4 in. (Pk of 2)
 Tactical: W 2 in. x L 36 in.



NUSTAT

OPTION	ITEM #
4 in. x 4 in.	30-0236
Tactical	30-0170





Emergency Trauma Dressing **ETD™**

- Premium cotton construction for max quality and absorption
- High quality elastic wrap with a sterile non-adherent pad
- Durable C-clasp securing device



SPECIFICATIONS

4 in. ETD: L 70 in. x W 4 in. Pad: H 6 in. x W 4 in. Weight: 2.6 oz Std Pkg: L 4.5 in. x W 2.5 in. x D 1.75 in. Flat Pkg: L 4.5 in. x W 3.5 in. x D 1.13 in.	Abd/Stump ETD: L 70 in. x W 6 in. Weight: 5oz Pad: H 12 in. x W 12 in. Pkg: L 5.25 in. x W 6 in. x D 1.51 in.
6 in. ETD: L 70 in. x W 6 in. Pad: H 6 in. x W 6 in. Weight: 4 oz Std Pkg: L 6 in. x W 2.5 in. x D 2.1 in. Flat Pkg: L 6 in. x W 3.5 in. x D 1 in.	8x10 inch ETD: L 70 in. x W 6 in. Weight: 4oz Pad: H 8 in. x W 10 in. Pkg: L 4.13 in. x W 7 in. x D 1 in.

EMERGENCY TRAUMA DRESSINGS

SIZE	STYLE	ITEM #	NSN #
4 inch	Standard	30-0031	6510-01-558-4108
4 inch	Flat	30-0093	----
6 inch	Standard	30-0032	6510-01-558-4114
6 inch	Flat	30-0094	----
Abd/Stump	----	30-0012	6510-01-541-8121
8x10 inch	----	30-0165	6510-01-654-6801



online video available
< scan to view

COMPRESSED GAUZE

- Low-cube, ultra-compact sterile gauze
- Premium rolled 100% cotton gauze for maximum absorption

SPECIFICATIONS



Gauze:
W 4.5 in. x L 4.1 yd
Packaged:
L 3 in. x W 2.5 in. x D 1 in.
Weight: 12 oz



COMPRESSED GAUZE

ITEM #	NSN #
30-0052	6510-01-503-2117

S-ROLLED GAUZE

- Unique dispenser, allows controlled linear feed application
- Premium 100% cotton

SPECIFICATIONS



Gauze:
W 4.5 in. x L 4.1 yd
Packaged:
L 3 in. x W 2.5 in. x D 1 in.
Weight: 12 oz



S-ROLLED GAUZE

ITEM #	NSN #
30-0003	6510-01-605-2234

CRAVAT (Triangular Bandage)

Bandage: 37 in. x 37 in. x 52 in.
Pkg: 3.5 in. x 2 in. x 0.63 in. Wgt: 1.6 oz



CRAVAT

ITEM #	NSN #
30-0089	6510-00-201-1755



TACTICAL COMPRESSION WRAP

Wrap: L 38 in. x W 5 in. Wgt: 6oz
Pkg: 5.4 in. x 2.7 in. x 2 in.



TAC-WRAP

ITEM #	NSN #
30-0016	6510-01-614-0938



BAC-PACK™ *Bougie-Aided Cricothyroidotomy Pack*

- All-in-one cricothyroidotomy kit for use in any environment
- Pre-inserted bougie in tube allows for a more rapid, easy procedure
- Tube meets TCCC recommendations for bougie-aided open surgical technique
- #10 Scalpel included for safety and ease of use
- Small, durable pouch keeps the items organized and ready to use
- Lightweight and ruggedized
- Components configured by order of use
- Scalpel, tube, syringe all provided sterile
- Patent pending

INCLUDED IN THE BAC-PACK:






SPECIFICATIONS


Packaged: H 8 in. x W 4 in. x D 1.75 in.
 Weight: 5.2 oz
 U.S. Patents: Pending

BAC-PACK

ITEM #
10-0052

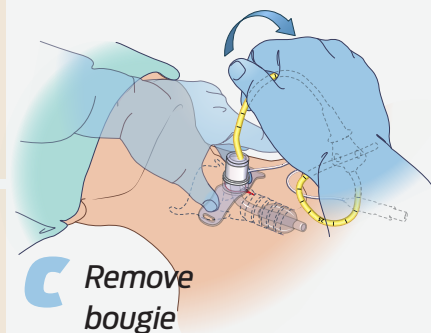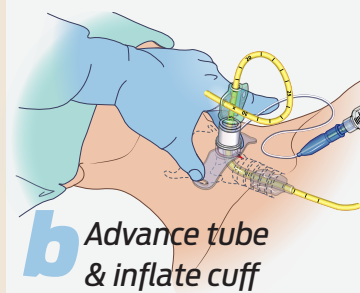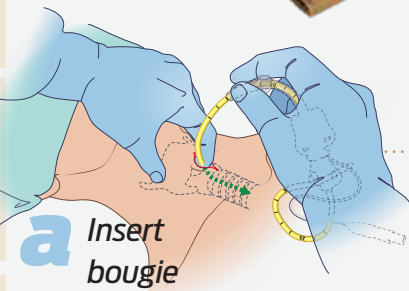

NAR Tracheostomy Kit with BOUGIE INTRODUCER

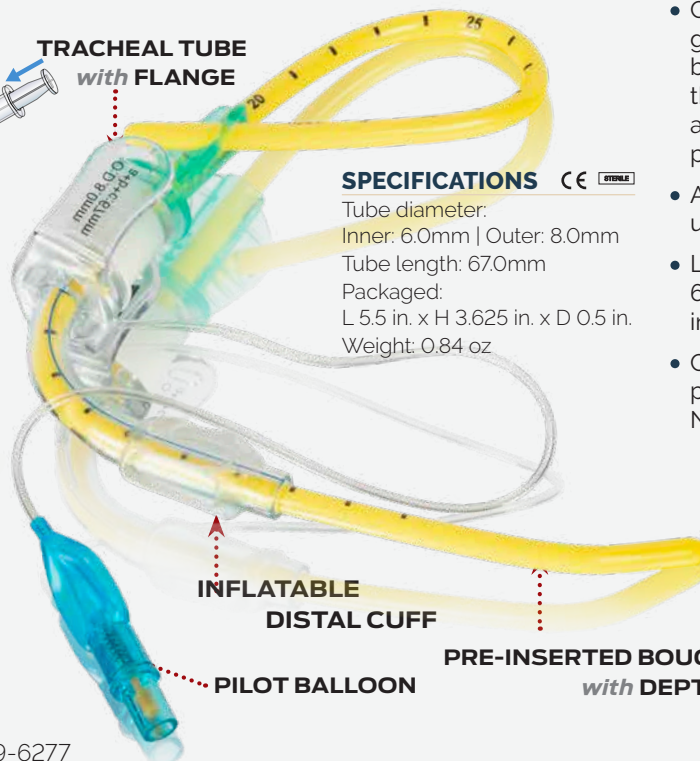

SPECIFICATIONS

Tube diameter:
 Inner: 6.0mm | Outer: 8.0mm
 Tube length: 67.0mm
 Packaged:
 L 5.5 in. x H 3.625 in. x D 0.5 in.
 Weight: 0.84 oz

- Convenient device that gives the user a pre-inserted bougie for introduction of the tube into the trachea and indication of correct placement prior to insertion
- Adjustable for maximum user preference
- Latex free tube features 6.0mm inner diameter and inflatable cuff
- Contains the same preferred flexible tube as NAR's Tactical CricKit®

TRACHEOSTOMY KIT WITH BOUGIE

ITEM #
10-0048

TACTICAL CRICKIT[®] *Tactical Surgical Cricothyroidotomy Kit*

- Contains all features required by CoTCCC for a surgical airway kit in the tactical setting
- Ideal length & size tracheal tube with securing device
- Configured in the sequence of use, ideal for high-stress, low-light, austere environments
- Lightweight, ruggedized and minimal in cube size
- Included rolled pouch provides a clean working surface



INCLUDED IN THE TACTICAL CRICKIT:



QTY **1** Nylon Pouch

QTY **2** Wipes, Antimicrobial

QTY **1** Gauze Pads 2"x2" (2pk)

QTY **1** #10 Scalpel

QTY **1** Tracheal Hook, NAR

QTY **1** Securing Strap

QTY **1** 6mm CircTube

QTY **1** Syringe Luer-Lock 10cc

SPECIFICATIONS

Packaged: H 6.5 in. x W 2.5 in. x D 1.25 in.
Unrolled: H 6.5 in. x W 7.5 in. x D 1.25 in.
Weight: 3.7 oz

TACTICAL CRICKIT



ITEM #	NSN #
10-0017	6515-01-540-7568



NAR TRACHEAL HOOK



NAR TRACHEAL HOOK

ITEM #	NSN #
10-0020	6515-01-562-9734

SPECIFICATIONS

L 4.2 in. x W 0.5 in. x D 0.8 in.
Weight: 0.5 oz

NAR TACTICAL SUCTION DEVICE



SPECIFICATIONS

L 9.75 in. x W 3.75 in. x D 3 in. | Wgt: 7.3 oz
Suction Capacity: 1000 mL
Vacuum Force: 100 mmHg

TACTICAL SUCTION

ITEM #	NSN #
10-0018	6515-01-540-7206

NASOPHARYNGEAL AIRWAY

with lubricant

SPECIFICATIONS

L 6 in. x W 16 in. x D 0.5 in.
Weight: 0.5 oz



NASO WITH LUBE

SIZE	ITEM #	NSN #
28F	ZZ-0034	6515-01-529-1187



22F-36F available at NARescue.com

KING LTS-D[™] *Supraglottic Airway Kit*

- Unique dual channel design, anatomical shaped distal tip/cuff & open second lumen
- Tapered ramp design provides additional ventilation outlet

SPECIFICATIONS

Folded: H 8 in. x W 11.43 in. x D 1.92 in.
Unfolded: H 8 in. x W 13.36 in. x D 1.22 in.
Weight: 5.7 oz

KING LTS-D

CLR	SIZE	ITEM #	NSN #
YEL	3	10-0025	----
RED	4	10-0026	6515-01-591-1578
PUR	5	10-0027	6515-01-591-0397

INCLUDED WITH KING LTS-D[™]



QTY **1** Sterile Lubricant

QTY **1** Gastric Tube, 18F

QTY **1** Syringe Luer-Lock 60cc





AIRWAY MANAGEMENT

INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

AMBU® RDIC *Mark III Resuscitator*

- Genuine Ambu Mark III reusable resuscitator equipped with impermeable elastomeric components on the compression bag and patient valve (butyl rubber)
- Specialized inlet filtration system, including C2 canister, for filtering contaminated air

SPECIFICATIONS

L 14 in. x W 7.5 in. x H 7 in.
Weight: 5.4 lb

AMBU RDIC

CE

ITEM #	NSN #
10-0047	6515-01-338-6602



CYCLONE® Improved Pocket BVM

- Unique collapsible design packaged 25% smaller than conventional BVMs
- Requires minimum pressure for compression of the bag to ventilate



SPECS

Diameter: 5.3 in. x H 2.9 in.
Weight: 17.4 oz

CYCLONE BVM

CE

ITEM #	NSN #
10-0019	6515-01-568-0193

AMBU® SPUR II BVM


SPECIFICATIONS

L 11.63 in. x W 5 in. x H 5 in.
Weight: 11.2 oz

AMBU SPUR BVM

CE

ITEM #	NSN #
10-0043	6515-01-519-3612

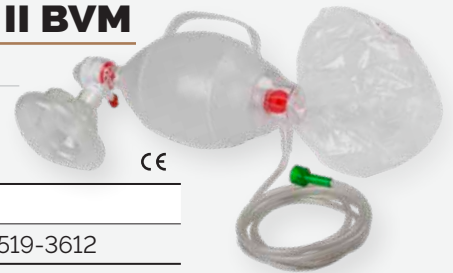

Special Operations

LARYNGOSCOPE SET

- Perform direct visual laryngoscopy even in the most severe environments
- Lightweight, compact kit to equip medical care providers in the field with tools to secure an airway
- Contains 2 single-use fiber optic blades, textured handle & removable light source all rolled up in a rugged nylon pouch

SPECIFICATIONS

Open: L 8 in. x W 5.5 in. x D 2 in.
Rolled: L 8 in. x W 2 in. x D 2 in.
Weight: 10 oz

LARYNGOSCOPE SET

ITEM #	NSN #
10-0001	6515-01-450-9790



INCLUDED IN THE LARYNGOSCOPE SET:



QTY 1 Pouch, Nylon



QTY 1 Light, Pocket Exam



QTY 1 #3 Miller Blade



QTY 1 #3 Mac Blade



QTY 1 Handle Layrngoscope



QTY 1 Spare Bulb




QTY 2 Size AA Batteries

NAR Fiber Optic LARYNGOSCOPE BLADES

LARYNGOSCOPE BLADES

CE

AVAILABLE SIZES	ITEM #	NSN# available on all sizes offered
Miller #1, #2, #3 & #4	See NARescue.com	See NARescue.com for details
Mac #2, #3 & #4	See NARescue.com	See NARescue.com for details

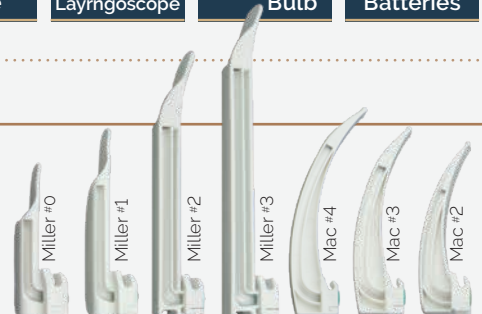




HYFIN[®] Chest Seals

Twin Pack

- Two individually packaged, vented/valve-type chest seals in one unique package designed for treatment of both entry/exit wounds or multiple penetrating injuries.
- Patented 3 channel pressure relief vents
- Ruggedized dual perforated pouch designed to allow the provider to open only one dressing at a time, as needed
- Advanced adhesive technology designed to provide a superior seal

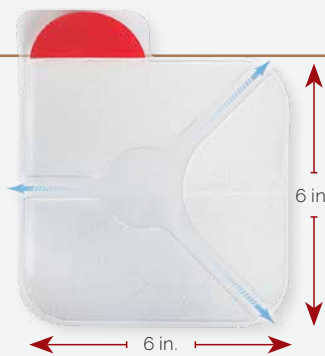

HYFIN CHEST SEAL TWIN PACK



ITEM #	NSN #
10-0037	6510-01-642-6212

SPECIFICATIONS

Device: 6 in. x 6 in. (excludes tab) Wgt: 2.5 oz
 Packaged: H 7.5 in. x W 4.5 in. x D 0.25 in.
 Unfolded: H 7.5 in. x W 9 in. x D 0.13 in.




online video available

< scan to view

HyFin[®] Vent Individual

HYFIN VENT INDIVIDUAL

ITEM #	NSN #
10-0029	6510-01-624-0840

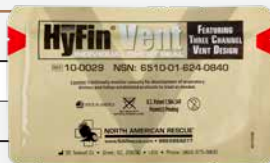

SPECIFICATIONS



Device: 6 in. x 6 in. (excludes tab) Wgt: 1.31 oz
 Packaged: H 4.63 in. x W 7.5 in. x D 0.1 in.

HyFin[®] Occlusive Individual

HYFIN INDIVIDUAL

ITEM #	NSN #
10-0015	6510-01-532-8019



SPECIFICATIONS



Device: 6 in. x 6 in. (excludes tab) Wgt: 1.31 oz
 Packaged: H 4.63 in. x W 7.5 in. x D 0.1 in.

USMC HyFin[®] Chest Seal COMBO PACK

- Versatile combo pack designed for the treatment of both entry/exit or multiple penetrating injuries to the chest
- Includes one HyFin[®] Occlusive & one HyFin[®] Vented Chest Seal packaged in a clear resealable pouch

SPECIFICATIONS

Packaged: L 9 in. x W 5 in. x D 0.2 in.
 Each Chest Seal (when deployed):
 H 6 in. x W 6 in. x D 0.1 in.
 Total Weight: 2.9 oz
 US Patent 7,504,549 &
 Patent(s) Pending



USMC COMBO PACK:



USMC HYFIN COMBO PACK

ITEM #	NSN #
85-1045	6510-01-654-6787





10 ga x 3.75 inch

SPEAR™ Simplified Pneumothorax Emergency Air Release

- Latest innovation in Needle Decompression technology
- Designed based on known gaps in the management of tension pneumothorax
- Enhanced needle tip for improved insertion control
- User friendly features like spin-lock connection, depth markings, and detachable low-pressure one-way valve
- Fenestrations along catheter for additional air flow
- Familiar, rugged, pen-style case for ease of storage and deployment
- Allows for either lateral decompression or anterior needle thoracostomy
- Catheter made of flexible kink-resistant material with a soft tip
- Single use only. Supplied sterile



online video available

< scan to view

SPECIFICATIONS

Packaged:
L 7.625 in. x Diameter 0.75 in.
Catheter: 10Ga x 3.75 in.
Weight: 1 oz
US Patents: Pending

SPEAR

GAUGE	ITEM #	NSN #
10 ga	10-0051	---



NOTE: Not for sale in Canada



Versatility that allows for both Lateral and Anterior Decompression Procedures...

DETACHABLE ONE-WAY VALVE


SPIN LOCK

10 gauge 3.75 in. FLEXIBLE CATHETER with DEPTH GAUGE


DISTAL FENESTRATIONS

SOFT SAFETY TIP

NEEDLE with MODIFIED TIP

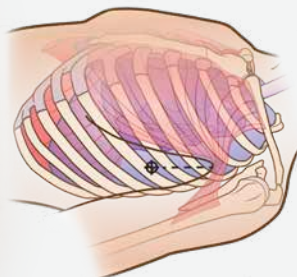

Lateral Approach

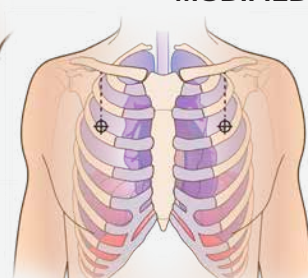

Anterior Approach

3.25 inch ARS® NEEDLE DECOMPRESSION KIT

- Penetrate the pleural space in 99% of casualties presenting the signs and symptoms of a tension pneumothorax
- Quality, reliable devices with 3.25 in. length needle/catheter
- Capless flash chamber provides immediate confirmation of correct placement
- Patented design packaged in an easy-to-open, crush-proof protective tube featuring easy-to-identify textured twist top with clip



14 ga

online video available



< scan to view

10 ga

SPECIFICATIONS

Pkg: H 5.7 in. x W 0.7 in. x D 0.8 in.
Weight: 0.6 oz

ARS NEEDLE DECOMPRESSION KIT

GAUGE	ITEM #	NSN #
14 ga	ZZ-0056	6515-01-541-0635
10 ga	ZZ-0298	6515-01-673-1701

U.S. Design Patents:
D584,409 S; D595,847S;
D584,410S; U.S. Patent
7,874,426B2; © North
American Rescue,®
LLC., 127820, 127822,
127823, 001002372-
0001, 001013940-0001,
001013940-0002



NAR SIMPLE THORACOSTOMY KIT



SPECIFICATIONS

Packaged: L 4.5 in. x W 7.5 in. x D 0.75 in.
Weight: 5.28 oz

- Compact kit that contains all necessary equipment to perform simple thoracostomy
- Includes 3-pack skin prep for performing the procedure bilaterally if needed
- HyFin Compact Vented Twin Pack Chest Seal to cover site if required


SIMPLE THORACOSTOMY KIT

ITEM #	UNIT	NSN#
83-0027	Single	-----

QTY 1 Hemostat 6.25" Rochester Pean Curved	QTY 1 Protected #10 Scalpel	QTY 2 Gauze Pads 4x4" (2pk)
QTY 1 HyFin® Vent Compact, Twin	QTY 1 Chloraprep Swabsticks, 3pk	QTY 1 Permanent Marker, Small

CHEST TUBE INSERTION KIT

- Designed to provide the necessary equipment needed to perform emergency chest tube thoracostomy
- Compact and lightweight kit packaged in an easy to open, vacuum sealed pouch
- Can be stored in your kit or crash cart

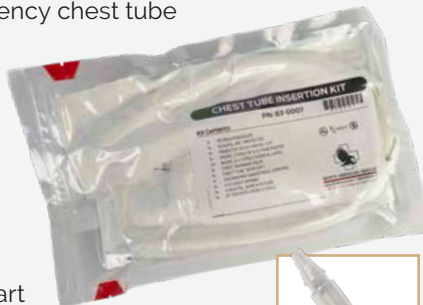

SPECIFICATIONS

Packaged: L 10 in. x W 6 in. x D 2.25 in.
Weight: 11.2 oz

CHEST TUBE INSERTION KIT


ITEM #	UNIT	NSN#
83-0007	Single	-----

QTY 1 Gauze, Petrolatum	QTY 1 Protected #10 Scalpel	QTY 1 5.5" Kelly Hemostat, Cvd	QTY 1 Drape, Fenestrated	QTY 1 Surgical Gloves, Latex
QTY 1 Valve Cook Chest Drain	QTY 1 36 French Chest Tube	QTY 1 Chloraprep Swabsticks, 3pk	QTY 1 4" x 4" Split Sponge	QTY 1 O-Silk Sutures, 30 in.
				QTY 1 8" Curved Roch-Peon



Meets all requirements for US Army Task: 081-831-1011

NAR SALINE LOCK KIT

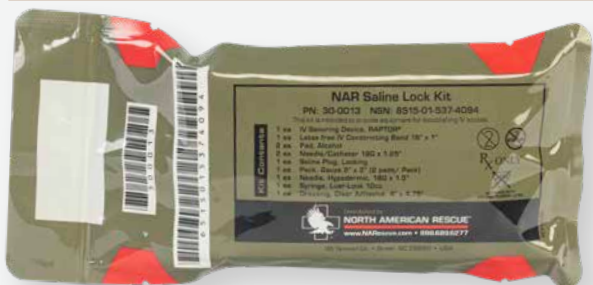

Essential components for establishing a saline lock for safer, quicker and more effective IV cannulation in the most demanding environments

NAR SALINE LOCK KIT

ITEM #	NSN #
30-0013	6515-01-537-4094

SPECIFICATIONS

Pkg: H 5.5 in. x W 3.25 in. x D 1 in.
Weight: 1.3 oz



SALINE LOCK KIT CONTENTS:



NAR NEEDLELESS SALINE LOCK KIT

IV safety catheter & needle-free valve port reduce or eliminate accidental needle-stick exposures

NAR NEEDLELESS SALINE LOCK

ITEM #	NSN #
30-0046	6515-01-607-1943

SPECIFICATIONS

Pkg: H 5.5 in. x W 3.25 in. x D 1 in.
Weight: 1.3 oz



NEEDLELESS SALINE LOCK KIT CONTENTS:



BOA® IV CONSTRICTING BAND

- Direct peripheral IV access utilizing a revolutionary rolling action for transfusing venous blood to the distal veins while still allowing arterial flow
- Quick Click connector & two-finger Quick Release improves "time on task"
- XL version stretches approximately 50% longer than the standard BOA®

SPECIFICATIONS

Standard: L 7.5 in. to 18 in. | XL: L 10 in. to 27 in.
Packaged: L 4 in. x W 3 in. x D 0.9 in. | Weight: 1.5 oz
US Patents: 5,304,202 & 5,607,448



online video available


< scan to view



BOA IV CONSTRICTING BAND



SIZE	ITEM #	NSN #
Standard	30-0009	6515-01-537-2611
XL	30-0071	----

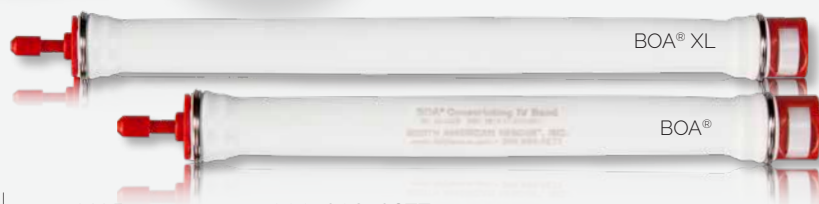

SHARPS SHUTTLE

SPECIFICATIONS

L 6.38 in. x Diameter 1.11 in.
Weight: 1.5 oz

SHARPS SHUTTLE

UNIT	ITEM #	NSN #
Single	30-0014	6515-01-518-8551
6 pack	30-00146	6515-01-619-5671



RAPTOR® IV Securing Device

SPECIFICATIONS

Extended: L 19 in. x W 1 in.
Weight: 19 oz

RAPTOR IV SECURING DEVICE CE

UNIT	ITEM #	NSN #
6 pack	30-0006	6515-01-523-3317



G3 POWER DRIVER by EZ-IO®

- Battery-powered solution for providing rapid vascular access to humerus & tibia in less than 10 seconds (Not for sternal use)
- Trigger guard prevents accidental activation
- Molded Hard Side Carrying Case (sold separately)

G3 POWER DRIVER

UNIT	ITEM #
Driver	30-0075
Case	30-0077

SPECIFICATIONS

H 4.5 in. x W 6.5 in. x D 2.5 in.
Weight: 15 oz



NOTE: EZ-IO Trainers also available - see page 108

MANUAL EXTREMITY NEEDLE SET


Manual I/O insertion at 6 extremity sites by EZ-IO®



SPECS

Device: L 6 in. x W 3.25 in. | Weight 1.4 oz

MANUAL NEEDLE SET:

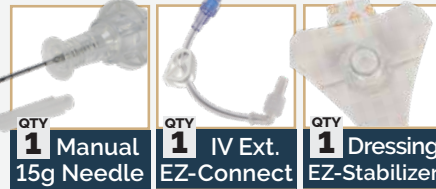

QTY 1 Manual 15g Needle
QTY 1 IV Ext. EZ-Connect
QTY 1 Dressing EZ-Stabilizer



QTY 1 Guide, I/O Reference
QTY 1 Sharps Block
QTY 1 EZ I/O Wristband

MANUAL NEEDLE

SIZE	ITEM #
38.5mm	30-0078

POWER DRIVER NEEDLE SET

for use with the EZ-IO® G3 Power Driver

POWER DRIVER NEEDLE SET

COLOR	SIZE	ITEM #
Pink	15mm/PD	30-0081
Blue	25mm/AD	30-0083
Yellow	45mm/LD	30-0085

SPECS

Packaged: H 4.75 in. W 4.25 in. D 0.7 in. Weight: 1 oz

POWER DRIVER NEEDLE SET:



QTY 1 Choice of Needle 15 ga
QTY 1 IV Ext. E-Z Connect
QTY 1 I/O Wristband
QTY 1 Sharps Block

T.A.L.O.N.™ MANUAL NEEDLE SET

Manual I/O insertion at 7 extremity sites - including the sternum

SPECIFICATIONS

Device: L 6 in. x W 3.25 in.
Weight 1.5 oz

T.A.L.O.N. MANUAL NEEDLE SET

SIZE	ITEM #
38.5mm	30-0132

T.A.L.O.N. MANUAL NEEDLE SET:




QTY 1 15g Manual Needle/Catheter
QTY 1 IV Ext. EZ-Connect
QTY 1 Sternal Locator
QTY 1 Guide, Quick Reference
QTY 1 Sharps Block
QTY 1 EZ I/O Wristband



CIRCULATION & Vascular Access

LIFEWARMER™ **Quantum™** IV Fluid & Blood Warming System

- **Lightweight** The Quantum™ has the lightest operational weight of any portable blood & fluid warming system (Total Operational Weight: 22 oz)
- **Powerful** Warms two units of 39°F refrigerated whole blood at 100mL/min+ on a single charge
- **Compact** Fewest components to carry (three)
- **Fast** Warm-up time of only 24 seconds
- **Intuitive** The thermal transfusion/infusion set looks, handles, and sets up/primes like any conventional IV set
- **Safe** Tested safe for blood, blood products & IV solutions
- **Effective** The system is plug-and-play for ease of setup under stressful conditions such as tactical field care (only two connection points, with no plug orientation or set-up order requirements – system status constantly visible at a glance)
- **Complete** Does not require additional standard transfusion or infusion set – all other warmers do
- **Tactical** Only system FDA cleared allowing easy disable of audible alerts to protect operators in combat

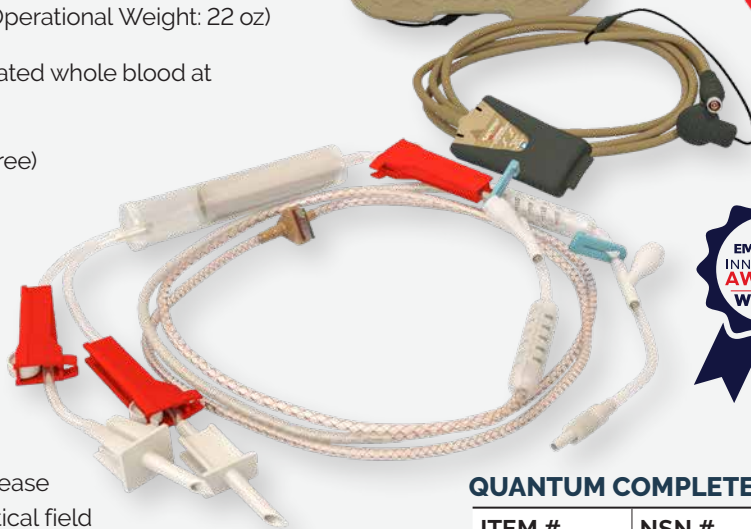

QUANTUM COMPLETE SYSTEM

ITEM #	NSN #
35-0006	6515-01-693-7810

OPERATIONAL CAPABILITIES

IV FLUIDS From a **single battery charge**, capable of warming **1700 mL** of **68°F IV fluids** to an output temperature of **100.4°F (38° +/- 2°C)**.
Warm up time: **24 seconds**.
Delivery Flow rate: **200mL/min**.

COMPLETE QUANTUM SYSTEM CONTENTS:

1 x Quantum Charger

1 x Quantum Battery

1 x Quantum Controller

2 x TTS-B Thermal Transfusion Set - Blood
(Blood, blood products & IV solutions)

SPECIFICATIONS

Total Operational Weight: 22 oz
(Quantum™ Controller, Battery & TTS-B)



Battery: 16 oz

Controller: 3 oz

Battery Charger: 22 oz

Controller Cover: 0.23 oz

Thermal Transfusion Set (TTS-B): 2.7 oz

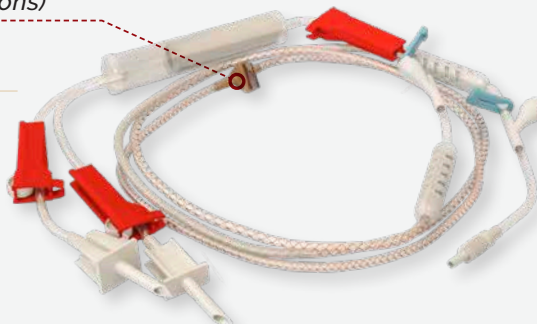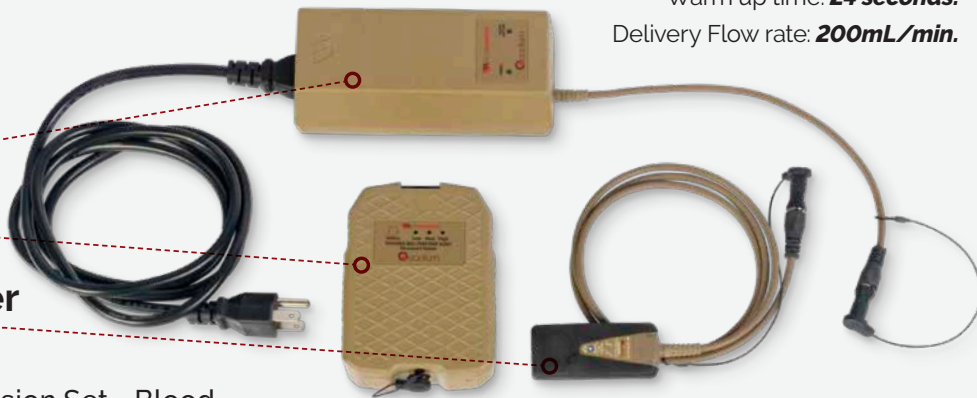

BLOOD or BLOOD PRODUCTS


From a **single battery charge**, capable of warming **2 units** of **39°F Whole Blood** to an output temperature of **100.4°F (38° +/- 2°C)**.
Warm up time: **24 seconds**.
Delivery Flow rate: **100mL/min**.

TTS-B Thermal Transfusion Set

- Compatible with whole blood, component therapy (including FDP), colloids, and isotonics
- Similar size and weight to traditional IV tubing
- Intelligent sensing & heating elements integrated into the walls of the tubing set
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller

SPECIFICATIONS

TTS-B package (each): 6.5 in. x 11 in.
 Tubing Length: 80 in.
 TTS-B Wgt (each): 2.6 oz
 Packaged Wgt (2 TTS-B & box): 11 oz

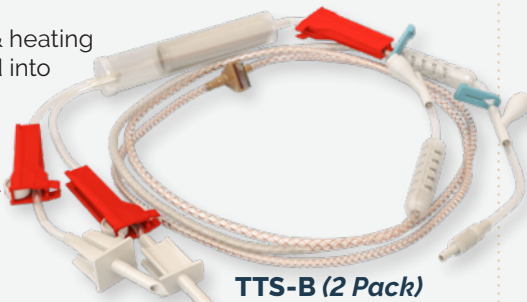

TTS-B (2 Pack)


ITEM #	NSN #
35-0005	6515-01-695-4444

TIS Thermal Infusion Set

- Similar size and weight to traditional IV tubing
- Intelligent sensing and heating elements integrated into the walls of the tubing set
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller
- Compatible with standard IV fluids such as colloids and isotonics

SPECIFICATIONS

TIS packaged (each): 6.5 in. x 11 in.
 Tubing Length: 80 in.
 Weight: TIS (each): 2.1 oz
 Packaged (2 TIS & box): 10 oz

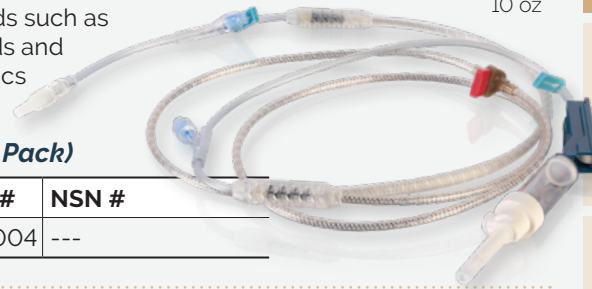

TIS (2 Pack)

ITEM #	NSN #
35-0004	---



Quantum BATTERY

- Lightweight at only 16 oz
- Designed for the tactical environment- water resistant, subdued color, and the option to disable the audible alarm
- 44V Lithium-Polymer Battery
- Intelligent, state-of-the-art battery technology
- Status button to check charge state at any time

QUANTUM BATTERY

ITEM #	NSN #
35-0001	---

SPECIFICATIONS

H 4.5 in. x W 3.25 in. x D 1.25 in.
 Weight: 16 oz

Quantum HARDCASE



SPECS

L 10.63 in.
 x W 9.63 in.
 x D 6.73 in.
 Weight: 3.22 lb

- Rugged hard case designed to protect your Quantum from dust, dirt, and moisture
- Features pressure relief valve for automatic pressure equalization
- Custom laser cut high density foam holds components securely in place
- Flexible storage capacity for multiple batteries, controllers and thermal tubing sets

QUANTUM HARDCASE

ITEM #	NSN #
35-0008	---

Quantum CONTROLLER

- Moderates the energy flow from the Quantum Battery to the Thermal Tubing
- Designed for ease of connecting and durability
- Colored LED lights constantly show during administration of blood/fluid to show status
- Lightweight & durable
- IP53 rated for dust & water resistance with jacket



Replacement cover is also available

- Minimum service life - 1000 insertions/removals

SPECIFICATIONS


Controller Body:
 W 1.25 in. x H 3 in. x D 0.75 in.
 Controller Total Length: 51.25 in.
 Weight (Controller only): 2.5 oz
 Packaged Weight: 5 oz

QUANTUM CONTROLLER

OPTION	ITEM #	NSN #
Controller	35-0003	---
Cover (Replacement)	35-0010	---

- Fully charge the Quantum Battery in <90 minutes
- Smart, safe charging
- Includes 6 ft AC power cord
- Colored LED lights to inform the user of charging and battery status
- IP 22 rated for dust and water resistance
- 100V to 264V AC input; 50Hz, 60Hz
- Minimum service life- 2500 insertions/removals

Quantum CHARGER

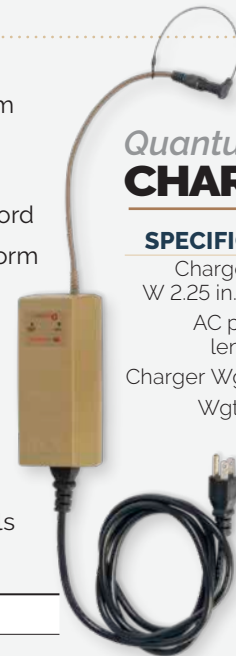

SPECIFICATIONS

Charger: H 6 in. x W 2.25 in. x D 1.75 in.
 AC power cord length: 6 feet
 Charger Wgt: 1 lb 4 oz
 Wgt with cord: 1 lb 10 oz

QUANTUM CHARGER

ITEM #	NSN #
35-0002	---



Quantum™ Field Blood Transfusion Kits

- Comprehensive kit designed for use with the Quantum Blood & Fluid Warming System allowing rapid, consistent normothermic field blood transfusion in far forward and near point of injury (POI) prehospital casualty management.
- Includes the Quantum Thermal Transfusion Set - (Blood) that incorporates an integrated heating system that does not disrupt the fluid path
- Equipped with Donor/Recipient Modules to verify blood types of the donor and recipient
- Also allows collecting of blood samples from the donor for additional surveillance testing

QUANTUM FIELD BLOOD TRANSFUSION KIT

ITEM #	NSN #
83-0028	---

SPECIFICATIONS

L 10 in. x W 7 in. x D 2.5 in.
Weight: 1.09 lb

QUANTUM FIELD BLOOD TRANSFUSION KIT CONTENTS:

 QTY 1 Single Blood Bag	 QTY 1 55in. Curved Kelly Hemostat	 QTY 1 9.5 in. 550 Cord	 QTY 1 Surgical Tape, 1 in.	 QTY 1 Quantum TTS-B Set
 QTY 1 Permanent Marker, Large	 QTY 2 Large Black Nitrile Gloves	 QTY 1 Beaded Cord, 11 in.	 QTY 1 DD572 Questionnaire	 QTY 1 Donor Module Kit
				 QTY 1 Recipient Module Kit

*Each Donor Module Kit contains:

 QTY 5 Alcohol Prep Pads	 QTY 2 Gauze Pads 4x4 (2pk)	 QTY 1 Constricting IV Band	 QTY 1 Adhesive Dressing, Clear	 QTY 1 Blood Bag Label	 QTY 1 Eldon Blood Type Kit
 QTY 1 Injection Site Adaptor	 QTY 1 18ga Needle	 QTY 1 10 cc Syringe	 QTY 1 Device Blood Transfer	 QTY 2 Vacutainer 5 ml Gold top	 QTY 2 Vacutainer 5 ml Pearl top
				 QTY 1 Vacutainer 5 ml Purple top	

**Each Recipient Module Kit contains:

 QTY 5 Alcohol Prep Pads	 QTY 2 Gauze Pads 4x4 (2pk)	 QTY 1 Constricting IV Band	 QTY 1 Adhesive Dressing, Clear	 QTY 1 18ga IV Catheter, 1.25 in.	 QTY 1 16ga IV Catheter, 1.25 in.	 QTY 1 3 Way Stopcock	 QTY 1 Injection Site Adaptor	 QTY 1 Eldon Blood Type Kit
--	---	---	---	---	---	--	---	---



All the needed supplies to allow the user to rapidly perform a field blood transfusion as outlined in the

Joint Theater Trauma System Clinical Practice Guidelines



FBTK Field Blood Transfusion Kits

- Comprehensive kit that provides all the needed supplies to allow the user to rapidly perform a field blood transfusion
- Equipped to allow for verifying the blood types of the donor and recipient and collecting blood samples from the donor for additional testing

FBTK

OPTIONS	ITEM #
FBTK Basic	83-0022
Blood Surveillance	83-0025
FBTK Complete*	83-0000

FBTK Basic:

- Donor Module:**
- 5 x Alcohol Prep Pad
 - 2 x 2-pack Gauze Pad - 4 in. x 4 in.
 - 1 x IV Constricting Band - 1 in. x 18 in.
 - 1 x Clear Adhesive Dressing - 4 in. x 4.75 in.
 - 1 x Blood Bag Label

Recipient Module:

- 5 x Alcohol Prep Pad
- 1 x IV Catheter - 18 ga x 1.25 in.
- 2 x 2-pack Gauze Pad - 4 in. x 4 in.
- 1 x IV Constricting Band - 1 in. x 18 in.
- 1 x Clear Adhesive Dressing - 4 in. x 4.75 in.
- 1 x IV Catheter - 16 ga x 1.25 in.

Additional Supplies:

- 1 x roll 1 in. Surgical Tape
- 2 x pair Black Nitrile Gloves, Large
- 1 x Black Permanent Marker
- 1 x Kelly Curved Hemostat - 5.5 in.
- 1 x Single Blood Bag (450ml)
- 1 x 550 Cord - 9.5 in.
- 1 x Beaded Security Tie - 11 in.
- 1 x IV Set Y-Type with Filter 170-260 UM
- 1 x Questionnaire Card - DD572 Modified

FBTK Blood Surveillance Module:

- Donor Module:**
- 1 x IV Hypodermic Needle - 18 ga x 1.5 in.
 - 1 x IV Syringe, 10cc
 - 1 x Eldon Card
 - 1 x Interlink Injection Site Needle
 - 1 x Blood Transfer Device
 - 2 x Vacutainer, Gold Top (5ml)
 - 1 x Vacutainer, Purple Top (5ml)
 - 2 x Vacutainer, Pearl Top (5ml)

Recipient Module:

- 1 x 3 Way Stopcock
- 1 x Eldon Card
- 1 x Interlink Injection Site Needle

*FBTK Complete:

Contains all components of *FBTK Basic* & *Blood Surveillance Module*



LifeFlow[®] FLUID INFUSER

- Hand-operated, rapid infusion device that allows health care providers to deliver fluids quickly and effectively
- Deliver 500ml of fluid in less than 2 minutes through a 20G IV catheter
- Designed to protect the syringe from bacterial contamination
- Can be operated with one hand


SPECIFICATIONS

12 in. x 7.5 in. x 2.75 in.
Weight: 0.95 lb

LIFEFLOW INFUSER

ITEM #	NSN #
30-0286	

LifeFlow[®] VOLUME RESUSCITATION. WHEN MINUTES MATTER.

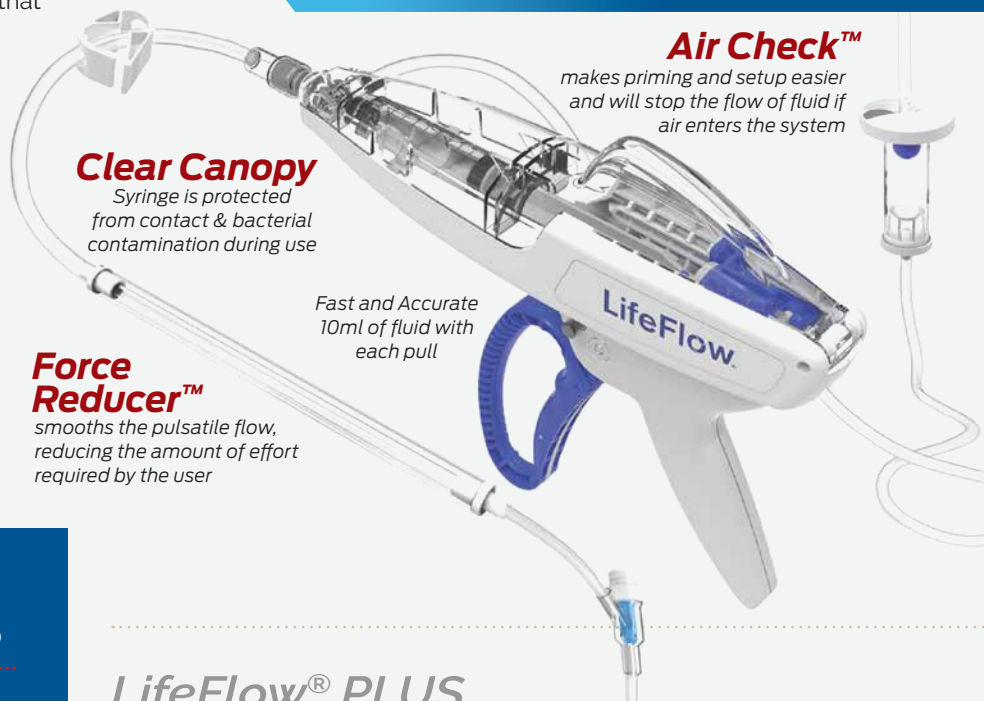

Air Check[™]

makes priming and setup easier and will stop the flow of fluid if air enters the system

Clear Canopy

Syringe is protected from contact & bacterial contamination during use

Fast and Accurate
10ml of fluid with each pull

Force Reducer[™]

smooths the pulsatile flow, reducing the amount of effort required by the user

Improve Resuscitation Workflow with LifeFlow[®]

> RESUSCITATE

500 ml in < 2 minutes
Stay with Your Patient

Deliver first bolus quickly to reverse hypotension in patients with shock

- Provide active resuscitation
- Maintain control over delivery

> RESTORE

Improve Perfusion / Reverse Shock

Assess using vitals & bedside markers of perfusion, decide if additional fluid bolus is needed.

Early fluids more rapidly reverse shock & **improve outcomes**.¹⁻³

> REASSESS

Rapid Decision / Disposition


Each bolus delivered provides an opportunity to evaluate and determine next steps.

Faster resuscitation may lead to:

- Less mechanical ventilation⁴
- More rapid disposition
- **Less donning/doffing of PPE**

LifeFlow[®] PLUS BLOOD & FLUID INFUSER

- Hand-operated, rapid infusion device that allows health care providers to deliver fluid & blood products quickly and effectively
- Deliver 500ml of fluid in less than 2 minutes through a 20G IV catheter
- Deliver 1 unit of whole blood in 2 minutes
- Goes from the field to the hospital
- Can be operated with one hand
- Handle and tubing packaged together for quick and easy set-up
- Designed to protect the syringe from bacterial contamination
- Average time to set up and prime: less than 3 minutes

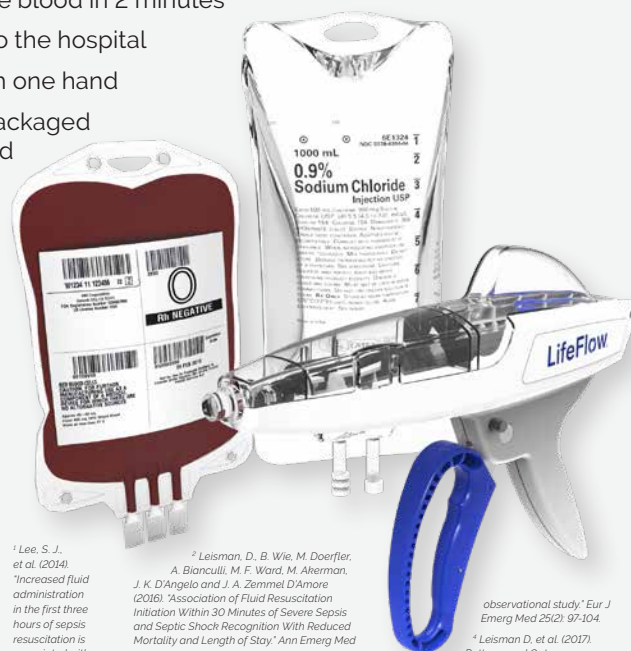

SPECIFICATIONS

12 in. x 7.5 in. x 2.75 in.
Weight: 0.95 lb

LIFEFLOW PLUS INFUSER

ITEM #	NSN #
30-0287	

¹ Lee, S. J. et al. (2014). "Increased fluid administration in the first three hours of sepsis resuscitation is associated with reduced mortality: a retrospective cohort study." Chest 146(4): 908-915.

² Leisman, D. B. Wie, M. Doerfler, A. Bianculli, M. F. Ward, M. Akerman, J. K. D'Angelo and J. A. Zemmel D'Amore (2016). "Association of Fluid Resuscitation Initiation Within 30 Minutes of Severe Sepsis and Septic Shock Recognition With Reduced Mortality and Length of Stay." Ann Emerg Med 68(3): 298-311.

³ Williams, J. M., J. H. Greenslade, C. A. Dymond, K. Chu, A. F. T. Brown and J. Lipman (2018). "Characteristics, treatment and outcomes for all emergency department patients fulfilling criteria for septic shock: a prospective

observational study." Eur J Emerg Med 25(2): 97-104.

⁴ Leisman D. et al. (2017). Patterns and Outcomes Associated With Timeliness of Initial Crystalloid Resuscitations in a Prospective Sepsis and Septic Shock Cohort. Crit Care Med. Oct.45(10):1596-1606.



HYPOTHERMIA Prevention & Management

HPMK[®] Hypothermia Prevention & Management Kit

- Includes reinforced Heat Reflective Shell (HRS™) of 4-ply, composite fabric with a protected non-conductive reflective layer
- Built-in hood, integrated fluid absorption pad & top-to-bottom taper providing excellent thermal insulation
- Continuous 1.5 in. Velcro® closures provide rapid 360° access to the casualty with minimal exposure to the elements
- Self-heating, 4-cell liner (Ready-Heat™ Blanket) provides 10 hours of continuous dry heat up to 125°F (52°C)
- Vacuum-packaged for a low-cube footprint



HPMK[®] CONTENTS:



HPMK

ITEM #	NSN #
80-0027	6515-01-532-8056

NOTE: Litter is sold separately

NOTE: Not for sale in Canada

SPECIFICATIONS



Shell Pkg: H 6.75 in. x W 10.5 in. x D 5.5 in.
 Shell Open: L 78 in. x W 43 in.
 US Patent: 7,766,950 B2
 Weight: 3 lb 8 oz

HRS[®] Heat Reflective Shell



As an individual unit, reinforced 4-ply, composite shell with built-in hood, weight tested to support 250+ lb.

HEAT REFLECTIVE SHELL

ITEM #	NSN #
ZZ-0051	6532-01-627-7888

SPECIFICATIONS



Shell Packaged: L 10 in. x W 6.75 in. x D 2.5 in.
 Open: L 78 in. x W 43 in. US Patent: 7,766,950 B2
 Weight: 1 lb 13 oz

READY HEAT™ BLANKET



As an individual unit, the self-heating, 4-cell liner provides 10 hours of continuous dry heat after reaching peak temperature approximately 15-20 minutes after opening

READY HEAT BLANKET

ITEM #
ZZ-0052

SPECIFICATIONS



Open: L 48 in. x W 36 in.
 Weight: 1 lb 11 oz

EMERGENCY SURVIVAL WRAP

SURVIVAL WRAP

COLOR	ITEM #
Orange/Silver	ZZ-0834
OD Green/Silver	ZZ-1390

- Designed to reflect heat back to the body
- Reversible Silver / high visibility Orange or OD Green

SPECIFICATIONS

Packaged: L 5 in. x W 3.75 in.
 Open Wrap: L 96 in. x W 60 in. | Weight: 3.7 oz



INTRO

M.A.R.C.H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX





KITS

COMBAT CASUALTY RESPONSE (CCRK®)


- Warfighter/Operators • Combat Lifesavers
- 1st Responders • Tactical Healthcare Professionals
- DoD Bleeding Control Response Kits **pages 24-52**



CUSTOMIZABLE SOLUTIONS

TRAINING

INDEX

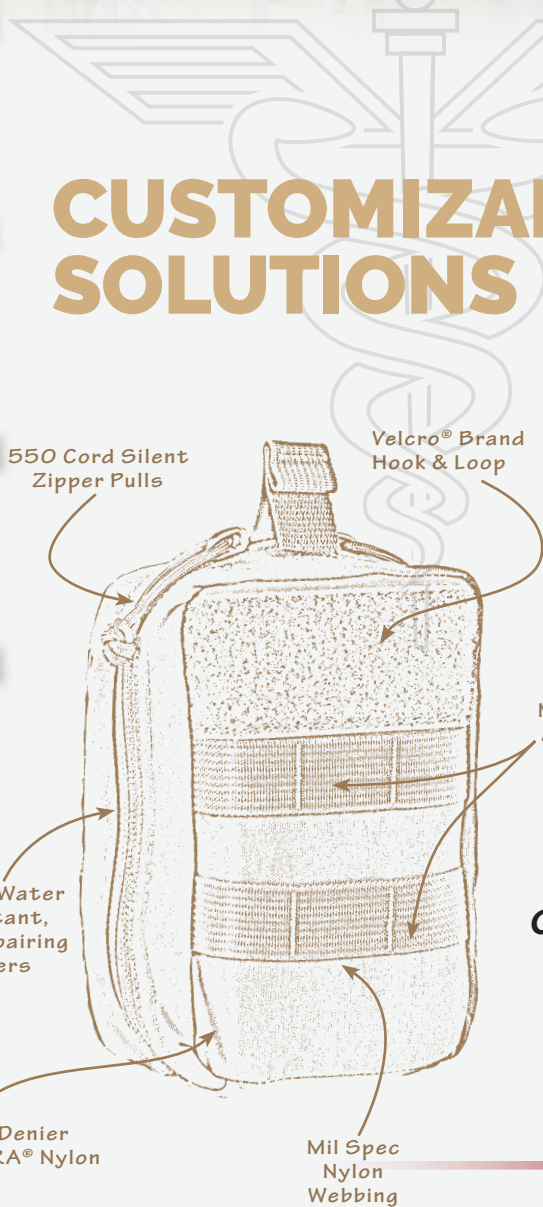


550 Cord Silent Zipper Pulls

Velcro® Brand Hook & Loop

Sand & Water Resistant, Self-Repairing Zippers

500 Denier CORDURA® Nylon

Mil Spec Nylon Webbing

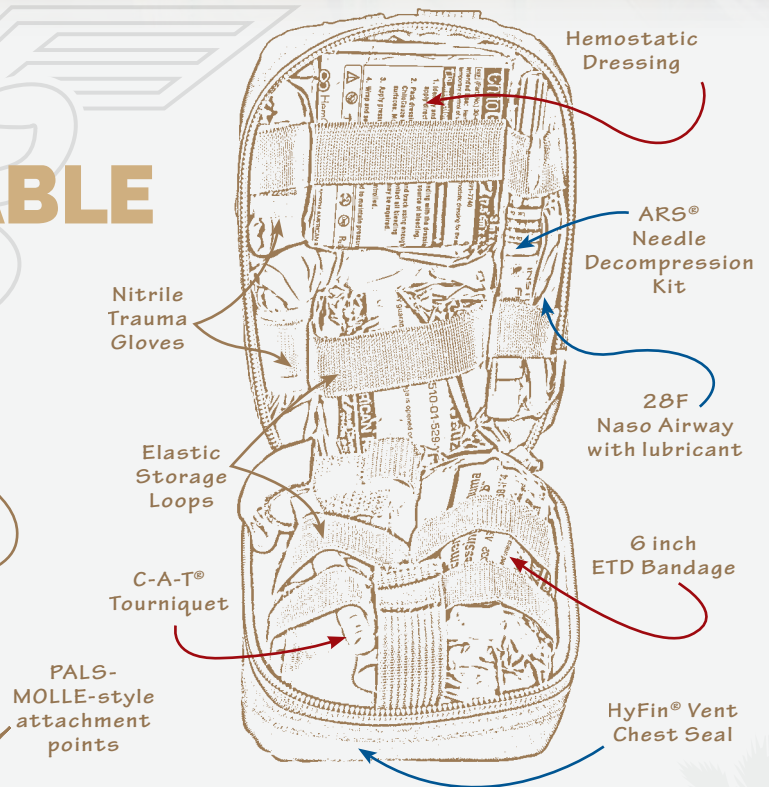

Nitrile Trauma Gloves

Elastic Storage Loops

C-A-T® Tourniquet

PALS-MOLLE-style attachment points

Hemostatic Dressing

ARS® Needle Decompression Kit

28F Naso Airway with lubricant

6 inch ETD Bandage

HyFin® Vent Chest Seal

...fulfilling YOUR complete Medical & Material logistics requirements!



North American Rescue® is **THE market leader** for designing, producing, integrating and assembling tactical field medical sets and kits **supporting the services of the Department of Defense, federal agencies** as well as **state and local law enforcement**.

Over the past several years, our ability to provide rapid response products and services has aided in the delivery of critical battlefield medical care when and where it's needed most with an **end result of saving more lives**.

A comprehensive list of supply chain services affords a unique opportunity to North American Rescue's customers in the procurement, assembly, storage and management of their **complete medical and material logistics** requirements. NAR maintains a 75,000 square foot warehouse facility which consists of 50,000 square feet for receiving and distribution staging and an additional 14,000 square feet solely dedicated to the assembly room where sets, kits and **custom medical solutions are built to order**.

Highly qualified logistics and warehouse staff ensure that orders, inventory and distribution processes **maintain the superior quality standards** associated with the North American Rescue brand. Our integral partnership with the South Carolina Vocational Rehabilitation **NISH Partner** is a pillar in the successful execution of these processes. ■



SIRK Supplemental IFAK ReSupply Kits

KIT CAPABILITIES INCLUDE:



SIRK Gen 2

SIRK GEN 2:



SPECIFICATIONS

H 7 in. x W 5 in. x D 2 in.
Weight: 14 oz

SIRK GEN 2

ARS OPTION	ITEM#
w/o ARS Needle	85-1738
with ARS Needle	85-1899

- Vacuum-sealed first aid kits ideal for re-stocking an existing IFAK that contains the essential items to treat leading causes of preventable combat death
- Designed for small cube space in rugged, clear package allowing contents to be easily visualized

SIRK Gen 1

SIRK GEN 1:



SPECIFICATIONS

H 7 in. x W 5 in. x D 2 in.
Weight: 14 oz

SIRK GEN 1

HEMOSTATIC	ITEM#
Combat Gauze	80-0291
ChitoGauze	80-0225**

IPOK Individual Patrol Officer Kit

KIT CAPABILITIES INCLUDE:



- Individual hemorrhage control kit small enough to be carried on your person
- Compact, vacuum-sealed & packaged to fit into a BDU pocket, vest pouch or individual bag
- Rugged, durable pouch features NAR's signature Red Tip Technology® ensuring quality & performance

IPOK:



IPOK

HEMOSTATIC/GAUZE	ITEM#
S-Rolled Gauze	80-0167**
Combat Gauze, LE	80-0168



SPECIFICATIONS

H 4 in. x W 6.5 in. x D 2.75 in.
Weight: 7.4 oz

M-FAK *Mini First Aid Kit*

KIT CAPABILITIES INCLUDE:



- Contains the essential point-of-wounding equipment needed for treating penetrating, blast and other traumatic injuries
- Super compact, rugged nylon platform that allows attachment both vertically (MOLLE backing) and horizontally (3 in. belt loop)
- Clam shell configuration utilizes (2) two main sleeves that open on both ends for easy access

M-FAK MINI FIRST AID KIT

COLOR	LEVEL	ITEM #
Black	Basic	80-0494
Coyote	Basic	80-0495
Multi-Cam	Basic	80-0496
OD Green	Basic	80-0497
Ranger Green	Basic	80-0996
Black	Combat Gauze LE	80-0498
Coyote	Combat Gauze LE	80-0499
Multi-Cam	Combat Gauze LE	80-0500
OD Green	Combat Gauze LE	80-0501
Ranger Green	Combat Gauze LE	80-0998
Black	Advanced	80-0502
Coyote	Advanced	80-0503
Multi-Cam	Advanced	80-0504
OD Green	Advanced	80-0505
Ranger Green	Advanced	80-0997



SPECIFICATIONS


H 6 in. x W 3 in. x D 3 in.
Weight: 13 oz



M-FAK BASIC:



M-FAK with COMBAT GAUZE:

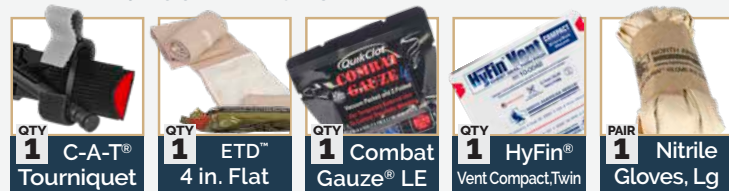

M-FAK ADVANCED:



INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

ROO M-FAK *Mini First Aid Kit*

KIT CAPABILITIES INCLUDE:



- Contains the essential point-of-wounding equipment needed for treating penetrating, blast and other traumatic injuries
- Unique front pouch that easily opens with the pull of a tab for rapid access to your tourniquet
- Super compact, rugged nylon platform that allows attachment both vertically (MOLLE backing) and horizontally (3 in. belt loop)
- Clam shell configuration utilizes (2) two main sleeves that open on both ends for easy access

Rapid one handed access to your tourniquet without ever having to open the full kit!



SPECIFICATIONS

H 6 in. x W 3 in. x D 3 in. | Weight: 13 oz

ROO M-FAK MINI FIRST AID KIT

COLOR	LEVEL	ITEM #
Black	Basic	80-1049
Coyote	Basic	80-1050
OD Green	Basic	80-1051
Ranger Green	Basic	80-1052
Black	Combat Gauze LE	80-1053
Coyote	Combat Gauze LE	80-1054
OD Green	Combat Gauze LE	80-1055
Ranger Green	Combat Gauze LE	80-1056
Black	Advanced	80-1057
Coyote	Advanced	80-1058
OD Green	Advanced	80-1059
Ranger Green	Advanced	80-1060

Bag only option also available at NARescue.com

ROO M-FAK BASIC:



ROO M-FAK with COMBAT GAUZE:



ROO M-FAK ADVANCED:



▲ NEW! ROO bag includes a Mini 1 in. PVC MED ID Patch



E-FAK Expeditionary First Aid Kit

KIT CAPABILITIES INCLUDE:



- Rugged nylon platform designed to meet the Expeditionary Military Forces First Responder requirements for critical point-of-wounding medical supplies
- Includes an additional TQ Holder and Tourniquet for primary placement separate from the EFAK, meeting the (2) Tourniquet Expeditionary Military requirement
- Allows attachment both vertically and horizontally
- Advanced Treatment pack out

SPECIFICATIONS

L 6 in. x W 3 in. x D 3 in.
Weight: 1 lb 3 oz

E-FAK

COLOR	ITEM #
Coyote	85-2533
Black	80-1044



OCHO IFAK

Compact and Low Visibility IFAK

KIT CAPABILITIES INCLUDE:



- Compact, low visibility, point-of-wounding IFAK designed for your duty/pistol belt
- Equipped with a 3 in. belt loop attachment to accommodate a duty/pistol belt
- Constructed with rip stop nylon, elastic loop organizational storage and easy open zipper
- Purposely designed to look like a non-tactical, accessory bag



Every OCHO IFAK includes these core items:



SPECIFICATIONS

L 7 in. x H 2.25 in.
x D 3.25 in.
Weight: 13 oz

OCHO IFAK

COLOR	LEVEL	ITEM #
Black/Gray	Basic	80-0374
Blue/Gray	Basic	80-0375
Red/Gray	Basic	80-0915
Black/Gray	Combat Gauze	80-0921
Blue/Gray	Combat Gauze	80-0922
Red/Gray	Combat Gauze	80-0917
Black/Gray	Advanced	80-0372
Blue/Gray	Advanced	80-0373
Red/Gray	Advanced	80-0916

PLUS+ BASIC

also includes:




COMBAT GAUZE

also includes:



ADVANCED also includes:

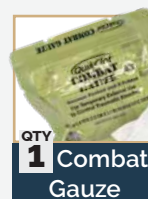

See OCHO Accessory Bag only on page 56

EAGLE IFAK

Compact and versatile Individual First Aid Kit

KIT CAPABILITIES INCLUDE:



- Can connect to the assault vest/ body armor or to a duty/pistol belt
- Contains the essential equipment needed for self-aid/buddy-aid to address the leading causes of potentially preventable death in a tactical environment
- Quick Pull Tab easily opens the kit with just one hand
- Unique, detachable shock-cord inner panel that can be adjusted to your components
- Inner panel is removable for easy access to your kit



Bag only option also available at NARescue.com

SPECIFICATIONS
H 8 in. x W 4 in. x D 3 in.
Weight: 1 lb 2 oz

EAGLE IFAK

COLOR	LEVEL	HEMOSTATIC	ITEM #
Black	Basic	None	85-0177
Coyote	Basic	None	85-0178
Multi-Cam	Basic	None	85-0179
OD Green	Basic	None	85-0180
Black	Advanced	Combat Gauze	85-0189
Coyote	Advanced	Combat Gauze	85-0190
Multi-Cam	Advanced	Combat Gauze	85-0191
OD Green	Advanced	Combat Gauze	85-0192
4 options online	Advanced	ChitoGauze	See NARescue.com

EAGLE IFAK BASIC:



EAGLE IFAK ADVANCED:



KIT CAPABILITIES INCLUDE:



RIG Series EAGLE IFAK

New, low profile IFAK with rapid access to your C-A-T® and other equipment



- Tourniquet attached for quick deployment by banding it to inside of front opening flap
- Protective hood with quick pull latch release helps protect IFAK components and includes sewn in embroidered medical cross identifier
- Detachable (hook & loop backed) shock-cord inner panel adjusts to secure internal components and is easily removable for rapid access

RIG EAGLE IFAK

COLOR	LEVEL	ITEM #
Black	Basic	80-0509
Coyote	Basic	80-0510
Multi-Cam	Basic	80-0511
Ranger Green	Basic	80-0991
Black	Combat Gauze LE	80-0517
Coyote	Combat Gauze LE	80-0518
Multi-Cam	Combat Gauze LE	80-0519
Ranger Green	Combat Gauze LE	80-0993
Black	Advanced	80-0513
Coyote	Advanced	80-0514
Multi-Cam	Advanced	80-0515
Ranger Green	Advanced	80-0992

SPECIFICATIONS

H 8 in. x W 4 in. x D 3 in.
 Kit Packed Weight: 118 lb
 Empty Bag Weight: 7 oz

Bag only option also available at NARescue.com

RIG SERIES EAGLE BASIC:



RIG SERIES EAGLE with COMBAT GAUZE LE:



RIG SERIES EAGLE ADVANCED:



INTRO
 M-A-R-C-H
 KITS
 T-A-C-E-V-A-C
 M-E-D-S-U-R-G
 R-E-S-C-U-E
 T-R-E-A-T-M-E-N-T
 S-U-P-P-O-R-T
 T-R-A-I-N-I-N-G
 I-N-D-E-X

KIT CAPABILITIES INCLUDE:



Tactical Operator Response Kit **T.O.R.K.**



- Streamlined, compact IFAK designed to MOLLE to the assault vest or body armor
- Includes essential life-saving equipment to be used at the point of wounding
- For self-aid/buddy-aid injuries caused by penetrating or blast trauma

SPECIFICATIONS

H 7 in. x W 4.5 in. x D 2.75 in.
Weight: 1 lb 5 oz

*Bag only options
also available at
NARescue.com*

T.O.R.K. TACTICAL OPERATOR RESPONSE KIT

COLOR	LEVEL/HEMOSTATIC	ITEM #
Black	Basic	85-0008
Coyote	Basic	85-0009
Digital Camo	Basic	85-2300
Multi-Cam	Basic	85-0012
OD Green	Basic	85-0054
Red	Basic	85-1805
Ranger Green	Basic	85-2301
Black	Basic w/Combat Gauze	80-1097
Coyote	Basic w/Combat Gauze	80-1098
Digital Camo	Basic w/Combat Gauze	80-1099
Multi-Cam	Basic w/Combat Gauze	80-1100
OD Green	Basic w/Combat Gauze	80-1101
Ranger Green	Basic w/Combat Gauze	80-1102
Red	Basic w/Combat Gauze	80-1108
Black	Advanced w/Combat Gauze	85-0048
Coyote	Advanced w/Combat Gauze	85-0049
Digital Camo	Advanced w/Combat Gauze	85-0050
Multi-Cam	Advanced w/Combat Gauze	85-0052
OD Green	Advanced w/Combat Gauze	85-0053
Ranger Green	Advanced w/Combat Gauze	85-2302
Red	Advanced w/Combat Gauze	85-1806
7 options online	Advanced w/ChitoGauze	NARescue.com



T.O.R.K. BASIC / BASIC with Combat Gauze:



T.O.R.K. ADVANCED:



USCG Mini IFAK

Designed to support the operations of the USCG Boarding Teams



- Improved first aid kit designed for USCG Boarding Teams providing for safe first aid with a minimal footprint design
- Featuring our modular system that can be worn by the operator, hooked directly to the duty belt, or attached to an assault vest or body armor with MOLLE-style connectors

SPECIFICATIONS

H 6.5 in. x W 3.5 in. x D 3.5 in. | Weight: 1 lb

USCG Mini IFAK:

 QTY 1 C-A-T® Tourniquet	 QTY 1 Combat Gauze Z-fold	 QTY 1 ETD™ 4 in. Flat	 QTY 1 HyFin® Occlusive, Sgl.
 PAIR 2 Nitrile Gloves, Lg	 QTY 1 Mini Duct Tape, 2 in.	 QTY 1 Model 7 Rescue Hook	 QTY 1 Permanent Marker, Large



USCG MINI IFAK

ITEM #	NSN #
80-1024	6545-01-682-2629

Bag only option
also available
at NARescue.com

USCG IFAK-Basic

Low Profile Basic IFAK designed for the US Coast Guard Boarding Teams



SPECIFICATIONS

H 8.25 in. x W 7.1 in. x D 5.5 in
Weight: 2 lb 10 oz

Versatile modular system can be worn by the operator as a drop-down leg rig, duty-belt add-on or assault vest/body armor attachment (via MOLLE-style connectors)

USCG IFAK BASIC

ITEM #	NSN #
80-1023	6545-01-682-2596

USCG IFAK Basic:

 QTY 1 C-A-T® Tourniquet	 QTY 1 Combat Gauze Z-fold	 QTY 1 ETD™ 4 in. Flat	 QTY 1 HyFin® Occlusive, Sgl.
 PAIR 2 Nitrile Gloves, Lg	 QTY 1 Mini Duct Tape, 2 in.	 QTY 1 Model 7 Rescue Hook	 QTY 1 Permanent Marker, Large

Bag only option also available at NARescue.com



USCG IFAK

Meets specific requirements of the US Coast Guard

KIT CAPABILITIES INCLUDE:



Versatile modular system can be worn by the operator as a drop-down leg rig, duty-belt add-on or assault vest/body armor attachment (via MOLLE-style connectors)

USCG IFAK

COLOR	ITEM #	NSN #
Coyote	80-0158	6545-01-587-8513
Black	80-0149	6545-01-572-7248
Resupply*	85-0104	----

USCG IFAK:



SPECIFICATIONS

H 8.25 in. x W 7.1 in. x D 5.5 in
Weight: 2 lb 10 oz

Bag only option also available at NARescue.com

*Resupply does not include the Nylon IFAK Bag



MARITIME ASSAULT KIT

Meets specific requirements of the Navy Expeditionary Combat Command

KIT CAPABILITIES INCLUDE:



Designed to meet the unique & specific requirements of the Navy Expeditionary Combat Command for a small, efficient IFAK

SPECIFICATIONS

H 7.4 in. x W 6.9 in. x D 5 in
Weight: 1 lb 13 oz

MARITIME ASSAULT KIT

COLOR	ITEM #	NSN #
Multi-Cam	80-0245	----
Coyote	80-0015	6545-01-529-4187
OD Green	80-0014	6545-01-531-6626
Black	80-0013	6545-01-529-4185
Digital Camo	80-0058	----

MARITIME ASSAULT KIT:



Bag only option also available at NARescue.com

INDIVIDUAL KIT

Designed for the Warfighter/Operator

KIT CAPABILITIES INCLUDE:



Designed with all the tools for the Warfighter/Operator to mitigate extremity bleeding, tension pneumothorax & airway obstruction

SPECIFICATIONS

H 8.5 in. x W 8.75 in. x D 4.5 in. Weight: 1 lb 12 oz

INDIVIDUAL KIT

COLOR	ITEM#	NSN#
Coyote	80-0005	----
OD Green	80-0004	----
Black	80-0003	6530-01-529-1134
Digital Camo	80-0091	----

Bag only option also available at NARescue.com



INDIVIDUAL KIT:



MEDIC LEG RIG

Rapid Access and Organization of Critical Equipment

KIT CAPABILITIES INCLUDE:



Designed for Tactical Healthcare Professionals configurable based on individual skill set & mission requirement

SPECIFICATIONS

H 9.5 in. x W 9 in. x D 7.75 in. Weight: 3 lb 9 oz

Bag only option also available at NARescue.com

MEDIC LEG RIG

COLOR	ITEM#	NSN#
Multi-Cam	80-0249	----
Coyote	80-0048	6545-01-543-3170
OD Green	80-0047	6545-01-543-3174
Black	80-0046	6545-01-543-3185
Digital Camo	80-0087	----



MEDIC LEG RIG:





SQUAD KIT

Enhance the survivability of the Assault Squad

KIT CAPABILITIES INCLUDE:



- Contoured profile, large central compartment, multiple attachment points & front trauma-shear storage
- Multiple-load carriage capabilities (fanny pack, over the shoulder, attachments for 3-day assault packs & MOLLE-style webbing)

Bag only option also available at NARescue.com

SQUAD KIT:



QTY 2 C-A-T® Tourniquet



QTY 4 ETD™ 6 in. Dressing



QTY 1 ETD™ Abdominal



QTY 4 Gauze S-Rolled



QTY 2 28F Naso Airway & Lube



QTY 1 Gauze, Petrolatum



QTY 1 HyFin® Occlusive, Sgl.



QTY 2 ARS® Needle, 14g



QTY 2 SAM® Splint II




PAIR 5 Nitrile Gloves, Lg



QTY 4 PES Eye Shield



QTY 1 Surgical Tape, 2 in.

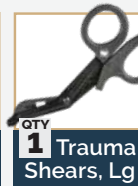

QTY 1 Trauma Shears, Lg



QTY 2 T2® Tag Triage Card



QTY 1 Casualty Reference Card

SPECIFICATIONS

H 7.75 in. x W 13 in. x D 9 in.
Weight: 4 lb 6.5 oz

SQUAD KIT

COLOR	ITEM #	NSN #
Multi-Cam	80-0251	----
Coyote	80-0039	6545-01-543-2956
OD Green	80-0038	6545-01-543-3159
Black	80-0037	6545-01-543-5721
Digital Camo	80-0093	6545-01-543-5733

KIT CAPABILITIES INCLUDE:



MARINER KIT

Meets requirements of EMAT



MARINER KIT

ITEM #	NSN #
80-0107	6545-01-543-5770

SPECIFICATIONS

H 7.75 in. x W 13 in. x D 9 in.
Weight: 4 lb 6.5 oz

Bag only option also available at NARescue.com

- Designed to meet the rigorous requirements of the Emergency Med Assist Team (EMAT) Response Bag
- Standard casualty response bag for the Submarine Forces
- Rugged rescue orange nylon bag for high visibility, easy transport & convenient storage
- Compact, contoured profile with large central compartment, multiple attachment points & front trauma-shear storage
- Designed to be worn as a belly bag or slung over the shoulder

MARINER KIT:



QTY 2 C-A-T® Tourniquet



QTY 4 ETD™ 6 in. Dressing



QTY 1 ETD™ Abdominal



QTY 4 Gauze S-Rolled



QTY 2 28F Naso Airway & Lube



QTY 1 Gauze, Petrolatum



QTY 1 HyFin® Vent. Twin Pk



QTY 1 BOA® IV Band



QTY 2 Saline Lock Kit



QTY 2 SAM® Splint II



QTY 1 Surgical Tape, 2 in.



QTY 4 PES Eye Shield



QTY 1 Trauma Shears, Lg



PAIR 5 Nitrile Gloves, Lg



QTY 2 T2® Tag Triage Card



QTY 1 Casualty Reference Card

NAR-4™ CHEST POUCH

Rapid access to critical medical equipment from your chest

KIT CAPABILITIES INCLUDE:



NAR-4 CHEST POUCH:



NAR-4 CHEST POUCH

COLOR	HEMOSTATIC	ITEM #
Multi-Cam	Combat Gauze	80-0247
Coyote	Combat Gauze	80-0175
OD Green	Combat Gauze	80-0289
Black	Combat Gauze	80-0173
Digital Camo	Combat Gauze	80-0176
5 options	Chito Gauze	Online

SPECIFICATIONS

H 7.75 in. x W 9.75 in. x D 4.5 in
Weight: 2 lb 4 oz

- Auxiliary NAR-4 companion kit streamlined to provide rapid access to critical basic & advanced equipment from the chest
- Designed to be secured to the front of the body armor using a MOLLE-style attachment system
- Ambidextrous accessibility to trauma shears



Bag only option also available at NARescue.com

KIT CAPABILITIES INCLUDE:



NAR-4™ AID KIT

Configurable for individual skill set and mission requirements



Bag only option also available at NARescue.com

- Low-profile, load carriage for Tactical Healthcare Professionals
- Modular system with 5 segregated removable equipment pouches to carry equipment consistent with TCCC Guidelines
- Designed to attach directly to the body armor with four rapid attachment straps or carried using shoulder straps

NAR-4 AID KIT

COLOR	HEMOSTATIC	ITEM #
Black	Combat Gauze	80-0181
Coyote	Combat Gauze	80-0183
Digital Camo	Combat Gauze	80-0184
Multi-Cam	Combat Gauze	80-0244
OD Green	Combat Gauze	80-0308
5 options	ChitoGauze	Online

SPECIFICATIONS

H 19 in. x W 15 in. x D 7 in. | Weight: 12 lb 8.3 oz

NAR-4 AID KIT:





MINI MEDIC KIT

Compact, multi-purpose trauma aid bag with vital medical gear

KIT CAPABILITIES INCLUDE:



- Compact, multi-purpose trauma aid bag with vital medical equipment for point of wounding treatment of multiple injuries
- Two external pockets with two external tourniquet holders
- One external pouch with additional drop-down pouch for use with breaching tools or storage for a hydration pack
- Expansion zipper to increase load capacity
- External MOLLE on sides and bottom of bag / Internal MOLLE rows on inside back of bag
- Three clear tear away pouches on the inside main compartment
- Two internal mesh pouches and one larger internal pouch on inner lid
- Internal tape roll strap

Inside Back Compartment for stowing extra equipment



Every MINI MEDIC KIT includes these core items:

QTY 2 C-A-T® TQ Holder	QTY 4 C-A-T® Tourniquet	QTY 4 Gauze S-Rolled	QTY 4 ETD™ 4 in. Flat	
QTY 1 Tactical Comp Wrap	QTY 2 28F Naso Airway & Lube	QTY 2 HyFin® Vent. Twin Pk	QTY 2 BurnTec 4"x4" Dressing	QTY 2 SAM® Splint II
QTY 1 NAR QuikLitter LITE	QTY 2 Surgical Tape 2" Flatfold	QTY 1 Mini, Tape Multipurpose	QTY 1 Trauma Shears, Lg	QTY 1 Scissor Leash™
PAIR 5 Nitrile Gloves, Lg	QTY 4 Combat Casualty Card	QTY 2 Permanent Marker, Large	QTY 1 Casualty Reference Card	QTY 1 Medical Cross Patch

PLUS+ BASIC also includes:

QTY 4 Gauze S-Rolled

COMBAT GAUZE also includes:

QTY 4 Combat Gauze Z-fold

ADVANCED also includes:

QTY 4 Combat Gauze Z-fold

QTY 1 BAC-Pack Bougie Aided Cric

QTY 2 ARS® Needle, 14g

QTY 1 Cyclone Pocket BVM

SPECIFICATIONS

H 14 in. x W 11 in. x D 5 in. (Depth expandable to 7.25 in.)
Weight: 12.5 lb

Bag only option also available at NARescue.com

MINI MEDIC KIT

COLOR	LEVEL	ITEM #
Black	Basic	80-0979
Coyote	Basic	80-0980
OD Green	Basic	80-0981
Ranger Green	Basic	80-1003
Black	Combat Gauze	80-0985
Coyote	Combat Gauze	80-0986
OD Green	Combat Gauze	80-0987
Ranger Green	Combat Gauze	80-1005
Black	Advanced	80-0982
Coyote	Advanced	80-0983
OD Green	Advanced	80-0984
Ranger Green	Advanced	80-1004

INTRO
M-A-R-C-H
KITS
TAC-EVAC
MED-SURG
RESCUE
TREATMENT SUPPORT
TRAINING
INDEX

KIT CAPABILITIES INCLUDE:



Slim Profile RUF KIT

- Slim profile that can be worn with single bandolier style sling strap and with a detachable waist strap
- Two tear away internal panels designed to hold standard IFAK contents
- Two additional tear away clear accessory pouches for additional equipment
- Internal shock cord area to hold compact, soft litter or additional equipment
- External MOLLE on sides of bag for external tourniquet pouches
- Large panel on exterior of kit for ID patches

RUF KIT

LEVEL	ITEM #
Basic	80-0974
Combat Gauze LE	80-0975
Advanced	80-0976

Bag only option also available at NARescue.com

SPECIFICATIONS

H 18 in. x W 8 in. x D 2 in.
Weight: 6.3 lb

Every RUF KIT includes these core items:



PLUS+ BASIC

also includes:




COMBAT GAUZE LE

also includes:



ADVANCED also includes:

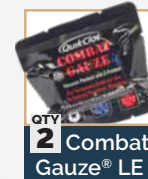

Can be worn with single bandolier style sling strap as well as with a detachable waist strap





Medic TRAUMA PACK

Comprehensive trauma management platform

KIT CAPABILITIES INCLUDE:



- 9 removable, color-coded pull bags and 2 removable side compartments provide modular organization
- Over 1,000 cubic inches of gear capacity and MOLLE-style attachment system compatible with removable side pouches



MEDIC TRAUMA PACK **NATO NSN

COLOR	ITEM #	NSN #
Black	80-0070	6545-01-572-7251
Coyote	80-0072	6545-01-588-1862
Digital Camo	80-0096	----
OD Green	80-0071	6545-01-594-5617**

SPECIFICATIONS

H 19 in. x W 19 in. x D 10 in
Weight: 18 lb 1 oz

MEDIC TRAUMA PACK:

 QTY 2 C-A-T® Tourniquet	 QTY 6 ETD™ 6 in. Dressing	 QTY 2 ETD™ Abdominal	 QTY 4 Gauze S-Rolled	 QTY 2 28F Naso Airway & Lube
 QTY 1 King LT-D Size 4	 QTY 1 Cyclone Pocket BVM	 QTY 1 Tactical Suction Device	 QTY 2 Gauze, Petrolatum	 QTY 2 ARS® Needle, 14g
 QTY 1 BOA® IV Band	 QTY 4 Saline Lock Kit	 QTY 1 Sharps Shuttle	 QTY 1 HRS, Heat Reflective Shell	 QTY 1 Tactical Traction Splint
 QTY 4 SAM® Splint II	 QTY 1 Trauma Shears, Lg	 QTY 2 Surgical Tape, 2 in.	 QTY 1 Armadillo Med Case	 QTY 1 Head Lamp
 PAIR 5 Nitrile Gloves, Lg	 QTY 2 PPE Kit	 QTY 6 PES Eye Shield	 QTY 4 T2® Tag Triage Card	 QTY 1 Casualty Reference Card

Bag only option also available at NARescue.com



9 removable, color-coded pull bags for modular organization

KIT CAPABILITIES INCLUDE:



Warrior Aid & Litter Kit **WALK**[®]



- Official vehicle life-saving kit of the U.S. Army
- For treatment of multiple trauma casualties and the immediate evacuation of one patient
- Includes full-size Talon[®] 90C quad-fold litter
- Pocket system holds self or buddy-aid medical supplies
- Back pack straps provide easy transport from the vehicle to the point of need
- Includes Hypothermia Prevention & Management Kit[™] to help stabilize patients during transport
- Interfaces with any standard or non-standard TACEVAC asset

WARRIOR AID & LITTER KIT

COLOR	ITEM #	NSN #
Black	80-0097	----
Coyote	80-0089	6545-01-587-1199
Digital Camo	80-0059	6545-01-532-4962
OD Green	80-0060	----

SPECIFICATIONS

H 22.75 in. x W 16.5 in. x D 11.5 in.
Weight: 29 lb 15 oz

WARRIOR AID & LITTER KIT:

 QTY 2 C-A-T [®] Tourniquet	 QTY 6 Trauma Dressing, 6"	 QTY 1 ETD [™] Abdominal	 QTY 4 Gauze S-Rolled	 QTY 2 28F Naso Airway & Lube
 QTY 2 Hyfin [®] Occlusive, Sgl.	 QTY 2 ARS [®] Needle, 14g	 QTY 1 HPMK [®] Kit	 QTY 2 SAM [®] Splint II	 QTY 1 Talon 90C Litter
 QTY 4 Litter Ret. Strap	 QTY 1 Casualty Equip. Bag	 QTY 1 Aviation Panel	 QTY 1 Trauma Shears, Lg	 PAIR 5 Nitrile Gloves, Lg
 QTY 1 Tape, 2 in. Gecko-Grip	 QTY 1 Surgical Tape, 2 in.	 QTY 6 PES Eye Shield	 QTY 2 Combat Casualty Card	 QTY 1 Casualty Reference Card

WALK ReSupply Kit (Item# 85-0020):

Contains WALK components for resupply - excluding nylon bag & litter



Bag only option also available at NARescue.com

INTRO

M.A.R.C.H

KITS

TACEVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX



NAVAL FIRST AID BOX

Response Kit Designed to stow away perfectly in the Naval First Aid Box

KIT CAPABILITIES INCLUDE:



- First responder rapid access and organization of critical life saving equipment for rendering on scene emergency care at sea
- Designed to meet the rigorous requirements of the U. S. Naval Ships, Military Sealift Command Ships and afloat commercial vessels
- Lightweight, durable slimline pouch is securable using anti-pilferage seals and features a clear window to readily view contents
- Pouch easily fits into the standard Naval First Aid Box (*not included*) and is ideal for restocking the Naval Gun Response Aid Kit as needed



Easily fits into the standard Naval First Aid Box (*not included*) and is ideal for restocking the Naval Gun Response Aid Kit (see below) as needed

Bag only option also available at NARescue.com

SPECIFICATIONS

H 15 in. x W 6 in. x D 4 in.
Weight: 1 lb 15 oz

NAVAL FIRST AID BOX RESPONSE KIT

OPTION	ITEM #	NSN #
Kit	80-0419	6545-01-648-2705
Pouch only	80-0357	8465-01-666-7328

NAVAL FIRST AID BOX RESPONSE KIT:

 QTY 1 C-A-T® Tourniquet	 QTY 1 ETD™ 6 in. Dressing	 QTY 1 ETD™ Abdominal	 QTY 1 28F Naso Airway	 QTY 1 34F Naso Airway	 QTY 1 Surgical Lubricant
 QTY 1 CPR Mask, Res-Cue Key	 QTY 1 HyFin® Vent. Twin PK	 QTY 2 SAM® Splint II	 QTY 1 Elastic Bandage, 3in.	 QTY 1 Trauma Shears, Lg	 PAIR 1 Nitrile Gloves, Lg
				 QTY 1 Permanent Marker, Large	 QTY 1 Fox Eye Shield w/garter



NAVAL GUN Response Aid Kit

Response Aid Kit Designed for the rigorous requirements of the US Navy

KIT CAPABILITIES INCLUDE:



- Designed for the rigorous requirements of the U.S. Navy for first responder emergency care equipment on board U.S. Naval Ships, Military Sealift Command Ships & afloat commercial vessels
- Ergo dynamic, low profile casualty response system features 2 removable equipment pouches that attach via hook-&-loop to a lightweight, durable outer bag
- Two internal pouches are identical & provide rapid access & organization of critical life saving equipment needed to render on-scene care at sea
- Both outer bag & 2 inner pouches are securable using anti-pilferage seals
- Slim-line, lightweight kit stores easily on a wall and can be carried over the shoulder or bandolier-style



SPECIFICATIONS

H 15 in. x W 6 in. x D 4 in.
Weight: 1 lb 15 oz

See Above for Contents of the Naval First Aid Box Response Kit



NAVAL GUN RESPONSE AID KIT

OPTION	ITEM #	NSN #
Kit	80-0418	6545-01-648-2690
Bags only	80-0417	8465-01-666-7169

Bag only option also available at NARescue.com

KIT CAPABILITIES INCLUDE:



NAR-5 SAR AID KIT

Specifically designed for U.S. Navy Search and Rescue Swimmers



Bag only option also available at NARescue.com

- The SAR Level A Medical Bag designed for U.S. Navy Search and Rescue (SAR) Swimmers
- Low profile & durable response bag featuring 6 internal pouches that attach via hook-and-loop allowing access to the specific pouch needed to provide on-scene medical care to the patient
- Various sized modules can be self labeled for the specific mission set; 4 smaller pouches for hemorrhage control products and 2 larger ones to accommodate airway & intravenous products
- Kit can be carried as a backpack using shoulder straps that stow away until needed

NAR 5 SAR KIT

NATO NSN*

OPTION	ITEM #	NSN #
Kit	80-0416	-----
Bags only	80-0399	8465-01-669-5796*

SPECIFICATIONS

H 24 in. x W 15 in. x D 7 in.
Weight: 14 lb 10 oz

NAR-5 AID KIT:

 QTY 4 C-A-T® Tourniquet	 QTY 4 ETD™ 4in. Dressing	 QTY 1 ETD™ Abdominal	 QTY 4 Gauze, Compressed	 QTY 1 CPR Mask
 QTY 4 PPE Kit	 QTY 2 1000 cc IV Solution	 QTY 1 Cyclone Pocket BVM	 QTY 4 HyFin® Vent. Twin Pk	 QTY 1 Survival Blanket
 QTY 2 SAM® Splint II	 QTY 4 Triangular Bandage	 QTY 1 Mini Maglight	 QTY 2 ACE Cervical Collar	 QTY 4 Elastic Bandage, 3 in.
 QTY 1 Trauma Shears, Lg	 QTY 1 Surgical Tape, 1 in.	 QTY 12 Flexible Fabric Bandages	 QTY 1 Sterile Burn Sheet	

Naval | USCG BOAT RESPONSE AID KITS

Provides Naval & USCG first responders with critical life saving equipment while at sea



For additional information about these products designed for emergency medical care on military & commercial afloat vessels, contact North American Rescue or visit www.NARescue.com

BOAT RESPONSE AID KITS

OPTION	ITEM #	NSN #
Naval Boat Response Aid Kit	80-0420	6545-01-648-3228
Naval Boat Response Aid Box Only	80-0354	8115-01-656-8316
USCG Boat Response Aid Kit	80-0353	6545-01-646-2623

NAVAL BOAT RESPONSE AID KIT:

 QTY 2 C-A-T® Tourniquet	 QTY 1 ETD™ 6 in. Dressing	 QTY 1 ETD™ Abdominal	 QTY 1 CPR Mask	 QTY 1 28F Naso Airway
 QTY 1 Surgical Lubricant	 QTY 1 HyFin® Vent. Twin Pk	 QTY 1 SAM® Splint II	 QTY 1 Dry Sterile Burn Dressing	 QTY 1 Elastic Bandage, 3 in.
 QTY 1 ACE Cervical Collar	 QTY 1 Blizzard Survival Blanket	 QTY 1 Fox Eye Shield w/garter	 QTY 1 Gauze Pads 4"x4" (2pk)	 QTY 1 Trauma Shears, Lg
 QTY 1 Biological PPE Kit				



USMC CLS KIT *U.S. Marine Corps Combat Life Saver Kit*

KIT CAPABILITIES INCLUDE:



The Official Combat Lifesaver Kit of the U.S. Marine Corps

- Constructed of durable 1000 denier Cordura® solution-dyed nylon
- Sand and water resistant zippers with silent pulls & ITW Nexus Hardware

USMC CLS ReSupply (#80-0424)

Contains items for resupply in a clear resealable plastic bag - *excludes nylon bag*



Bag only #ZZ-0252 also available at NARescue.com

SPECIFICATIONS

Packaged: L 12 in. x W 10 in. x D 8 in.
Open Length: 37.5 in.
Weight: 6 lb 14 oz

USMC CLS KIT

ITEM #	NSN #
80-0159	6545-01-571-4470

USMC CLS KIT:

 QTY 3 C-A-T® Tourniquet	 QTY 1 Combat Gauze XL	 QTY 3 Gauze, Compressed	 QTY 3 ETD™ 6 in. Flat Dressing	 QTY 1 ETD™ 8"x10" Dressing			
 QTY 1 22F Naso Airway	 QTY 1 Guedel AirWay, 80 mm	 QTY 1 HyFin® Combo Pack	 QTY 4 ARS® Needle, 14g	 QTY 1 SAM® Splint II		 QTY 5 Combat Casualty Card	
 QTY 1 BurnTec 4"x16" Dressing	 QTY 3 Dry Sterile Burn Dressing	 QTY 2 PES Eye Shield	 QTY 1 Trauma Shears, Lg	 QTY 1 Model 8 Rescue Hook			 QTY 8 Alcohol Prep Pads
 QTY 6 Gauze Pads 4"x4" (2pk)	 QTY 2 Surgical Lubricant, 5g	 PAIR 6 Nitrile Gloves, Lg	 QTY 1 Mini Duct Tape, 2 in.				



US Army CLS RESUPPLY KIT

KIT CAPABILITIES INCLUDE:



US ARMY CLS RESUPPLY KIT

HEMOSTATIC	ITEM #
Combat Gauze	85-0197
ChitoGauze	85-0196

- Convenient kit containing the essential items to resupply the US Army CLS Combat Life Saver Kit delivered in a durable and easy-to-open resealable plastic bag
- Available with 2 Hemostatic options: Combat Gauze or ChitoGauze

USMC CLS RESUPPLY:

 QTY 2 C-A-T® Tourniquet	 QTY 2 ETD™ 6 in. Dressing	 QTY 1 ETD™ Abdominal	 QTY 2 Gauze, Compressed	 QTY 2 Choice of Hemostatic	 QTY 1 28F Naso Airway & Lube	 QTY 2 HyFin® Occlusive, Sgl	 QTY 2 ARS® Needle, 14g	 QTY 1 Survival Blanket
 QTY 1 Blanket Ready Heat	 QTY 1 SAM® Splint II	 QTY 3 Triangular Bandage	 QTY 1 Elastic Wrap 6 in. Bandage	 QTY 1 Trauma Shears, Lg	 QTY 1 Surgical Tape, 2 in.	 QTY 4 PES Eye Shield	 PAIR 4 Nitrile Gloves, Lg	 QTY 6 Alcohol Prep Pads
								 QTY 1 Permanent Marker, Small

KIT CAPABILITIES INCLUDE:



SQUAD RESPONDER KIT

Medical equipment Kit for the US Army Combat Life Saver (CLS)



- The NAR alternative to the US Army Combat Lifesaver (CLS) Kit
- Equipped with optimal medical gear to treat 2 to 4 casualties
- Constructed of Durable 1000 Denier Cordura® Solution-Dyed Nylon
- Large central compartment and multiple attachment points
- Sand and water resistant zippers with silent pulls
- ITW Nexus Hardware



SQUAD RESPONDER KIT

ITEM #	NSN #
85-1764	----

Bag only #80-0932 also available at NARescue.com

US ARMY CLS KIT:

 QTY 2 C-A-T® Tourniquet	 QTY 2 Trauma Dressing, 6in.	 QTY 2 Gauze, Compressed	 QTY 2 Combat Gauze Z-fold	 QTY 1 28F Naso Airway & Lube					
 QTY 2 HyFin® Vent. Twin Pk	 QTY 2 ARS® Needle, 14g	 QTY 1 Blizzard Survival Blanket	 QTY 1 Blanket Ready Heat	 QTY 1 SAM® Splint II	 QTY 3 Dry Sterile Burn Dressing	 QTY 1 Elastic Bandage, 6 in.			
 QTY 1 Surgical Tape, 2 in.	 QTY 4 Fox Eye Shield w/garter	 PAIR 4 Nitrile Gloves, Lg	 QTY 6 Alcohol Prep Pads	 QTY 1 Trauma Shears, Lg	 QTY 1 Scissor Leash™	 QTY 1 Model 8 Rescue Hook	 QTY 1 TORQ Tactical Light	 QTY 2 Combat Casualty Card	 QTY 2 Permanent Marker, Large

SPECIFICATIONS

H 14 in. x W 14 in. x D 5 in
Weight: 4 lb 14 oz

SQUAD RESPONDER RESUPPLY

KIT CAPABILITIES INCLUDE:



Items to resupply the Squad Responder Kit, NAR's alternative to the US Army Combat Lifesaver (CLS) Kit

US ARMY CLS RESUPPLY

ITEM #	NSN #
80-2495	----

US ARMY CLS RESUPPLY:

 QTY 2 C-A-T® Tourniquet	 QTY 2 Trauma Dressing, 6in.	 QTY 2 Gauze, Compressed	 QTY 2 Combat Gauze Z-fold	 QTY 1 28F Naso Airway & Lube	 QTY 2 HyFin® Vent. Twin Pk	 QTY 2 ARS® Needle, 14g	 QTY 1 Blizzard Survival Blanket	 QTY 1 Blanket Ready Heat	 QTY 1 SAM® Splint II
 QTY 1 Elastic Bandage, 6 in.	 QTY 1 Surgical Tape, 2 in.	 QTY 4 Fox Eye Shield w/garter	 PAIR 4 Nitrile Gloves, Lg	 QTY 6 Alcohol Prep Pads	 QTY 1 Trauma Shears, Lg	 QTY 1 Scissor Leash™	 QTY 1 Model 8 Rescue Hook	 QTY 2 Combat Casualty Card	 QTY 2 Permanent Marker, Large

KIT CAPABILITIES INCLUDE:



Expeditionary Casualty Response Kit **ECRK**

- Constructed of 500D Cordura® solution dyed nylon
- Sand- and water-resistant, self repairing zippers with silent 550 pull cords
- Separate lid compartment for the R.E.E.L. Non-Rigid Litter
- Six segregated removable pouches secured by hook and loop:
 - > 4 x Internal Bleed Pouches with Red pull handle
 - > 1 x Airway Pouch with Blue pull handle
 - > 1 x Accessories Pouch with Coyote pull handle
- 2 internal lid pouches to secure additional medical products
- 2 external side pouches on either side of the kit for additional medical products
- 2 external straps with dual adjust MILSPEC IR reduced buckles to secure the kit in a vehicle or boat

SPECIFICATIONS

Case closed:
H 17 in. x W 16.5 in. x D 10 in.
Weight: 18 lb

ECRK

ITEM #	NSN #
80-0889	----

Bag only **#80-0888**
also available at NARescue.com



Lightweight platform to rapidly access and organize all your essential trauma management equipment

EXPEDITIONARY CASUALTY RESPONSE KIT:

 QTY 6 C-A-T® Tourniquet	 QTY 4 Combat Gauze Z-fold	 QTY 10 Gauze, Compressed	 QTY 6 ETD™ 6 in. Dressing	 QTY 1 ETD™ Abdominal
 QTY 4 28F Naso Airway & Lube	 QTY 4 HyFin® Chest Seal Combo Pk	 QTY 4 ARS® Needle, 14g	 QTY 1 Blizzard Survival Blanket	 QTY 1 R.E.E.L. Litter Black
 QTY 1 BurnTec Minor Burn Kit	 QTY 3 Elastic Bandage, 6 in.	 QTY 2 SAM® Splint II	 QTY 2 Triangular Bandage	
 QTY 6 PES Eye Shield	 QTY 2 Surgical Tape, 1 in.	 QTY 10 Antiseptic Wipes	 PAIR 6 Nitrile Gloves, Lg	 QTY 1 Trauma Shears, Lg
 QTY 1 Model 7 Rescue Hook	 QTY 2 Mini Duct Tape, 2 in.	 QTY 6 Combat Casualty Card	 QTY 2 Permanent Marker, Small	 QTY 1 Casualty Reference Card

INTRO
 M-A-R-C-H
 KITS
 TAC-EVAC
 MED-SURG
 RESCUE
 TREATMENT SUPPORT
 TRAINING
 INDEX

RANGE TRAUMA KIT

Addresses injuries common to shooting and training ranges

KIT CAPABILITIES INCLUDE:



- Essential equipment for medical and non-medical personnel to address injuries common to shooting and training ranges
- Treatment capabilities include penetrating injuries, chest wounds, burns, eye injuries, minor lacerations, sprains and fractures

RANGE TRAUMA KIT

*NATO NSN#

HEMOSTATIC	COLOR	ITEM #	NSN #
Combat Gauze	Orange	80-0299	----
Chito Gauze	Orange	80-0213	6545-01-610-5411*

SPECIFICATIONS

H 7.75 in. x W 13 in. x D 9 in
Weight: 5 lb 10 oz



Bag only option also available at NARescue.com

RANGE TRAUMA KIT:

 QTY 2 C-A-T® Tourniquet	 QTY 4 ETD™ 6 in. Dressing	 QTY 4 Gauze S-Rolled	 QTY 1 Choice of Hemostatic	 QTY 1 CPR Mask	 QTY 1 HyFin® Vent. Twin Pk	 QTY 2 SAM® Splint II	 QTY 2 Triangular Bandage	 QTY 2 Elastic Bandage, 4 in.	 PAIR 5 Nitrile Gloves, Lg	 QTY 4 PES Eye Shield	 QTY 1 Surgical Tape, 2 in.	 QTY 1 Trauma Shears, Lg	 QTY 1 Forceps, Splinter, 4.5 in.	 QTY 5 Gauze Pads 4"x4" (2pk)	 QTY 1 Eye Wash, 4 oz	 QTY 50 Flexible Fabric Bandages	 QTY 2 BurnTec 5"x10" Dressing
-----------------------------	-------------------------------	--------------------------	--------------------------------	--------------------	--------------------------------	--------------------------	------------------------------	----------------------------------	-------------------------------	--------------------------	--------------------------------	-----------------------------	--------------------------------------	----------------------------------	--------------------------	-------------------------------------	-----------------------------------

RANGE TRAUMA AID KIT

in Hard Case



- Water-tight case maintains the integrity of your medical products when you need them most
- Kit contents arranged in a modular system featuring rugged lid storage organizer and removable nylon pouches, in the main compartment, with easy-to-read treatment labels

RANGE TRAUMA AID KIT

COLOR	ITEM #
Orange	85-0889

SPECIFICATIONS

H 9.5 in. x W 13.5 in. x D 6 in.
Weight: 8 lb 12 oz

RANGE TRAUMA AID KIT HARD CASE:

 QTY 2 C-A-T® Tourniquet	 QTY 2 ETD™ 6 in. Dressing	 QTY 2 Gauze, Compressed	 QTY 2 Combat Gauze® LE	 QTY 1 CPR Mask	 QTY 1 HyFin® Vent. Twin Pk	 QTY 2 SAM® Splint II	 QTY 2 Triangular Bandage	 QTY 1 Permanent Marker, Large	 QTY 2 Elastic Bandage, 4 in.	 PAIR 4 Responder Nitrile Gloves, Lg	 QTY 2 PES Eye Shield	 QTY 1 Surgical Tape, 2 in.	 QTY 1 Responder Trauma Shears, Lg	 QTY 1 Forceps, Splinter, 4.5 in.	 QTY 5 Gauze Pads 4"x4" (2pk)	 QTY 1 Eye Wash, 4 oz	 QTY 50 Flexible Fabric Bandages	 QTY 1 BurnTec 5"x5" Dressing
-----------------------------	-------------------------------	-----------------------------	----------------------------	--------------------	--------------------------------	--------------------------	------------------------------	-----------------------------------	----------------------------------	---	--------------------------	--------------------------------	---------------------------------------	--------------------------------------	----------------------------------	--------------------------	-------------------------------------	----------------------------------

KIT CAPABILITIES INCLUDE:



K-9 HANDLER ROO M-FAK



Easy access to C-A-T tourniquet

- One of NAR's most compact, versatile Mini First Aid Kits featuring rapid one handed access to your tourniquet from an easy-to-open front pouch
- Contains the First Responder's most requested critical point-of-wounding medical equipment for treating injuries typical to K-9s and their handlers operating in tactical or patrol environments
- Super compact, rugged nylon platform that allows attachment both vertically (MOLLE backing) and horizontally (3 in. belt loop)
- Multiple elastic loops for secure gear organization
- Includes 1 in. Mini PVC Med Cross ID Patch



K-9 HANDLER ROO M-FAK

COLOR	LEVEL	ITEM #
Black	Basic	80-1067
Coyote	Basic	80-1068
OD Green	Basic	80-1069
Ranger Green	Basic	80-1070
Black	Combat Gauze LE	80-1071
Coyote	Combat Gauze LE	80-1072
OD Green	Combat Gauze LE	80-1073
Ranger Green	Combat Gauze LE	80-1074

Bag only option also available at NARescue.com

SPECIFICATIONS

L 7 in. x W 3.75 in. x D 4 in.
Weight: 13 oz

K-9 HANDLER ROO M-FAK:



Rapid one handed access to your tourniquet without ever having to open the full kit!



KIT CAPABILITIES INCLUDE:



NAR K-9 LITTER & Carrying Strap Accessories



SPECIFICATIONS

Litter: W 35 in. x H 48 in.
Adjustable Carry Strap:
L 9 to 54 in. x W 2 in.
Carabiner open gate clearance: 27 mm
Packaged:
W 5.5 in. x H 10 in. x D 2.5 in.
Weight: 1.38 lb

- K-9 Litter can be carried in cargo pocket or pouch
- Leg cut outs to safety carry the canine in the litter
- Six large and four small carry handles that allow for one to four person carries
- Optional K-9 Carry Strap & 2 Carabiners for over the shoulder carry of injured K-9 also available

K-9 LITTER

OPTIONS	ITEM#
K-9 Litter Only	60-0074W
Carry Strap & Carabiners	ZZ-1004
K-9 Litter with Carry Strap & Carabiners	60-0076



K-9 HANDLER IFAK

Designed especially for K-9 and Handler

KIT CAPABILITIES INCLUDE:



K-9 HANDLER IFAK:



- Designed for immediate lifesaving treatment to a wounded handler, operator or K-9
- Treatment capabilities range from major trauma typical to operating in tactical or patrol environments to K-9 emergencies

SPECIFICATIONS

H 7 in. x W 4 in. x D 2.75 in | Weight: 1 lb 3 oz

K-9 HANDLER IFAK

COLOR	HEMOSTATIC	ITEM #
Black	----	80-0439
Coyote	----	80-0440
Multi-Cam	----	80-0441
OD Green	----	80-0442
Black	Combat Gauze	80-0443
Coyote	Combat Gauze	80-0444
Multi-Cam	Combat Gauze	80-0445
OD Green	Combat Gauze	80-0446

Bag only option also available at NARescue.com

KIT CAPABILITIES INCLUDE:




K-9 TACTICAL FIELD KIT

Kit is designed specifically to treat wounded K-9s



K-9 TACTICAL FIELD KIT:

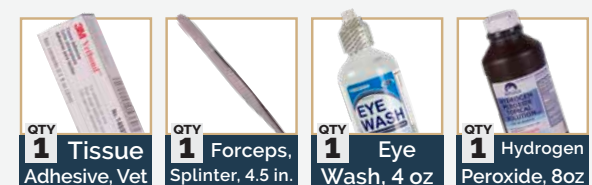

- Designed for the handler or medic to deliver lifesaving treatment to a wounded K-9 at the point of injury
- K-9 treatment capabilities from major trauma typical to operating in tactical or patrol environments to more common K-9 emergencies
- Provides the essential equipment to render care to the injured K9 when veterinary care is not immediately available

SPECIFICATIONS

H 7.75 in. x W 13 in. x D 9 in.
Weight: 5 lb 2 oz

K-9 TACTICAL FIELD KIT

COLOR	HEMOSTATIC	ITEM #
Black	Combat Gauze	80-0300
Coyote	Combat Gauze	80-0301
OD Green	Combat Gauze	80-0303
3 options online	ChitoGauze	See NARescue.com

Bag only option also available at NARescue.com

INTRO

MAR·CH

KITS

TAC·EVAC

MED·SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX


K-9 MINI MEDIC KIT

Versatile, Multi-purpose Platform for K9 Trauma in the field

KIT CAPABILITIES INCLUDE:



- Compact, multi-purpose trauma aid bag with life-saving trauma equipment and K9 specific gear for point of wounding treatment of multiple traumatic injuries
- K9 specific gear enables handlers to render care to an injured canine when veterinary care is not immediately available
- Versatile platform features two external pockets with two external tourniquet holders, one external pouch with additional drop-down pouch for extra storage
- External MOLLE on sides and bottom of bag, internal MOLLE rows on inside back of bag, three clear tear away pouches | two internal mesh pouches on the inside main compartment, and one larger internal pouch on inner lid
- Internal tape roll strap also included

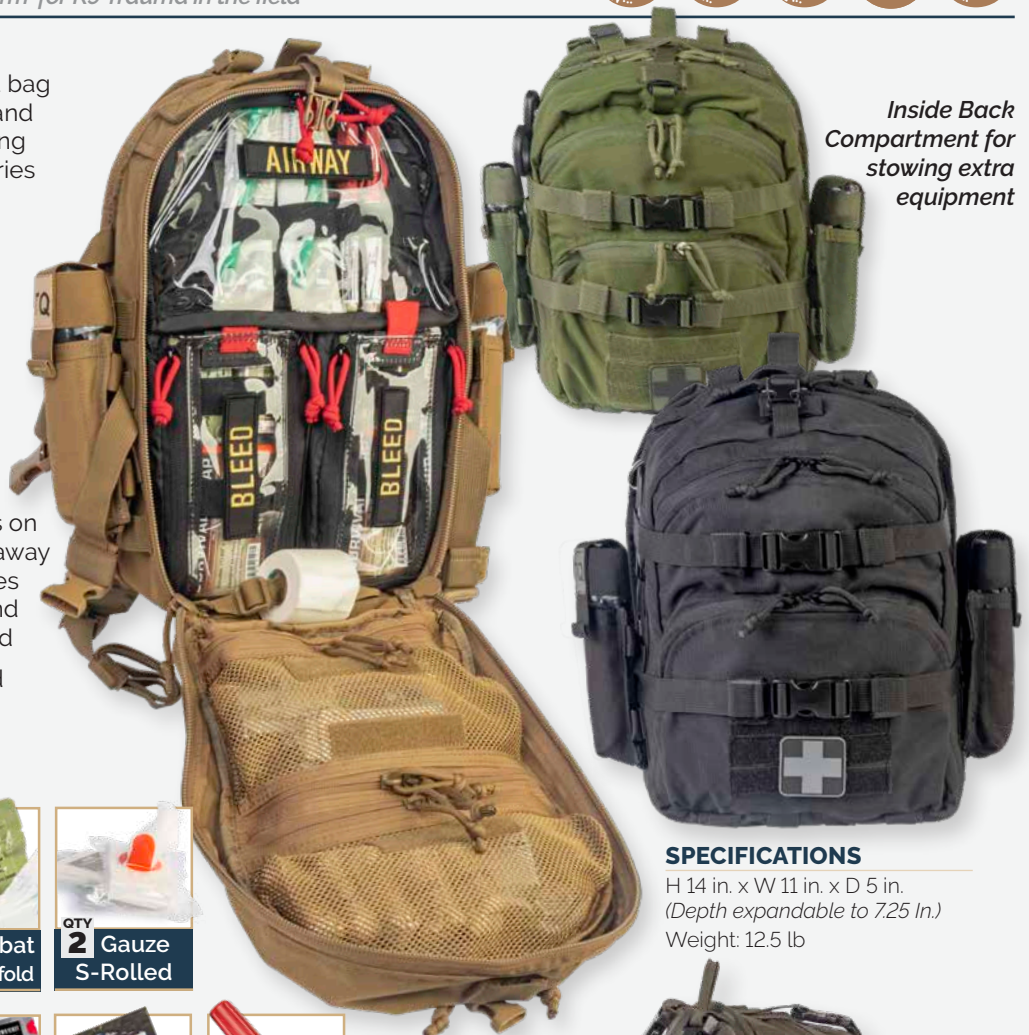

SPECIFICATIONS

H 14 in. x W 11 in. x D 5 in.
(Depth expandable to 7.25 in.)
Weight: 12.5 lb

K-9 MINI MEDIC KIT:

 QTY 4 C-A-T® Tourniquet	 QTY 2 C-A-T® TQ Holder	 QTY 2 Combat Gauze Z-fold	 QTY 2 Gauze S-Rolled
 QTY 4 ETD™ 4 in. Flat	 QTY 2 28F Naso Airway & Lube	 QTY 2 HyFin® Vent. Twin Pk	 QTY 4 Gauze, Petrolatum
 QTY 2 ARS® Needle, 14g	 QTY 2 Survival Blanket, Org	 QTY 1 Tactical Muzzle	 QTY 1 NAR K-9 Litter
 QTY 4 Green ChemLight™	 QTY 2 (6 Pk) IR Square	 QTY 1 Head Lamp	 QTY 10 Flexible Fabric Bandages
 QTY 1 Surgical Tape, 2 in.	 QTY 1 Trauma Shears, Lg	 QTY 1 Veterinary Thermometer	 PAIR 5 Nitrile Gloves, Lg
 QTY 2 Permanent Marker, Large	 QTY 1 Medical Cross Patch		



K-9 MINI MEDIC KIT

COLOR	ITEM #
Black	80-1075
Coyote	80-1076
OD Green	80-1077
Ranger Green	80-1078

Bag only option also available at NARescue.com



K-9 RUF KIT



Slimline Aid Pack ideal for K9 Handler or Tracking Team




- Slim profile trauma aid pack designed for the K9 Handler or Tracking Team member
- Carries a minimum of two complete sets of IFAK contents along with K9 medical supplies
- Carriage options include single bandolier style sling strap as well as detachable waist strap
- Two tear away internal panels hold standard IFAK contents
- Two additional tear away clear accessory pouches to store additional equipment
- Internal shock cord area to hold K9 Litter
- External MOLLE on sides of bag for external tourniquet pouches
- Large exterior loop panel for ID patches

Bag only option also available at NARescue.com

SPECIFICATIONS


H 18 in. x W 8 in. x D 2 in.
Weight: 6.3 lb

K-9 RUF KIT:

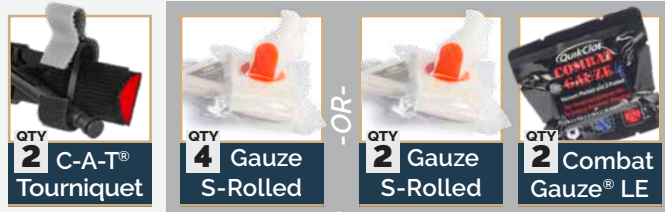

K-9 RUF KIT

OPTION	ITEM #
S-Rolled Gauze	80-1027
with Combat Gauze LE	80-1026

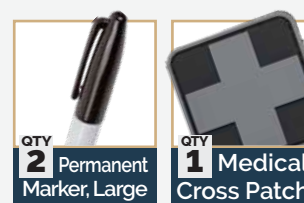

Can be worn with single bandolier style sling strap as well as with a detachable waist strap



MULTI-MISSION TRAUMA PACK

Customizable modular solution ideal for personalized gear organization spanning multiple Mission profiles

- A 3200-cubic inch backpack with (4) compartments
- Multi-mission belt attached for different configurations
- Internal Molle straps covered in hook and loop
- Detachable compartments
- Converts to trauma sheet for vehicles, boats, or aircrafts
- Internal channeling for electrical communications wiring
- Dedicated loops, providing rappelling or lowering of the bag

SPECIFICATIONS

H 25.5 in. x W 11 in. x D 9 in.
Weight: 7 lb 8 oz

MULTI-MISSION TRAUMA PACK *NATO NSN

COLOR	ITEM #	NSN #
Black	80-0950	6545-01-695-5545*
Coyote	80-0949	-----



Tiered Trauma Medical Bag **TTMB**



A. ReSupply Bag
Item# 80-0271
(MTC only)



B. NAR-4 Aid Bag
Item# 80-0241 (MTC)

Other colors available
see NARescue.com



C. NAR-4 Chest Pouch
Item# 80-0248 (MTC)
Other colors available



D. Leg Rig Conversion Panel
Item# 80-0272 (MTC only)

E. MOLLE Directional Conversion Panel
Item# 80-0273
(MTC only)

- 5-piece modular ensemble targeting the military medic's "Soldier-as-a-System" requisite
- Includes low-profile backpack-style aid bag, ample re-supply bag, medical chest pouch, leg rig conversion panel & directional MOLLE converter
- Interchangeable internal pouches with removable labels compatible with the NAR-4 system or any hook and loop attachment
- Leg rig converter & directional change MOLLE panel allow the chest pouch to be mounted & carried using a drop-down, leg-attachment style platform

SPECIFICATIONS

- (A) ReSupply Bag: H 20 in. x W 15 in. x D 5 in.
 (B) NAR-4 Aid Bag: H 20 in. x W 15 in. x D 5 in.
 (C) NAR-4 Chest Pouch: H 7.75 in. x W 9.75 in. x D 4.5 in.
 (D) Leg Rig Conversion Panel: L 14 in. x W 6 in.
 (E) MOLLE Conversion Panel: L 7.5 in. x W 6 in.

TIERED TRAUMA MEDICAL BAG (TTMB)

BAG	COLOR	ITEM #
TTMB Complete System	Multi-Cam	80-0274
ReSupply Bag (A)	Multi-Cam	80-0271
NAR-4 Aid Bag (B)	Multi-Cam	80-0241
NAR-4 Chest Pouch (C)	Multi-Cam	80-0248
Leg Rig Panel (D)	Multi-Cam	80-0272
MOLLE Conv. Panel (E)	Multi-Cam	80-0273

NOTE: Bags only all contents sold separately

INTRO

M.A.R.C.H

KITS

TAC-EVAC

MED-SURV

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX



OCHO BAGS

Low Visibility | Customizable Modular Solutions

- Low visibility aid bags purposely designed to resemble standard off-the-shelf bags to minimize the threat of theft and maximize low visibility discretion
- Featuring hidden compartments with interchangeable pouches backed with hook and loop fields
- Compatible with CrossBreed Holsters (not included)



OCHO MESH AID POUCH

SPECIFICATIONS

H 7 in. x W 4.5 in. x D 3 in.
Weight: 187 oz

COLOR	ITEM #
Black/ODG	80-0259

OCHO ACCESSORY BAG

SPECIFICATIONS

H 7 in. x W 2.25 in. x D 3.25 in.
Weight: 2.75 oz

COLOR	ITEM #
Blue/Gray	80-0371
Gray/Black	80-0370
Red/Gray	80-0571



OCHO Accessory Bag does NOT include Mesh Aid Pouches

OCHO SHOCK CORD PANEL

SPECIFICATIONS

H 6 in. x W 3 in. x D 0.25 in.
Weight: 116 oz

COLOR	ITEM #
Black	80-0260



Gear sold separately

MEDIC TRAUMA SHEET



- Versatile, rugged, lightweight equipment access platform for managing casualties during TACEVAC operations
- Oversized pockets & large grommet attachment points facilitate rapid set-up and equipment retrieval
- Quick attach and release straps allow for secure attachment to the inside wall of an aircraft or vehicle

SPECIFICATIONS

Open: H 41.5 in. x W 21.5 in. x D 2.5 in.
Folded: H 20.25 in. x W 21.5 in. x D 5 in.
Weight: 4 lb 12 oz

Bag only; contents are sold separately

MEDIC TRAUMA SHEET

COLOR	ITEM #	NSN #
Black	80-0067	6515-01-540-7612
Coyote	80-0069	6515-01-540-7603
Digital Camo	80-0094	6515-01-540-7579
OD Green	80-0068	6515-01-540-7585





Horizontal & Vertical Mounting capability available with the 4.5 x 6.5 in. Panel



- Allows the operator the fastest access to their kit during the critical moments after wounding
- Allows any Molle/PALS kit to be detachable
- Gives user faster access to Molle/PALS equipped kits
- Mil-Spec hook & loop and heavy duty snaps to ensure kit security
- Designed to be compatible with all of NAR IFAK's, including a **NEW horizontal mount option** in the medium version (W 4.5 x H 6.5 in.)

Small QUICK DETACH PANEL (3 in. x 6 in.)

- Utilizes only 3 rows by 2 columns on user vest/bag, vertical mount configuration
- Designed specifically for NAR's most compact M-FAK kits

SPECIFICATIONS

W 3 in. x H 6 in. x D 0.5 in.
Weight: 2.8 oz

SMALL PANEL

SIZE	COLOR	ITEM #
3 x 6 in.	Black	80-0893
3 x 6 in.	Coyote	80-0894



Medium DUAL DIRECTION QUICK DETACH PANEL (4.5 in. x 6.5 in.)

- Utilizes 3 rows by 3 columns vertical mount or 2 rows by 4 columns horizontal mount
- Designed for compatibility with all NAR IFAKs such as horizontal mount for M-FAK and vertical or horizontal mount of the TORK

SPECIFICATIONS

W 4.5 in. x H 6.5 in. x D 0.5 in.
Weight: 4.5 oz

DUAL DIRECTION PANEL

SIZE	COLOR	ITEM #
4.5 x 6.5 in.	Black	80-0895
4.5 x 6.5 in.	Coyote	80-0896



Large QUICK DETACH PANEL (9 in. x 6.5 in.)

- Utilizes 3 rows by 7 columns on user vest/bag, vertical mount configuration
- Designed specifically for use with NAR 4 Chest Pouch & Emergency Evacuation Litter (EEL)

SPECIFICATIONS

W 9 in. x H 6.5 in. x D 0.5 in.
Weight: 5.6 oz

LARGE PANEL

SIZE	COLOR	ITEM #
9 x 6.5 in.	Black	80-0897
9 x 6.5 in.	Coyote	80-0898





D-BCRK *Department of Defense* Individual BLEEDING CONTROL RESPONSE KITS

- Individual Bleeding Control Response Kit designed to provide essential equipment that empowers DOD Forces and Federal Employee responders in public areas to take immediate action in stopping life threatening bleeding.
- Contained in a ruggedized, nylon carry bag or vacuum-sealed pouch
- Provides intuitive and easy-to-use tools that are proven to help to save lives
- Available configurations: Basic, Intermediate, & Advanced

Nylon Bag Option:

- Easy to restock in cases of use
- Convenient to carry and easy to open in stressful situations
- Opens flat for ease of use and immediate identification of enclosed devices
- Easily recognizable red nylon bag with ample hook & loop field for attaching an ID patch

Vacuum Sealed Option:

- Cost-effective alternative to Nylon
- Compact and rugged Vacuum-Sealed pouch
- Large Identifying Card Vacuum-Sealed inside the pouch for ease of identification

SPECIFICATIONS

Nylon Bag: H 7 in. x W 5 in. x D 2.75 in.

Vacuum Sealed: H 7 in. x W 5 in. x D 2.75 in.

D-BCRK INDIVIDUAL KITS

LEVEL	BAG	WEIGHT	ITEM #
Basic	Red Nylon	15 oz	85-2057
Intermediate	Red Nylon	16.5 oz	85-2058
Advanced	Red Nylon	17.4 oz	85-2059
Basic	Vacuum Sealed	11.7 oz	85-1749
Intermediate	Vacuum Sealed	13.25 oz	85-2066
Advanced	Vacuum Sealed	14 oz	85-2069

Also available at NARescue.com:

#ZZ-1000 Individual Clear Polycarbonate Wall Cabinet



Every D-BCRK INDIVIDUAL BLEEDING CONTROL RESPONSE KIT:



PLUS these Individual Components based upon configuration of choice:



D-BCRK 8-Pack Department of Defense Bleeding Control Response Kits (Nylon or Vacuum-Sealed)



- Eight (8) Nylon or Vacuum-sealed Individual Bleeding Control Response Kits inside an easy-to-open nylon carrying case
- Available in Basic, Intermediate, and Advanced configurations

SPECIFICATIONS

Each 8-Pack: H 14 in. x W 11 in. x D 8 in.
 Individual Bleeding Control Kit:
 H 7 in. x W 5 in. x D 2.5 in.



D-BCRK 8-PACKS

LEVEL	BAG	WEIGHT	ITEM #
Basic	Red Nylon	8.5 lb	85-2072
Intermediate	Red Nylon	9.3 lb	85-2073
Advanced	Red Nylon	9.7 lb	85-2074
Basic	Vacuum Sealed	8.5 lb	85-1750
Intermediate	Vacuum Sealed	9.3 lb	85-2067
Advanced	Vacuum Sealed	9.7 lb	85-2070

D-BCRK 8-PACK - Vac



D-BCRK 8-PACK - Nylon



See page left for contents of the Individual Kits

D-BCRK Wall Station Department of Defense Bleeding Control Response Kits (Nylon or Vacuum-Sealed)

- Wall Station contains One D-BCRK 8-Pack each containing Eight Nylon or Vacuum-sealed Individual Bleeding Control Response Kits
- Available in Basic, Intermediate, and Advanced configurations
- Includes a durable, clear polycarbonate cabinet with rounded corners provides both physical and environmental protection for medical equipment
- Includes a wire shelf to hold the 8-pack & an external eyelet that allows for a secondary safety seal to easily identify if the station has been opened or tampered with
- External mounting tabs and hardware enable easy installation on any wall
- Clear design allows visibility at any angle

SPECIFICATIONS

Clear Wall Case: H 18 in. x W 15.5 in. x D 10 in.
 Each 8-Pack: H 14 in. x W 11 in. x D 8 in.
 Individual Bleeding Control Kit: H 7 in. x W 5 in. x D 2.5 in.

D-BCRK WALL STATIONS

LEVEL	BAG	WEIGHT	ITEM #
Basic	Red Nylon	23.5 lb	85-2075
Intermediate	Red Nylon	24.3 lb	85-2076
Advanced	Red Nylon	24.7 lb	85-2077
Basic	Vacuum Sealed	23.5 lb	85-1751
Intermediate	Vacuum Sealed	24.3 lb	85-2068
Advanced	Vacuum Sealed	24.7 lb	85-2071

D-BCRK CLEAR WALL STATION:



Also available at NARescue.com:

#ZZ-0921 Semi-Recessed Metal Wall Cabinet

Cases are water resistant. To best protect your critical lifesaving products, cases are best mounted under cover and out of the direct elements.



TACEVAC

CASUALTY EVACUATION

- FRAMED LITTERS • ACCESSORIES
- NON-RIGID LITTERS
- SPINE BOARDS

see pages 54-59





Rigid Litters | **TACEVAC** 

TALON® II 90C
Collapsible Handle Litter

Advanced, innovative combat evacuation system features portable quad-fold design, retractable handles, six IV-attachment points and two patient securing straps



- Standard 90 in. NATO-compatible evacuation platform
- Developed to U.S. Army requirements
- Certified Safe-to-Fly U.S. Army Rotary Wing
- Retractable, ergonomic handles and automatic locking hinges
- Slip-resistant, chemical-resistant rip-stop fabric tested & approved by Federal & private labs using live warfare agents

SPECIFICATIONS CE
 Folded w/handles collapsed: H 20.5 in. x W 8.5 in. x D 9.5 in.
 Open w/handles collapsed: L 78 in. x W 22.5 in. x D 6 in.
 Open w/handles extended: L 90 in. x W 22.5 in. x D 6 in.
 Weight: 15 lb 6 oz (6.97 kg)
 Static Payload: 1,800 lb (816 kg)
 Ground Clearance: 1.5 in.
 U.S. Patents: 6,842,923 & 5,598,592, Patent(s) Pending

TALON II 90C LITTER

ITEM#	NSN#
60-0002	6530-01-504-9051

RAVEN® 90C
Bi-fold Litter

- The official litter of the U.S. Navy and certified safe-to-fly by the U.S. Air Force
- Bifold design for reliable long-term storage
- Ergonomically designed handles
- Slip-resistant, chemical-resistant, rip-stop fabric tested & approved by Federal & private labs using live warfare agents
- Includes 6 IV-attachment points and 2 patient-securing straps



SPECIFICATIONS CE
 Folded: H 45.5 in. x W 7.5 in. x D 8 in.
 Open: 90 in. x W 22.5 in. x D 6 in.
 Weight: 16 lb 8 oz (7.5 kg)
 Static Payload: 2,000 lb (909 kg)
 Ground Clearance: 1.5 in.
 U.S. Patents: 6,842,923 5,598,592, Patent(s) Pending

RAVEN 90C LITTER

ITEM#	NSN#
60-0001	6530-01-432-5114

T3™ LITTER
Lightweight, Assault

- Ultra light, ultra strong litter weighing less than 12 lb
- Specifically designed for Special Operations Forces
- Provides maximum maneuverability during TACEVAC in urban terrains
- Overall length of 81 in. enables floor loading on a UH-60 helicopter
- 2 straps for securing patients
- Ergonomically designed handles
- Quad-fold design for maximum cube space storage



SPECIFICATIONS CE
 Folded: H 21 in. x W 8.75 in. x D 7.25 in.
 Open: L 81 in. x W 22.5 in. x D 6 in.
 Weight: 11 lb 13 oz (5.34 kg)
 Static Payload: 1,200 lb (544.3 kg)
 Ground Clearance: 1.5 in.

T3 LITTER

ITEM#
60-0022



Wheeled LITTER CARRIER

- Lightest wheeled litter cart on the market for single person transport to and from MedEVAC vehicles
- Compatible with the Talon® II 90C, Raven® 90C, MedEvac & T3, in addition to all NATO litters
- Quick and easy to deploy or disassemble for storage
- Locking clamp system keeps the litter firmly in place and allows for different litter widths
- Performs as litter stand in the field due to its unique wheel locking system
- Also available with optional Carrying Case

SPECIFICATIONS

Deployed:
L 56.29 in. x W 27.5 in. x H 36.2 in.
Folded:
L 42.1 in. x W 27.5 in. x H 9.84 in.
Case: L 48 in. x W 29.5 in. x H 12.5 in.
Weight:
Carrier only: 38 lb 8 oz (17.5 kg)
Carrier with Case: 88 lb (39.9 kg)
Payload: 550 lb (250 kg)



WHEELED LITTER CARRIER *NATO NSN

OPTION	ITEM #	NSN #
Carrier with Case	60-0063	-----
Carrier only	60-0041	6530-01-663-2308*

LITTER STANDS

33 in. Collapsible

- Set of 2 lightweight, compact litter stanchions designed to collapse for convenient storage in a standard issue military carrying case
- Quickly locks open
- Allows the capacity to store a quad-folded Talon®II 90C (sold separately)



SPECIFICATIONS

H 17.95 in. x W 25.5 in. x D 6.92 in.
Extended: 33 in.
Base Width: 25.5 in.
Weight: 10 lb 12 oz

COLLAPSIBLE 33 in. LITTER STANDS

Unit	ITEM#
Set of 2	60-0064

BACKREST TALON® or RAVEN®

- Allows patients on a Talon® 90C or Raven® 90C litter to be placed at 90° or 45° angles
- Constructed of a lightweight aluminum frame with Nomex/Kevlar-blend fabric covering

SPECIFICATIONS

H 31.5 in. x W 29.25 in. x D 1.6 in.
Weight: 2 lb 13 oz

LITTER BACKREST

ITEM#	NSN#
60-0027	6530-01-594-7444



LITTER STANDS

33 in. Standard

- Set of 2 lightweight, foldable stanchions for standard sized litters
- Constructed of a lightweight aluminum frame that conveniently folds flat for ease in transport & storage
- Quickly locks open to safely secure the litter in place



SPECIFICATIONS

H 34 in. x W 25.5 in. x D 6.92 in.
Extended: 33 in.
Base Width: 23.5 in.
Weight: 8 lb 8 oz

33 in. LITTER STANDS

UNIT	ITEM#	NSN#
Set of 2	60-0014	6530-01-523-3314

IV POLE

Adjustable

Exclusively designed to fit the pre-molded holes of our framed litters

SPECIFICATIONS

H 19-33 in. x W 8.75 in. x D 0.6 in.
Weight: 4.8 oz

IV POLE

ITEM#	NSN#
60-0008	6530-01-504-9054





RIG SERIES CARRIER

Assault Litter Carrier

- Enhanced carrier with additional strap over the top of the litter to hold it in place during drastic movement and FRIES (*Fast Rope Insertion Extraction System*) operations
- Upgraded shoulder straps provide comfort for long wear
- Includes pocket in the top flap for storage and three additional handles facilitate packing and unloading of the litter
- Accommodates all NAR quad-fold litter platforms

- Equipment attaching webbing (MOLLE) along the outside of the carrier allows attachment of additional equipment

SPECIFICATIONS

H 22 in. x W 9 in. x D 9 in.
Weight: 3 lb (1.36 kg)

RIG SERIES LITTER CARRIER

COLOR	ITEM #	NSN #
Multi-Cam	60-0057	6530-01-647-8571*
Coyote	60-0056	6530-01-643-9138*

*NATO NSN

**SPECIFICATIONS**

Collapsed: H 20 in. x W 9 in. x D 2 in.
Expanded: H 20 in. x W 9 in. x D 10 in.
Weight: 1 lb 7 oz

RETENTION STRAP

Litter Support Bar

Designed to prevent the litter support bar from collapsing

RETENSION STRAP**ITEM#**

60-0026

SPECS

L 104 in. x W 1 in.
Weight: 2 oz



LITTER BUCKLE

*Replacement Set***SPECIFICATIONS**

L 3.75 in. x W 2.38 in. x D 0.5 in.
Weight: 1 oz

ORDERING INFO

ITEM#	NSN#
ZZ-0032	6530-01-435-7337



TALON® CARRIER

Assault Litter Carrier

- Transport any NAR quad-fold litter (Talon, MedEvac & T3)
- MOLLE-style attachment points
- Includes thin back plates, drain hole, non-slip, adjustable shoulder straps & top-carrying handle

TALON LITTER CARRIER

*NATO NSN

COLOR	ITEM #	NSN #
Coyote	60-0011	6530-01-531-6854
OD Green	60-0010	6530-01-504-9055
Black	60-0009	6530-01-531-6907
Digital Camo	60-0025	6530-01-643-9134*

EQUIPMENT BAG

Store a casualty's personal items

- Roomy nylon bag with reinforced shoulder straps and info tag to bag & tag casualty's personal items during evacuation
- Lightweight & vacuum-packed for low cube, easy transport

SPECIFICATIONS

H 5 in. x W 5 in. x D 2.5 in.
Weight: 6.25 oz

EQUIPMENT BAG

ITEM#	NSN#
80-0132	8465-01-575-6644





MEDEVAC LITTER



- The ultimate, cost efficient casualty evacuation stretcher for First Responders operating in Active Shooter, Mass Casualty and Disaster incidents
- Compact, collapsible, quad-fold litter is easy to deploy and carry
- Two patient securing buckles, six IV-attachment points, ergonomically designed handles & patented joint configuration with locking leg crossbar

SPECIFICATIONS

Folded: H 21.5 in. x W 8.5 in. x D 7.5 in.
Open: H 84.5 in. x W 22.5 in. x D 4 in.
Weight: 13 lb 8 oz
Static Payload: 1,200 lb (544.3 kg)
Ground Clearance: 1.5 in.




MEDEVAC LITTER

BED COLOR	FRAME COLOR	ITEM #
Black	OD Green	60-0051
Orange	Black	60-0052

OSL™ LITTER *OverSized Litter*

- Safely evacuate larger patients with 30% larger surface area compared to standard-sized evacuation platforms
- Bi-fold configuration enables lower cube storage capability
- 6 litter stirrups, 6 padded carrying loops & 3 patient securing straps for firmer patient support as well as better weight distribution on the floor
- Flame-retardant cover & attached storage pouch made of a durable Nomex/Kevlar-blend
- Certified safe-to-fly by the U.S. Air Force

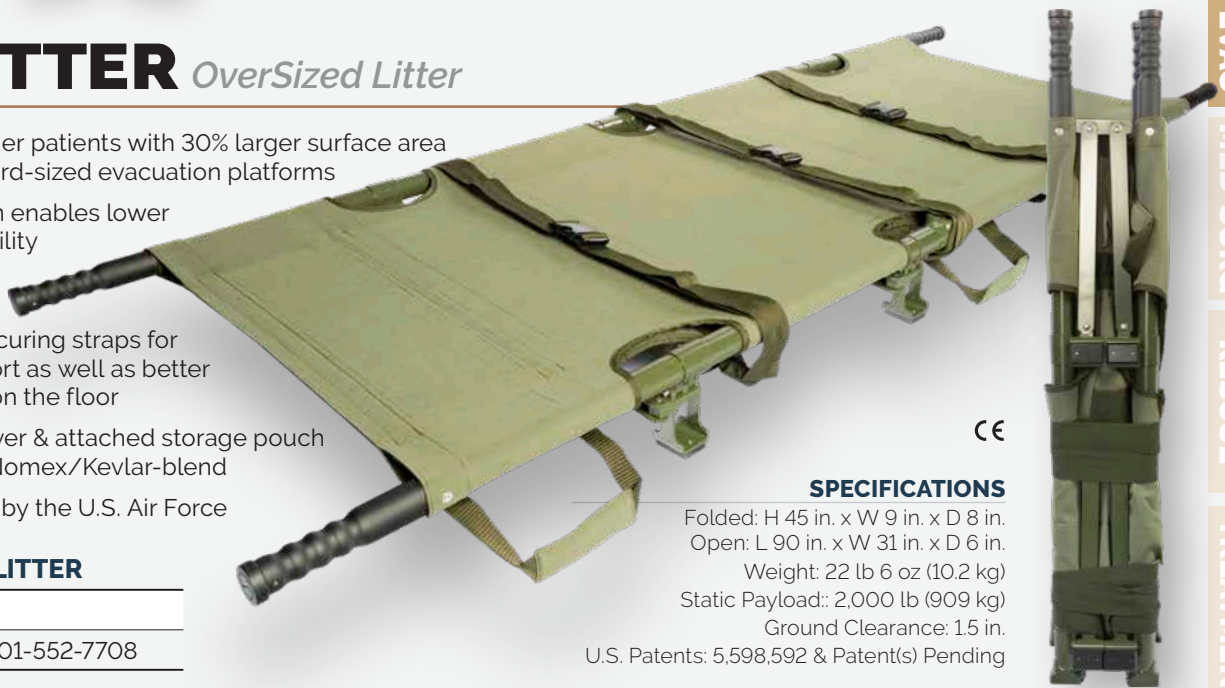

SPECIFICATIONS

Folded: H 45 in. x W 9 in. x D 8 in.
Open: L 90 in. x W 31 in. x D 6 in.
Weight: 22 lb 6 oz (10.2 kg)
Static Payload: 2,000 lb (909 kg)
Ground Clearance: 1.5 in.

U.S. Patents: 5,598,592 & Patent(s) Pending

OSL OVERSIZED LITTER

ITEM#	NSN#
60-0023	6530-01-552-7708

OSL™ LITTER STANDS

for the OverSized Litter™

- Set of 3 lightweight aluminum stanchions accommodate our oversized OSL™ Litter
- Quickly opens and securely locks the litter in place
- Easily folds for transport & storage



SPECIFICATIONS

Folded: H 34 in. x W 25 in. x D 2.5 in.
Extended Height: 32 in.
Base Width: 23.5 in.
Weight: 28 lb 8 oz (12.93 kg)



OSL LITTER STANDS

UNIT	ITEM#
Set of 3	60-0029

OSL™ BACKREST

for the OverSized Litter™

- Position patient at 45° & 90° angles for maximized treatment options & increased comfort
- Constructed of a lightweight aluminum frame with durable Nomex /Kevlar-blend fabric cover
- Compatible with our OSL™ litter



SPECIFICATIONS

H 31.5 in. x W 30.25 in. x D 1.6 in.
Weight: 2 lb 15 oz (1.33 kg)



OSL BACKREST

ITEM#
60-0024



TED *Tactical Extrication Device*

- Designed to allow two people to carry a casualty to safety
- Lighter, more compact non-rigid evacuation platform designed specially for use in high-threat environments & limited spaces
- Constructed of heavy duty vinyl with 4 sturdy carrying handles
- Packaged in a sturdy, transport-friendly bag with adjustable shoulder strap

**SPECIFICATIONS**

Folded: H 5 in. x W 10.25 in. x D 3.5 in.
Open: L 43.5 in. x W 19 in.
Weight: 1 lb 8 oz (0.68 kg)
Payload: 600 lb (272 kg)

TACTICAL EXTRICATION DEVICE

COLOR	ITEM #	NSN #
OD Green	60-0006	6530-01-492-4114
Black	60-0005	6530-01-531-6808



R.E.E.L. *Rapid Expeditionary Evacuation Litter*

- Cost-effective, non-rigid casualty evacuation platform designed for quick and effective casualty evacuation
- Rugged, lightweight, nylon construction
- 8 carry handles accommodate mission teams ranging from 2 to 8 rescuers
- Packaged in a compact, small cube, vacuum sealed pouch

SPECIFICATIONS

Packaged - L 13 in. x W 10 in. x D 0.75 in.
Litter (Open) - L 78 in. x W 30 in. | Weight: 1 lb 3 oz
Carrying Capacity: 700 lb

R.E.E.L. LITTER**ITEM#**

60-0075

NAR **K-9 LITTER** *& Carrying Strap Accessories*

- Lightweight, compact litter for the movement of an injured canine
- Litter can be carried in cargo pocket or pouch
- Leg cut outs to safety carry the canine in the litter
- Six large and four small carry handles that allow for one to four person carries
- Optional K-9 Carry Strap & 2 Carabiners for over the shoulder carry of injured K-9 also available

**SPECIFICATIONS**

Litter: W 35 in. x H 48 in.
Adjustable Carry Strap:
L 9 to 54 in. x W 2 in.
Carabiner open gate
clearance: 27 mm
Packaged:
W 5.5 in. x H 10 in. x D 2.5 in.
Weight: 1.38 lb

K-9 LITTER

OPTIONS	ITEM#
K-9 Litter Only	60-0074W
Carry Strap & Carabiners	ZZ-1004
K-9 Litter with Carry Strap & Carabiners	60-0076

STINGRAY® *Poleless Litter*

- Ideal for decontamination roller systems & moving wet patients
- Non-skid, flame-retardant monofilament polypropylene fabric that is resistant to mildew, acids, alkalis and penetrating liquids
- Tested by Federal & private labs using live warfare agents
- Includes six heavy-duty carrying handles & 2 casualty securing straps
- Easy to store and clean

SPECIFICATIONS

Folded: H 20 in. x W 12 in. x D 2 in.
Open: L 78 in. x W 24 in.
Weight: 13 oz (0.37 kg)
Payload: 1,000 lb (453.6 kg)

**STINGRAY POLELESS LITTER**

ITEM #	NSN #
60-0004	6530-01-492-4079



RESCUE SLEEVE II

- The industry's ultimate patient securing system refined for enhanced stability & ease of use
- Ideal for rescues requiring vertical and/or horizontal hoisting capabilities

SPECIFICATIONS

Folded: L 25.5 in. x W 13 in. x D 8 in.
Open: L 73 in. x W 18 in. x D 2 in.
Weight: 16 lb (7.26 kg)
Payload: 1,000 lb (453.6 kg)



RESCUE SLEEVE II

COLOR	ITEM #	NSN #
Black	50-0011	6530-01-531-6649
Orange	50-0012	6530-01-477-8525

NOTE: Spine board sold separately. Hoisting straps & carabiners not included


NAR CASUALTY IMMOBILIZATION SYSTEM

- Constructed of black polypropylene webbing with hook & loop closures
- Featuring 6 attachment points: 2 "suspender" type, over-the-shoulder straps & 4 adjustable torso/limb straps with adjustable center strap
- No choking "V" configuration by the patient's neck



CASUALTY IMMOBILIZATION SYSTEM

ITEM #	NSN #
50-0031	6530-01-586-2411

SPECIFICATIONS 
Pkg: L 16.5 in. x W 4 in. x D 2.75 in.
Weight: 1 lb 1 oz

NAR CASUALTY RESTRAINT STRAPS

Heavy-duty nylon webbing & convenient seat belt buckle utility

- A. 2 piece Plastic Buckle**
L 5 ft x W 2 in. • Weight: 4 oz
Black: 50-0023 • Orange: 50-0024
- B. 2 piece Plastic Buckle/Speed Clip**
L 5 ft x W 2 in. • Weight: 7.8 oz
Black: 50-0026 • Orange: 50-0025
- C. 2 piece Metal Buckle**
L 5 ft x W 2 in. • Weight: 8 oz
Black: 50-0027 • Orange: 50-0028
- D. 2 piece Metal Buckle/Speed Clip**
L 5 ft x W 2 in. • Weight: 11.7 oz
Black: 50-0029 • Orange: 50-0030



SPINE BOARD

- Thin & lightweight with tapered design
- Surface texture to prevent slippage & 14 raised, oversized hand-holds
- Black spine board includes ten rapid attachment pins

SPECIFICATIONS

L 72 in. x W 16 in. x D 1.75 in.
Weight: 14 lb (6.35 kg) Payload: 1,100 lb (498.96 kg)
Float Capacity: 250 lb (113.4 kg)

SPINE BOARD

COLOR	ITEM #	NSN #
OD Green	50-0013	6530-01-490-2487
Orange	50-0014	6530-01-531-6771
Black w/Pins	50-0021	----



NAR HEAD IMMOBILIZER

- Quick, effective head block constructed of soft, closed cell foam with hook & loop closures and padded straps
- Ears and head are accessible without removing the device



SPECIFICATIONS

H 10 in. x W 16 in. x D 5.75 in.
Weight: 2 lb 3 oz

HEAD IMMOBILIZER

ITEM #	NSN #
50-0020	6530-01-619-1777

CERVICAL COLLAR ACE®

- Adjustable neck support with foam lining
- Extended back to ensure complete
- Adjusts to fit 16 adult neck sizes

SPECIFICATIONS

H 6.75 in. x W 20.5 in. x D 0.75 in.
Weight: 5.2 oz

ACE CERVICAL COLLAR

COLOR	ITEM #	NSN #
OD Green	50-0010	6515-01-541-8147
White	50-1007	6515-01-650-0036





MED•SURG

MEDICAL | SURGICAL PRODUCTS

- DIAGNOSTICS
- FIELD HOSPITAL
- SURGICAL SUPPLIES

see pages 62-73





VSCAN™ with DUAL PROBE Pocket-sized Ultrasound Tool from GE

- Directional Color Flow, Voice Recorder, Auto Image Optimization
- Intuitive thumb-wheel navigation with system presets helps speed diagnostic decisions and enables efficient patient triage
- Access to Vscan™ Web Portal / Vscan™ Gateway Software for PC
- Storage capacity: 4GB upgradable to 32GB on microSD or microSDHC cards
- Data storage formats: .jpg for still frames | .mpg for loops | .wav for voice

• Clinical applications:

- Cardiac
- Abdominal
- Renal
- Obstetrics | gynecology
- Urology
- Fetal
- Evaluation of presence of fluid
- Peripheral vascular
- Procedure guidance for arterial and venous vessels
- Imaging guidance for needle/catheter placement
- Thoracic/Lung
- Thyroid and other small organs
- Long bone
- Hip and knee joints
- Pediatric



VSCAN ULTRASOUND

ITEM

20-0092



KIT CONTENTS:

- 1 x GE Vscan™ with Dual Probe Pocket-sized Ultrasound Tool
- 1 x Docking Station with Global AC Adaptor & USB cable
- 1 x Rechargeable Battery
- 1 x 4 GB microSD card
- 1 x Soft Case
- 1 x Vscan™ with Dual Probe User's Manual



SPECIFICATIONS

- Display unit :
H 5.32 in. x W 2.87 in. x D 1.1 in.
(135 x 73 x 28 mm)
- Probe: H 5.08 in. x W 1.5 in. x D 1.1 in.
(129 x 39 x 28 mm)
- Weight (unit and probe):
0.96 lb (436 g)
- Display: 3.5 in.
- Image Resolution:
240 x 320 pixels



MIGHTYSAT™ Rx

Fingertip Pulse Oximeter from Masimo

features Masimo SET® Measure-through Motion and Low Perfusion™ Pulse Oximetry which overcomes the challenges of low perfusion and motion artifact that limit conventional pulse oximetry.

This Signal Extraction Technology® (SET®) uses five parallel signal processing engines that measure through challenging conditions allowing users to accurately measure oxygen saturation and pulse rate even during movement or low blood flow.

- **High Resolution Display**
Rotating and adjustable color screen for real-time display of all measurements
- **Plethysmographic Waveform**
Provides high resolutions display of the patient's pleth waveform
- **Touchpad**
Provides easy access of settings and device features
- **Durable, Lightweight Design**
Small and portable for use in clinical and home-care environments
- **Flexible Pad**
Easily conforms to finger shape for enhanced comfort & secure application to finger to minimize external interference
- **Signal IQ (SIQ™)**
An indicator of the SpO₂ measurement confidence
- Bluetooth Models (20-0094 & 20-0095) are compatible with free downloadable Masimo Professional Health App

KIT CONTENTS:

- 1 x MightySat™ Rx Device
- 1 x Lanyard
- 2 x AAA Batteries
- 1 x Carrying Case

MIGHTYSAT PULSE OXIMETER

OPTION	ITEM #
Standard	20-0093
Bluetooth LE	20-0094
Bluetooth LE/RRp/PVi	20-0095

SPECIFICATIONS

Device Dimensions :
H 2.9 in. x W 1.6 in. x D 1.2 in.
(7.4 cm x 4.1 cm x 3.0 cm)
Weight (with batteries):
2.56 oz (73 g)
Display: 3.3 cm



FEATURES	20-0093	20-0094	20-0095
Masimo SET® Measure-through Motion and Low Profusion™ Technology	✓	✓	✓
Measurements			
SpO ₂ (Oxygen Saturation)	✓	✓	✓
PR (Pulse Rate)	✓	✓	✓
PI (Perfusion Index)	✓	✓	✓
PVI (Pleth Variability Index)			✓
Display			
Pleth Waveform	✓	✓	✓
3.3 cm Color Display	✓	✓	✓
Rotational Screen	✓	✓	✓
Bluetooth LE		✓	✓
Compatible with Masimo Professional Health App*		✓	✓
Four Year Limited Warranty	✓	✓	✓

INTRO
M.A.R.C.H.
KITS
T.A.C.E.V.A.C.
M.E.D.-S.U.R.G.
R.E.S.C.U.E.
T.R.E.A.T.M.E.N.T.S.U.P.P.O.R.T.
T.R.A.I.N.I.N.G.
I.N.D.E.X.



FIELD CORPSMAN KITS

Portable ENT diagnostics

- Lightweight, durable ENT diagnostic sets for portable field care of the eyes, ears & nose
- A wide selection of high intensity lights & magnifiers for enhancing both day and night examinations

FIELD CORPSMAN KITS:



Additional items in Deluxe Kit only



SPECIFICATIONS

H 8 in. x W 4.25 in. x D 2.25 in.
 Weight: Basic Kit: 6.2 oz
 Deluxe Kit: 11.6 oz
 Ophthalmoscope only: 3.7 oz

FIELD CORPSMAN KITS

KIT OPTIONS	ITEM #	NSN #
Basic	20-0001	6515-01-458-6057
Deluxe (with Heine)	20-0002	6515-01-458-6178
Ophthalmoscope only	20-0007	6540-01-491-6651

BP/STETHOSCOPE COMBO KIT *from NAR*

- Ultra-compact, multi-pocketed BP kit that can be stored in most trauma kits
- Blood-pressure cuffs in 4 sizes (adult, large adult, child, infant).
- ADC stethoscope with extra-large bell & ultra-sensitive diaphragm
- Ambidextrous palm gauge

BP/STETHOSCOPE COMBO KIT

ITEM #	NSN #
20-0057	6515-01-542-7723

SPECIFICATIONS

H 6.75 in. x W 9 in. x D 5 in.
 Weight: 2 lb



HEINE MINI 3000 *Ophthalmoscope*

- High quality optics for examining the interior of the eye - in a pocket-size
- Enhanced XHL Xenon Halogen Technology; concentrated white light for perfect illumination & crisp retinal detail



HEINE MINI 3000

ITEM #	NSN #
20-0007	6540-01-491-6651

SPECIFICATIONS

H 6.75 in. x W 1.25 in. x D 0.9 in.
 Weight: 4 oz

STETHOSCOPE *from ADC*


- Best full-range stethoscope for auscultation & comfort
- Superior acoustics for accurate auscultation of heart, lung & blood-pressure sounds

SPECIFICATIONS

Device Length: 30.5 in.
 Weight: 13.5 oz

STETHOSCOPE

ITEM #
20-0060

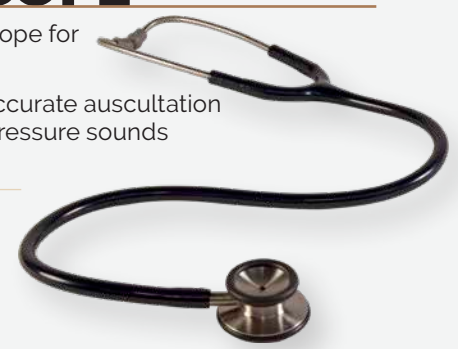



ARMADILLO Medication & Storage Case

- Compact solution for keeping critical medications safe, secure and easily accessible
- Features Automatic Pressure Equalization Valve
- Hardware constructed of Marine grade 316 SS
- Case is water-resistant and crush-proof
- Hook & loop attachment system allows storage compartment to be tailored to hundreds of custom configurations

SPECIFICATIONS

Ext: L 7.95 in. x W 5.1 in. x D 2.14 in.
Int: L 6.45 in. x W 3.87 in. x D 1.75 in.
Weight: 10.5 oz

ARMADILLO

* NATO NSN

LID STYLE	ITEM #	NSN #
Black	ZZ-0081	8145-01-573-2533
Clear	ZZ-0259	6545-01-594-9399*

also available: Half Inch Hook Discs with pressure sensitive adhesive Set of 250 • Item# ZZ-0258 • NSN# 8145-01-617-9075*



Clear Lid
ZZ-0259 only
includes Loop Liner
on the bottom

I&A CASE Injury & Ailment Case

- Superbly engineered medical equipment organization system
- Lid's clear pockets and hook & loop fasteners allow secure storage as well as superior inventory management
- Polymer O-ring provides dust and water-proof seal

- Automatic Gore-TEX® pressure equalization valve stops moisture from entering the case – also prevents vacuum-lock
- Double-throw latches engineered to tolerate up to 400 lb of pressure
- All-weather, sure-grip handle and reinforced padlock protectors

SPECIFICATIONS

H 8.5 in. x W 20.5 in. x D 17 in.
Wgt: 13 lb 12 oz



I&A CASE

ITEM

ZZ-0222

M.U.L.E. Multi-Use Luggage & Equipment Carrier

- Heavy-duty, water-resistant nylon platform with large compartment and mesh pockets on the inside lid for maximum storage of critical items
- Outside end pockets feature open channels for expanded storage capability
- Suitcase-style carrying handles on top and side; Telescoping handle and wheels for easy transport
- Hidden rucksack shoulder pads



SPECIFICATIONS

H 30.5 in. x W 18.5 in. x D 11 in.
Wgt: 11 lb 10 oz | Payload: 70 lb

M.U.L.E.

ITEM

ZZ-0212

M.U.L.E. II Multi-Use Luggage & Equipment Carrier II

- Extra large, heavy-duty, water-resistant 1000 denier nylon platform featuring 9,000 cubic in. of storage with large reinforced open compartment and multiple zippered mesh pockets for maximum storage
- Top & side web carry handles and 3 durable 3.75 in. diameter wheels for easy transport



Available in Black & Coyote

M.U.L.E. II

COLOR	ITEM #
Black	ZZ-1157T
Coyote	ZZ-1158T

SPECS

Exterior: H 40 in. x W 21 in. x D 15 in.
Interior: H 40 in. x W 15 in. x D 15 in.
Weight: 15.4 lb
Capacity: 9,000 cubic in.



INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

PAXLIGHT® LED Field Medical Light

- Compact, lightweight precision lighting solution with the rugged engineering required for use in field hospitals
- LED Technology with intense light fifteen times brighter than traditional halogen field lights
- Patented 3-pod head eliminates shadow interference common with single point light sources
- Energy efficient, cool white color with an exceptionally long life
- Stays cool to the touch, projecting no heat on the patient
- Locking joints make it fully adaptable and secure with adjustable mounting clamp for attachment to hospital bed, litter or litter stand
- Includes a rugged storage case with wheels for easy transport



SPECIFICATIONS

Light in case:
L 47.2 in. x W 16.5 in. x D 6.7 in.

Light Head: 12.8 in. Diameter

Light Arm/Post:
L 78 in. x 2.8 in. Diameter

Weight: Case & Light: 37 lb 8 oz
Light with Arm/Post:
11 lb 8 oz

Light Output: 50,000 lux
(4590 foot candles)

PAXLIGHT FIELD LIGHT

ITEM #	NSN #
60-0042	6530-01-549-1089

SPECIFICATIONS

Base: W 38 in.
Assembled: H 38 in.
Weight: Base: 13.5 lb
Assembled: 22 lb

WHEELED BASE

ITEM #
60-0045



SPECIFICATIONS

L 4.5 in. x W 1 in. x D 1 in.
Shield Diameter: 4.25 in.

EZ LIGHT HANDLES

UNIT	ITEM #
Pack of 10	60-0046



MARK IV Field Hospital Bed

- Lightweight solution designed for patient comfort in the combat support or field hospital/treatment tent
- Durable & easy-to-open, supporting up to 850 lb of static weight
- Designed to be rapidly deployed with minimal training in less than 60 seconds by a single person
- Constructed of a bi-fold aluminum frame that is adjustable to 6 positions
- 3 in. polyurethane foam pad with water-resistant decontaminable cover easily cleans with 10% bleach solution or other standard hospital cleaning products
- Includes IV Pole and a removable table platform that rests on its adjustable bed rails

SPECIFICATIONS

Folded: L 42 in. x W 17 in. x D 32 in.
Opened: L 79.25 in. x W 33 in. x H 35.25 in.

Weight:
Mattress and Frame: 54 lb
Mattress, Frame & Accessories: 68 lb
U.S. Patents Pending

MARK IV FIELD HOSPITAL BED

ITEM #	NSN #
66-0019	6530-01-624-7626



microMEND™ Skin Closure Devices

- Easy-to-use alternative to sutures or surgical staples when skin closure is required
- Safe, secure, and rapid skin closure
- Equivalent holding strength to sutures
- Patent pending microstaples anchor the device to the skin
- Three times faster to apply than sutures
- More comfortable and easier to remove than standard invasive methods of skin closure
- Can be used on a variety of lacerations and incisions and with wounds of various lengths, as the product can be applied right next to each other

SPECIFICATIONS

Large

Device:
43mm x W 12mm
Box Wt.: (24) 3.5 oz

Small

Device:
L 32mm x W 12mm
Box Wt.: (24) 3.5 oz
(48) 7.1 oz

X-Small

Device:
L 28mm x W 7mm
Box Wt.: (24) 3.5 oz

Wide

Device:
L 38mm x W 28mm
Box Wt.: (24) 6.9 oz



X-SMALL: Ideal for small areas, such as the fingers, toes, ears, and nose



SMALL: Broadest application for closing most wounds. It is used for closing most lacerations (face, neck, and hands)



LARGE: For closing larger wounds on the extremities, back, chest, and abdomen. Great for closing hard to repair skin tears

microMEND™ Wide

- Three times wider for easier and faster closure of longer incisions with the same holding strength as other microMend sizes
- Ultra-flexible and elastic making it ideal for use in joint surgeries, Cesarean sections, and other procedures in areas of increased movement and varying skin architecture
- Available with apertures (Wide A) to allow drainage, or without apertures (Wide S) for a solid seal to keep the wound clean from contaminants




STEP 1:




Remove device from packaging by pulling **TAB** upwards

STEP 2:

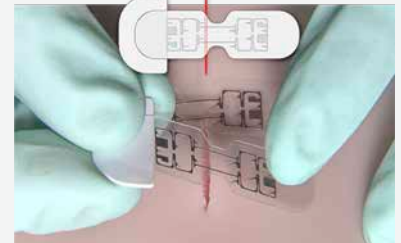

Position device so the edge of narrow section on the **TAB**-side is over wound. Press **MICROSTAPLE**-side onto skin.

STEP 3:

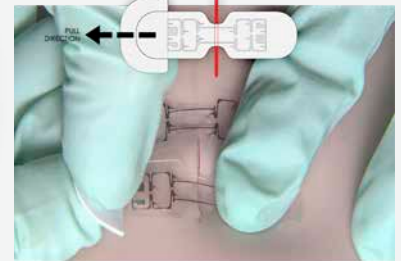

Pull **TAB**-side across the wound. Press **MICROSTAPLE** of tab-side down onto skin. Remove the tab with a twist.

microMEND DEVICES

DEVICE	UNIT (per box)	ITEM #
Large	24	30-0203
Small	24	30-0230
Small	48	30-0249
X-Small	24	30-0252
Wide A	24	30-0295
Wide S	24	30-0296

U.S. Patents Pending

Available with Aperture or Solid Seal



R-SICK Rapid Surgical Intervention Convenience Kits

- Flexible austere and prolonged field care surgical instrument storage and recovery system tray
- For both major and minor surgical procedures in less than ideal or highly mobile environments
- Sequentially set up for rapid "incision, spreading, damage control, and closure"
- Ideal for complex, emergent procedures such as surgical airways, hemorrhage control, field amputation, and REBOA



Instruments are available in Floor & OR Grade

SPECIFICATIONS

Packaged: H 8 in. x W 7 in. x D 3.5 in.
 Fully deployed (flat): L 24 in. x W 43 in. x D 0.1 in.
 Weight: Kits 3 lb • Bag Only 11 oz

R-SICK CONTENTS:

- 2 x Weitlaner Retractor, 3 x 4 Blunt (5.5 in.)
- 1 x Trousseau Tracheal Dilator (5.25 in.)
- 2 x Adson Forceps, 1 x 2 (4.75 in.)
- 1 x Ewald Forceps 1 x 2
- 1 x Plain Splinter Forceps (4.5 in.) Serr. Handle
- 1 x Kelly Forcep Heavy
- 1 x Converse Needle Holder (4.24 in.) Long
- 1 x Olsen-Hegar Needle Holder (7.25 in.) Serr
- 1 x Stevens Tenotomy Scissors Trbl/Bl
- 1 x Stevens Tenotomy Scissors Cvdbl/Bl
- 1 x Scissors, Mayo (5.5 in.) Str
- 1 x Mayo Dissecting Scissors, Crv (5.5 in.)
- 1 x Knife Handle (#3)
- 2 x Gigli Saw Blade (20 in.)
- 2 x Gigli Saw Handle (3.5 in.) Loop
- 1 x Hiupp Tracheal Hook, (6.5 in.)
- 2 x Forcep, Rochester-Pean (7.25 in.) Curv
- 6 x Scalpel #10 w/Handles

- 6 x Scalpel #11
- 6 x Scalpel #15
- 6 x Blade, #10 Pluss/S(Tefl)
- 10 x Blade, #15 Pluss/S(Tefl)
- 2 x Suture Ethilon 4-0 Blk Mono Fs-218
- 2 x Suture 3-0 Silk
- 2 x Black Braided Silk Suture, 2-0/2-60 in.
- 2 x 2-0 Suture Ct-2
- 4 x 4-0 Silk Suture Blk (30 in.)
- 1 x Sharps Shuttle
- 1 x Farabeuf Retractor



R-SICK

INSTRUMENTS	ITEM #
Floor Grade	83-0016
OR Grade	83-0017
Bag Only	80-0871



FBTK *Field Blood Transfusion Kits*

FBTK Basic:

Donor Module:
 5 x Alcohol Prep Pad
 2 x 2-pack Gauze Pad 4 in. x 4 in.
 1 x IV Constricting Band - 1 in. x 18 in.
 1 x Clear Adhesive Dressing - 4 in. x 4.75 in.
 1 x Blood Bag Label

Recipient Module:

5 x Alcohol Prep Pad
 1 x IV Catheter - 18 ga x 1.25 in.
 2 x 2-pack Gauze Pad - 4 in. x 4 in.
 1 x IV Constricting Band - 1 in. x 18 in.
 1 x Clear Adhesive Dressing - 4 in. x 4.75 in.
 1 x IV Catheter - 16 ga x 1.25 in.

Additional Supplies:

1 x roll 1 in. Surgical Tape
 2 x pair Black Nitrile Gloves, Large
 1 x Black Permanent Marker
 1 x Kelly Curved Hemostat - 5.5 in.
 1 x Single Blood Bag (450ml)
 1 x 550 Cord - 9.5 in.
 1 x Beaded Security Tie - 11 in.
 1 x IV Set Y-Type with Filter 170-260 UM
 1 x Questionnaire Card - DD572 Modified

FBTK Blood Surveillance Module:

Donor Module:

1 x IV Hypodermic Needle - 18 ga x 1.5 in.
 1 x IV Syringe, 10cc
 1 x Eldon Card
 1 x Interlink Injection Site Needle
 1 x Blood Transfer Device
 2 x Vacutainer, Gold Top (5ml)
 1 x Vacutainer, Purple Top (5ml)
 2 x Vacutainer, Pearl Top (5ml)

Recipient Module:

1 x 3 Way Stopcock
 1 x Eldon Card
 1 x Interlink Injection Site Needle

*FBTK Complete:

Contains all components of
FBTK Basic & Blood Surveillance Module

FIELD BLOOD TRANSFUSION


OPTIONS	ITEM #
FBTK Basic	83-0022
Blood Surveillance	83-0025
FBTK Complete*	83-0000



MINOR SURGERY *Sets*

MINOR SURGERY Contents:

1 x Tray, Shallow
 2 x Kelly Forceps, Straight 5.5 in.
 2 x Kelly Forcep, Curved 5.5 in.
 1 x Mayo Hegar Needle Holder, 6 in.
 1 x Adson Forceps, Serr 4.75 in.
 1 x Forceps Tissue, 5 in.
 1 x Collier Needle Holder, Straight 5 in.
 1 x Handle Surgical Knife #3
 1 x Probe W. Eye 5 in.
 1 x O.R. Scissors, Sb St 5.5 in.
 1 x Wrap, Kimguard 24 x 24 in.

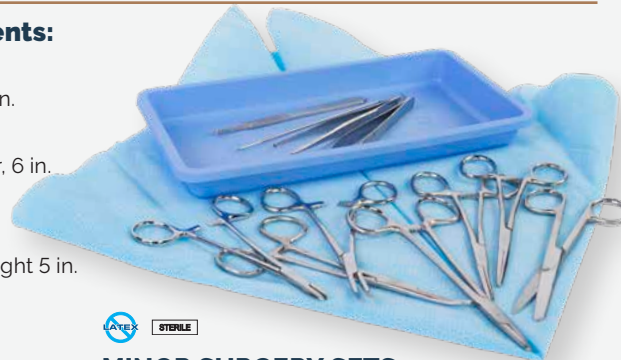

MINOR SURGERY SETS

ITEM #	UNIT	NSN#
ZZ-0960	Single	-----
ZZ-0959	Pack of 20	-----

SPECIFICATIONS

Packaged: L 9 in. x W 5 in. x D 1 in.
 Weight: 5 oz

SURGICAL SUTURE *Kits*

SURGICAL SUTURE Contents:

1 x Tray, Shallow
 1 x Adson Forceps w/Teeth
 1 x Webster Needle Holder
 1 x Iris Scissors, Straight
 1 x Mosquito Hemostat, Curved
 1 x Poly Bag w/Twist Tie
 1 x Poly-backed Towel
 1 x Fenestrated Drape
 1 x Alcohol Prep Swab
 1 x PVP Prep Swab
 1 x Cotton Filled Sponge.



Sutures
 not included



SURGICAL SUTURE KITS

ITEM #	UNIT	NSN#
ZZ-0958	Single	-----
ZZ-0961	Pack of 20	6515-01-311-0360

SPECIFICATIONS

Packaged: L 9 in. x W 5 in. x D 1 in.
 Weight: 6 oz

TOENAIL REMOVAL *Trays*

TOENAIL REMOVAL Contents:

1 x Tray, Half Deep
 1 x CSR Wrap, Kimguard 24 in. x 24 in.
 1 x Kelly Forcep, Straight 5.5 in.
 1 x Kelly Forcep, Curved 5.5 in.
 1 x Adson Forceps, 1 x 2 Teeth, 4.75 in.
 1 x Iris Scissors, Straight 4.5 in.
 1 x Cuticle Nipper 4 in.
 1 x Probe w/Eye 5.5 in.
 10 x 12-Ply Gauze, 4 in. x 4 in.
 1 x Kwiktear Header Bag, 10 in. x 15 in.



TOENAIL REMOVAL TRAYS

ITEM #	UNIT	NSN#
ZZ-0956	Single	-----
ZZ-0957	Pack of 10	6515-01-630-8050

SPECIFICATIONS

Packaged: L 9.72 in. x W 5.4 in. x D 2 in.
 Weight: 11 oz



INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

SUTURE REMOVAL KIT #1



QTY 1 KC 100 Kimguard



QTY 1 Littauer Scissors



QTY 1 Adson Forceps



QTY 1 Gauze Pad, 3 x 3

ITEM #

ZZ-0664

SPECIFICATIONS

STERILE

Pkg: L 6 in. x W 3.5 in. x D 0.38 in.
Weight: 2 oz

SUTURE REMOVAL KIT #2



QTY 1 KC 100 Kimguard



QTY 1 Crvd Iris Scissors




QTY 1 Adson Forceps



QTY 1 Gauze Pad, 3 x 3



QTY 1 Alcohol Prep Pad

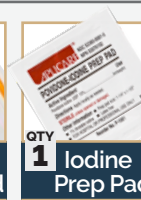

QTY 1 Iodine Prep Pad

ITEM #

ZZ-0665

SPECIFICATIONS

STERILE

Pkg: L 6 in. x W 3.5 in. x D 0.38 in.
Weight: 2 oz

SUTURE REMOVAL KIT #3



QTY 1 KC 100 Kimguard



QTY 1 Str Iris Scissors



QTY 1 Adson Forceps



QTY 1 Gauze Pad, 3 x 3



QTY 1 Alcohol Prep Pad

ITEM #

ZZ-0666

SPECIFICATIONS

STERILE

Pkg: L 6 in. x W 3.5 in. x D 0.38 in.
Weight: 2 oz

SUTURE REMOVAL KIT #4



QTY 1 KC 100 Kimguard



QTY 1 Littauer Scissors



QTY 1 Adson Forceps



QTY 1 Alcohol Prep Pad



QTY 1 Iodine Prep Pad

ITEM #

ZZ-0667

SPECIFICATIONS

STERILE

Pkg: L 6 in. x W 3.5 in. x D 0.38 in.
Weight: 2 oz

SURGICAL INSTRUMENTS

Precision standard surgical instruments available individually

STERILE

ADSON FORCEPS

FP, 1x2 • 4.75 in.




ITEM #

ZZ-0705

HALSTEAD MOSQUITO

Curved • 5 in.

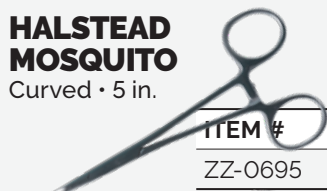

ITEM #

ZZ-0695

HALSTEAD MOSQUITO

Straight • 5 in.




ITEM #

ZZ-0694

IRIS SCISSORS

Curved • 4.5 in.

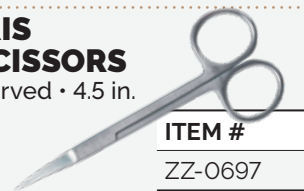


ITEM #

ZZ-0697

IRIS SCISSORS

Straight • 4.5 in.

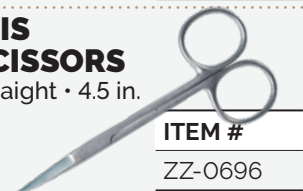

ITEM #

ZZ-0696

KELLY HEMOSTAT

Curved • 5.5 in.



ITEM #

ZZ-0701

KELLY HEMOSTAT

Straight • 5.5 in.



ITEM #

ZZ-0700

LISTER SCISSORS

5.5 in.



ITEM #

ZZ-0698

O.R. SCISSORS

SB Straight 5.5 in.



ITEM #

ZZ-0699

ROCHESTER PEAN

Hemostat, Curved 6.25 in.



ITEM #

ZZ-0702

WEBSTER NEEDLE HOLDER

Holder, 5 in.




ITEM #

ZZ-0703




FACIAL LACERATION TRAY

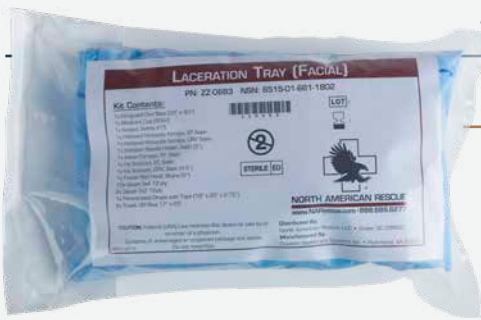

- Cost effective, sterile, single use kit for facial lacerations with a 3 year shelf life
- Conveniently packaged in a peel-to-open pouch for easy access


ITEM #	NSN#
ZZ-0683	6515-01-661-1802

SPECIFICATIONS STERILE
 Pkg: L 9 in. x W 5 in. x D 2.75 in.
 Weight: 12 oz

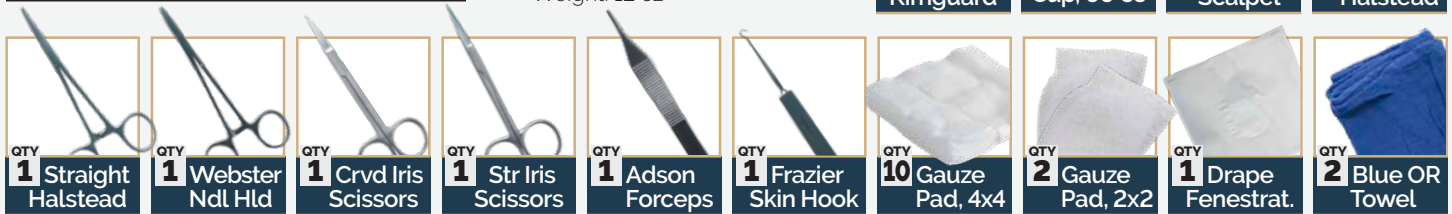

MINOR LACERATION TRAY

- Cost effective, sterile, single use kit for facial lacerations with a 3 year shelf life
- Conveniently packaged in a peel-to-open pouch for easy access

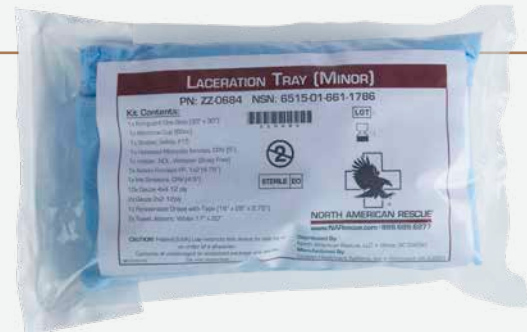

ITEM #	NSN#
ZZ-0684	6515-01-661-1786


SPECIFICATIONS STERILE
 Pkg: L 9 in. x W 5 in. x D 2.75 in.
 Weight: 12 oz



- Cost effective, sterile, single use kit for facial lacerations with a 3 year shelf life
- Conveniently packaged in a peel-to-open pouch for easy access to components

ITEM #
ZZ-0685

SPECIFICATIONS STERILE
 Pkg: L 9 in. x W 5 in. x D 2.75 in.
 Weight: 12 oz

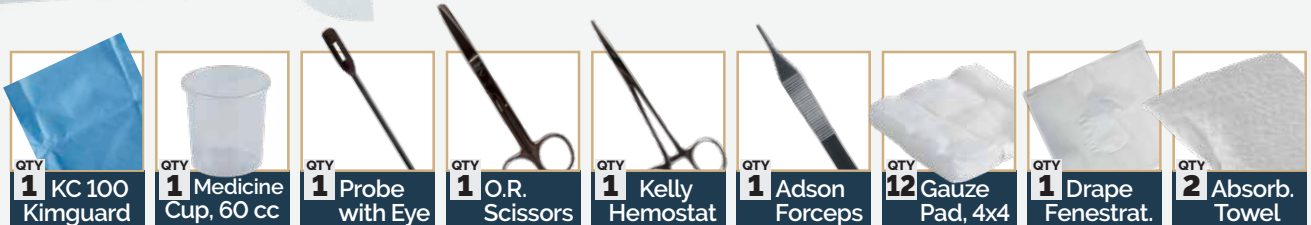

Custom Procedure Trays/Kits available • Contact North American Rescue for information






RESCUE

EQUIPMENT & ACCESSORIES

see pages 74-85

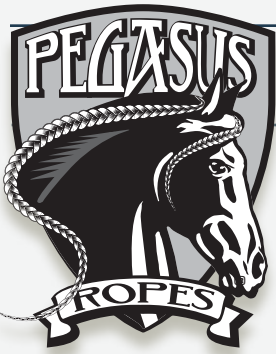

PEGASUS ROPE®



8MM HAVOC Assault Rope

Innovative Synergy of 2 unique Fibers engineered for unmatched Confidence & Resilience from Ropes in Challenging Environments

North American Rescue & Edelrid, the industry leader in rope manufacturing with 150 years of experience, have joined forces with multi-disciplinary vertical specialists to innovate the **most tenacious small diameter rope on the market today** – engineered to counter the most severe environmental pathologies found throughout full-spectrum operations... whether at work or at play.

The unique material engineering of the 8mm HAVOC core, specifically it's **Innegra™ Braided Fiber Core Integration**, provides **superior durability, low fatigue** properties, an **extremely low specific gravity**, and **superior resistance to the harshest contaminants** (acids, bases, and organic solvents). The novel tensioning and fiber refinement techniques utilized within the **12-braid Innegra™ core amplifies its overall robustness** and resilience providing an operational rope which acts as a counter-measure to anticipated environmental challenges.

This specialized core of the 8mm HAVOC is protected by a higher density, precision double braided **Technora sheath** to synergistically interface with the core for **optimized strength, elongation, decreased creep, heat resistance**, and **operational safety**. To reinforce the **abrasion resistance, cut resistance**, and **overall durability**, the 8mm HAVOC's sheath was bulked up to a **52% sheath ratio**.

EDELRID

in partnership
with



NORTH AMERICAN RESCUE®

Engineered to counter the pathologies innate to the following environments and/or mission sets:

- **Mountain/Rough Terrain Movement & Evacuation**
- **Urban Access, Extraction, & Rescue**
- **Waterborne/Maritime Access & Rescue**
- **Cave/Subterranean Response**
- **Snow/Ice Rescue**
- **Shipboard Response, Rescue, & Evacuation**
- **Canyon Movement & Rescue (Class B & C)**
- **Confined Space Rescue**
- **Wilderness SAR**
- **Wildland REMS**



Wreak HAVOC with Confidence!



INNEGRA
CORE



**52%
Technora®
Sheath**

CORE

- **Innegra™** Lightest Weight Structural Fiber on the market with a very Low Specific Gravity
- Hydrophobic
- Resistant to Harsh Contaminants
- Low Flex Fatigue
- Superior Durability
- Meets EN 564

8MM HAVOC ASSAULT ROPE

LENGTH	ITEM #
75 ft	90-0403
100 ft	90-1402
150 ft	90-0404
200 ft	90-2402
300 ft	90-3402
600 ft	90-6402

Specifications:

- Diameter: 8 mm
- Weight: 2.9 lb/100 ft
- Elongation: 300 lb: 3.58%
- 3-Sigma MBS: 5,615 lb (25 kN)
- Fibers (Core/Sheath): Innegra™/Technora®
- Sheath Ratio (by wt.): 52.5%

SHEATH

- **52% Technora Sheath** Ratio (by weight)
The most robust on the market!
- Resistant to Cuts, Abrasion & Heat
- Optimized Strength
- Decreased Creep
- Increased Durability & Safety

INTRO

M.A.R.C.H

KITS

TAC·EVAC

MED·SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX



WARHORSE™ SERIES

100% Innegra™ fiber braided core

- High quality, multipurpose rope
- Static rope provides reliability, high performance and great handling
- Excellent knot-ability with high strength-to-weight ratio
- Low stretch/sheath slippage – wet, dry or in any other climatic condition
- Ultra-lightweight, strong Innegra™ core for reliable strength, durability & flotation
- Properties remain unchanged even when wet, dry or contaminated with chemicals
- Low dielectric and UV stable



SPECIFICATIONS

Diameter	9.1 mm	9.5 mm
Weight:	4.88 lb/100 ft	3.34 lb/100 ft
Tensile Strength:	5,670 lbf (22.7 kN)	5,100 lbf (22.7 kN)
Fibers (core/sheath):	Innegra™/Polyester	Innegra™/Polyester
Sheath Ratio:	63% Polyester	28% Polyester
Ratings:	NFPA 1983/2012 Throwline (Pending)	NFPA 1983/2012 Throwline (Pending)

ORDERING INFORMATION

DIAMETER	COLOR	RUNNER	ITEM #
9.1 mm	Black	Gray	90-0207
9.5 mm	Black	Gray Reflective	90-0206
9.5 mm	Yellow	Blue Reflective	90-0205
9.5 mm	White	Blue	90-0208





OLYMPUS™ SERIES

Engineered for rough terrain

- Ideal for use in challenging and abusive environments
- Durable, lightweight and great handling offers reliability in all weather conditions
- Technora sheath adds strength and abrasion resistance
- Buoyant even with a 62% Technora sheath (9.5mm)

ORDERING INFORMATION

DIAMETER	COLOR	RUNNER	ITEM #
9.5 mm	Gold	Black/Red	90-0401
11 mm	Gold	Black/Blue	90-0400

SPECIFICATIONS

Diameter	9.5 mm	11 mm
Weight:	3.82 lb/100 ft	4.99 lb/100 ft
Tensile Strength:	5,495 lbf (24.4 kN)	6,834 lbf (30.4 kN)
Fibers (core/sheath):	Innega™/Technora	Innega™/Technora
Sheath Ratio:	62% Technora	60% Technora



STALLION™ SERIES

Low visibility - high performance

- 100 ft of this 9.5 mm rope offers 5,100 lb of strength under 3 lb – wet or dry
- Chemically resistant to acids, bases & common industrial chemicals
- Low sheath slippage, excellent strength-to-weight ratio and hydrophobic

ORDERING INFORMATION

DIAMETER	COLOR	RUNNER	ITEM #
9.5 mm	Black	Blue/Reflective	90-0103
9.5 mm	Black	Blue	90-0104

SPECIFICATIONS

Diameter	9.5 mm
Weight:	2.9 lb/100 ft
Tensile Strength:	5,100 lbf (22.7 kN)
Fibers (core/sheath):	Innega™/Innega™
Sheath Ratio:	24% Innega™

TRITON™ SERIES

Specially designed for water rescue

- 10% lighter & 250% stronger than polypropylene alone
- 100% Innega™ fiber braided core, the lightest structure fiber available, provides additional buoyancy, low elongation and consistent strength

SPECIFICATIONS

Diameter	7.5 mm
Weight:	2 lb/100 ft
Tensile Strength:	3,050 lbf (13.6 kN)
Fibers (core/sheath):	Innega™/Polypropylene
Sheath Ratio:	46% Polypropylene
Ratings:	NFPA 1983/2012 Throwline (Pending)

ORDERING INFORMATION


DIAMETER	COLOR	RUNNER	ITEM #
7.5 mm	Yellow	Blue Reflective	90-0303
7.5 mm	Red	Yellow	90-0305





R.A.C.E. SYSTEM Rapid Access Casualty Extraction

- Lightweight, multiple use casualty evacuation system enabling first responders/medical providers to ascend, descend, hoist & lower both rescuer & casualty to safety
- Combines 150 ft of lightweight durable 9.5 mm Olympus Pegasus Rope® with the latest in lightweight multipurpose hardware
- Capable of rigging 2:1, 3:1, 4:1, 5:1, 6:1, and 9:1 systems in different configurations for a wide variety of rescuer options with fully illustrated instruction booklet
- Front pocket allows for optional storage
- Carriage options include 4 Point, D-ring configurable modified side and/or backpack carry-style adjustable strap

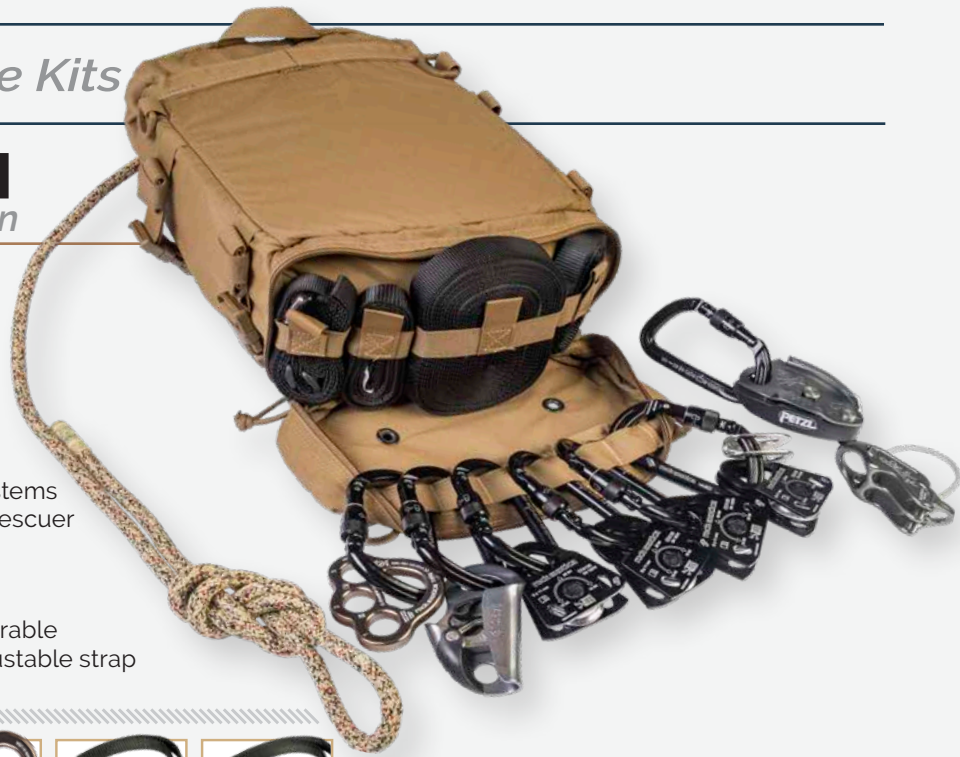

R.A.C.E. SYSTEM CONTENTS:



R.A.C.E. SYSTEM

COLOR	ITEM #
Multi-Cam	90-0063
Coyote	90-0077

SPECIFICATIONS

H 15 in. x W 10 in. x D 7.5 in.
Weight: 10 lb 14 oz

Water Rescue BELT THROW BAG

- Multipurpose water rescue solution designed for easy deployment and retrieval for rescuers
- Quick release reverse tapered throw bag for rapid packing of Triton rope and a slim profile for users
- 2 in. adjustable belt with quick release strap, 1.5 in. plastic cam buckle, and red toggle ball for easy ID



BELT THROW

ITEM #
90-0065

SPECIFICATIONS

Rope Bag: L 11 in. x W 4.5 in. x D 4 in.
Gear Pouch: L 5.5 in. x W 5 in. x D 2 in.
Total Weight: 2 lb 2 oz

50 Foot THROW BAG

- Must-have necessity for water rescue
- Durable, easy-to-use rescue device with highly visible orange nylon bag and mesh panels that drain quickly for use in multiple rescue operations
- Packed with 50 ft of 7.5 mm Triton Pegasus Rope® rated at 3,050 lbf
- Triton rope is strong enough for multiple uses: as a throw line, tension diagonal, tow line and for pinned kayak recovery



50 FT THROW BAG

ITEM #
90-0064

SPECIFICATIONS

Rope Bag: L 11 in. x W 4.5 in. x D 4 in.
Gear Pouch: L 5.5 in. x W 5 in. x D 2 in.
Total Weight: 2 lb 2 oz



SPECIFICATIONS

L 11 in. x W 7 in. x D 7 in.
Weight: 5 lb 13 oz

HELO TAGLINE


ITEM #

90-0069

HELO TAGLINE KIT

- **Help control litter movement during hoist operations** - a requirement on all Army MEDEVAC helicopters
- Contains 250 ft of 7.5 mm Triton Pegasus Rope® featuring a highly visible polypropylene sheath with an Innegra™ core

HELO TAGLINE KIT CONTENTS:

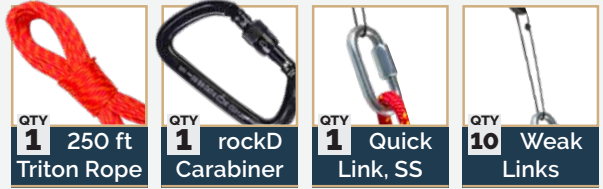


A-VRR *Adaptable Vertical Response Rack Kit*

The Adaptable Vertical Response Rack was designed meet and exceed full-spectrum operational vertical requirements for urban through mountain / canyon access, rescue, and evacuation variables. The kit capabilities were driven utilizing NFPA 1006 (Operations & Technician Level) JPR's and the KPP's (Key Performance Parameters) within the USSOCOM SOF CASEVAC Program framework. The A-VRR Kit's organic versatility allows it to be scaled up to meet all relevant job performance requirements (JPR's) for complex team access and/or rescue or be scaled down for a "fast & light" mission specific assault configuration.

P-A-C-E methodology was engineered throughout the kit design, enabling multiple options for progressive capture and rope grab capabilities depending on the environmental pathology.

- Haul / Mechanical Advantage Capability – Multiple simple, compound and complex configurations ranging from a 2:1 through 9:1
- Full spectrum Friction Configurations for Rappels, Bail-outs, And One or Two Person Lowering Systems
- Traditional and MA-assisted Ascending Techniques
- Multiple Basic and Anchor Configurations to Include Contingency and Releasable Systems
- Load Transfer Capabilities to include Knot Passing on Rappel and within a system and casualty Pick-Offs

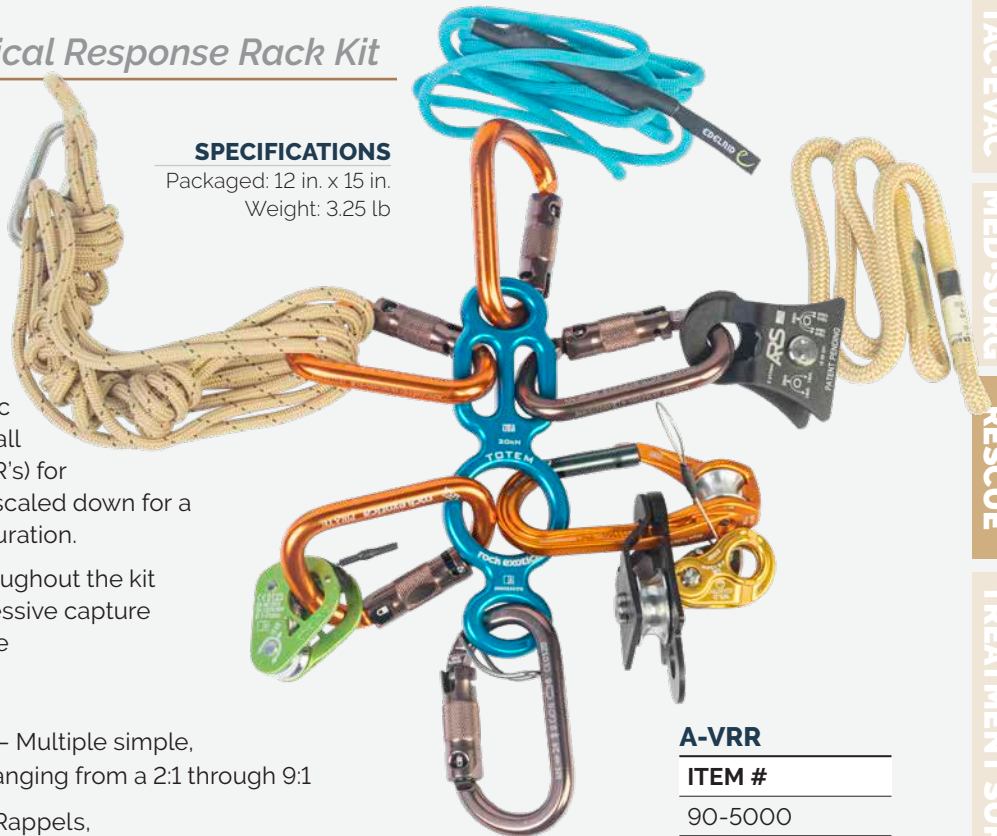

SPECIFICATIONS

Packaged: 12 in. x 15 in.
Weight: 3.25 lb

A-VRR

ITEM #

90-5000

A-VRR KIT CONTENTS:





R.A.C.E. SUB-T *Rapid Access Casualty Extraction Subterranean Kit*

Specifically designed for full-spectrum vertical operations within harsh environments

- Ultra-Lightweight Components
- Hydrophobic Software
- Innate requisite variety (redundancy)
- Easy sustainment
- HRO (High Reliability Organizing) Engineered
- Scalable based on mission context, organic assets, & end-user
- Resilient to environmental pathology
- Capability focused - up to Technician Level JPR's (Job Performance Requirements)
- Meets and exceeds operational requirements of SOF CASEVAC Program



**Subterranean
Deep Hydro | Cave
Mountain | Rough Terrain Access | Evac
Class B & C Canyon | Slot
Mountain | Jungle**

R.A.C.E. SUB-T CONTENTS:

 QTY 1 R.A.C.E. Nylon Bag - MTC	 QTY 1 Rope Edge Guard	 QTY 1 100ft 9.5mm Olympus Rope	 QTY 1 Totem Decent Control Device	 QTY 1 Oval Auto-lock Carabiner	 QTY 1 Oval Auto-lock Carabiner w/Wire Eye	 QTY 3 HMS Auto-lock Carabiner
 QTY 1 MagnaPulley Small Pulley Pair with 7/16 in. w & w/o Beckett	 QTY 1 Pulley with Rope Grab	 QTY 1 H-Frame Auto-lock Carabiner w/Pulley	 QTY 1 Ascender / Rope Clamp	 QTY 1 16ft 6mm XTEC Cord	 QTY 1 8mm VT Prusik	 QTY 1 120cm Sling w/Aramid Core
	 QTY 1 22ft 1in. Nylon Tubular Webbing					 QTY 1 7mm Screw Link

SPECIFICATIONS

H 15 in. x W 10 in. x D 7.5 in.
Weight: 10 lb 8 oz

R.A.C.E. SUB-T

ITEM

90-5002

FLIGHT RISK KIT

- Ultra-lightweight, compact rescue, and emergency egress / bail-out modular system that easily fits within pocket, bag, or harness
- Auto-Locking F3 Descent device engineered with a double-action lever to give a brake assist in the event of a panic grab
- Gives the end-user a rapid, smooth, and safe egress capability
- Counter-balance casualty evacuation capability

FLIGHT RISK KIT

ITEM

90-5003

FLIGHT RISK KIT CONTENTS:

 QTY 1 IPOK Nylon Holder-BLK	 QTY 1 Oval Auto-lock Carabiner w/Wire Eye	 QTY 1 HMS Auto-lock Carabiner	 QTY 1 120cm Sling w/Aramid Core	 QTY 1 F3 Anti-Panic Descent Device	 QTY 1 75ft 6.8mm RIT 900 Cord
--	---	---	---	--	---



SPECIFICATIONS

H 6.5 in. x W 6 in. x D 3.5 in.
Weight: 3.75 lb



For Multiple Operational Mission Sets:

- **Austere** | Remote Vehicle Recovery
- **Heavy Lifting in Confined Space**
- **Subterranean**
- **Tunnel** | Mine Environments

SMAS-H KIT

Scalable Mechanical Advantage System-HEAVY

- Scalable Mechanical Advantage System designed for specialized heavy rigging situations constrained by limited amount of personnel
- Simple, rapid, and adaptable rigging configurations from a 3:1 up to a 23:1
- Developed and tested within Mobility, Vehicle Recovery, Subterranean, Confined Space, CBRN, Recovery, and Large Animal movement, stabilization, and rescue environments

SPECIFICATIONS
 H 7.75 in. x W 13 in. x D 9 in.
 Weight: 5.7 lb

SMAS-H KIT	
ITEM #	
90-5004	

SMAS-H CONTENTS:

 QTY 1 Squad Nylon Bag - ODG	 QTY 1 75ft 9.5mm Olympus Rope	 QTY 2 Oval Auto-lock Carabiner	 QTY 3 HMS Auto-lock Carabiner	 QTY 1 MagnaPulley Small Pulley Pair 7/16 in. w & w/o Beckett
 QTY 1 Pulley with Rope Grab	 QTY 1 H-Frame Auto-lock Carabiner w/ Pulley	 QTY 2 Ascender / Rope Clamp	 QTY 1 8mm VT Prusik	 QTY 1 120cm Sling w/ Aramid Core
 QTY 1 JAG Traxion Double Pulley	 QTY 1 50ft 6mm XTEC Cord	 QTY 1 5mm Screw Link	 QTY 1 22ft 1in. Nylon Tubular Webbing	 QTY 2 Lg. Rigging Rubber Band



iMAS KIT

Individual Mechanical Advantage System

- Ultra-compact, lightweight - Grab & Go "Piggy-Back" Mechanical Advantage System engineered for full spectrum operational requirements
- Simple to use yet expandable to the level of expertise of the rescuer and the complexity of the rescue
- Easily carried or packed/attached onto an aid bag, mountain pack, cargo pocket, harness gear loop, or fast-rope bag
- Traditional hauling operations (3:1, 4:1, 5:1, and 9:1)

SPECIFICATIONS
 H 7 in. x W 4.5 in. x D 2.75 in.
 Weight: 4.25 lb

iMAS KIT	
ITEM #	
90-5005	

iMAS CONTENTS:

 QTY 1 TORK Nylon Bag - BLK	 QTY 2 Oval Auto-lock Carabiner	 QTY 1 HMS Auto-lock Carabiner	 QTY 1 MagnaPulley Small Pulley Pair 7/16 in. w & w/o Beckett
 QTY 1 Ascender / Rope Clamp	 QTY 1 8mm VT Prusik	 QTY 1 120cm Sling w/ Aramid Core	 QTY 1 JAG Traxion Double Pulley
 QTY 1 30ft 6mm XTEC Cord	 QTY 1 5mm Screw Link	 QTY 1 Lg. Rigging Rubber Band	



HIGH THREAT EXTRACTION KIT



SPECIFICATIONS

L 8.5 in. x W 8.75 in. x D 4.5 in.
Weight: 3 lb 6 oz



- Facilitates immediate casualty extraction so that Care Under Fire casualty treatment may be performed as soon as possible
- Self-contained kit with all of the essential tools for tactical operators and rescue personnel
- Includes a C-A-T® tourniquet to immediately address extremity hemorrhage at the point of wounding
- Rugged nylon bag attaches to the duty belt as a drop down leg rig or to the assault vest/body armor with PALS/MOLLE-style connectors

HIGH THREAT EXTRACTION KIT

COLOR	ITEM #
Black	80-0139
OD Green	80-0140
Coyote	80-0141
Digital Camo	80-0142

HIGH THREAT EXTRACTION KIT:



QTY 1 C-A-T Tourniquet

QTY 1 RAT Strap

QTY 1 Dragon Handle System

QTY 1 Hasty Harness

QTY 1 Carabiner Non-locking

QTY 1 Auto-lock Carabiner

QTY 1 Pulley Carabiner

QTY 2 Red ChemLight



RESCUE | Rescue Gear

IRUPU HARNESS

- Multi-disciplinary rugged harness for use in operational environments (*Mountain evac, Sub-T, Canyons, & Caves*)
- 2 robust, floating gear loops
- Enhanced stitching and replaceable seat provide uniform protection in environments with abrasion hazards
- Stitched clip-in buckles with tension straps for easy on & off
- Left/ right leg loop markings make it easier to put the harness on correctly
- Large, textile tie-in point



SPECS

UIAA: Yes

Certification: EN 12277 Type C

Waist: 65-110 cm | Leg Loops Max: 70 cm

Weight: 21.16 oz (600 gm)

IRUPU HARNESS

ITEM #
90-4020



MADILLO HELMET

The first & only foldable climbing helmet on the market

EDELRID

SPECS

CE marking: CE 0123

Weight: 13.76 oz (390 gm)

MADILLO HELMET

ITEM #
90-4019

- Robust and impact-resistant injection-moulded ABS shell
- Innovative combination of 3 different foam types (EPS, EPP and EVA) deliver optimum safety and comfort
- Fully adjustable chinstrap for optimal fit
- 4 robust headlamp clips
- Detachable, washable padding
- Side/front/back impact tested to EN 12492



RAT STRAP

Rescue Assault Tether

- For head/feet first linear drags during shield-assisted, open-field/indoor-corridor rescues or single-rescuer extractions
- Innovative Daisy Chain configuration offers multiple reinforced attachment points



SPECIFICATIONS

Pkg: H 1 in. x W 4 in. x D 5.5 in.
 Extended: L 62 in.
 Minimum Break Strength: 2,702 lb
 Weight: 12.7 oz

RAT STRAP

ITEM #	NSN #
90-0002	6530-01-632-7872

DRAGON HARNESS

Under the Vest Harness

- Built-in handle system for immediate high threat extraction
- Operator harness to aid in the rescue of downed personnel without the need to holster or sling weapons
- Prevents the "ride up" of the casualty's ballistic vest, reducing exposed, unprotected body regions



SPECIFICATIONS

Height: 20 in. from waist to non-deployed handle
 Length: 14 in. from harness attachment point to top of handle
 Weight: 1 lb 1.6 oz.
 Waist Size: 30 in. to 45 in.

DRAGON HARNESS

ITEM #
90-0003

EX-TRACT STRAP

Point of Wounding Extraction

- Point of wounding extraction
- Back pack carry
- Hasty harness for traditional and confined space lowering and hauling



SPECIFICATIONS

Pkg: H 18 in. x W 3.5 in. x D 1 in.
 Extended: L 58 in.
 Weight: 9.3 oz

EX-TRACT STRAP

ITEM #
90-0082

DRAGON HANDLE

Tool offering head/feet first drags

- Maintain an anatomically correct posture for rapidly extracting greater weight
- Decreases manpower needs in "zeroed in" locations
- Black Omega, non-locking "D" carabiner rated for 31kN (6,968.8 lb).



SPECIFICATIONS

H 6.3 in. x W 4.1 in. x D 2.25 in.
 Weight: 11.3 oz
 Device Length: L 42 in.
 US Patent(s): Pending

DRAGON HANDLE

ITEM #	NSN #
90-0004	6530-01-598-9595

HASTY HARNESS *Multi-purpose Harness*

- Indispensable gear for every webbing equipment toolbox
- Resourceful, multi-purpose harness for unlimited applications & maximum versatility
- Constructed of 22 ft loop of 1 in. Tech Tape (MBS: 4,400 lb) with opposing color thread to aid in visual safety inspection

SPECIFICATIONS

Loop: L 10 ft 11 in. x W 1 in.
 Circumference - 21 ft 10 in.
 Min. Break Str: Tech Tape (4,400 lb)
 Sewn Sling (6,000 lb)



HASTY HARNESS NATO NSN*

COLOR	ITEM #	NSN #
Red	90-0044	----
Black	90-0001	6530-01-674-9774*

INTRO
 M-A-R-C-H
 KITS
 TAC-EVAC
 MED-SURG
 RESCUE
 TREATMENT SUPPORT
 TRAINING
 INDEX



INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

TIBLOC

Ascender/Rope Grab by Petzl®

Ultra-light Compact Ascender/Progress Capture Device for rope ascents, hauling systems or friction knot replacement in self-rescue situations

SPECS

Rope Size: 8-11 mm
Weight: 1.23 oz (35 g)

ITEM

ZZ-0547



ROPE GRAB PULLEY-SPOC by Edelrid®

Ultra-light pulley designed for lifting loads, crevasse rescue operations, or as an emergency ascender

SPECS

Rope Size: 7-11 mm
Weight: 2.12 oz (60 g)
Diameter Pulley: 20 mm

ITEM

90-4005



ROLLCLIP PULLEY Carabiner

Pulley-carabiner with a non-locking gate opening

SPECS

Rope: 7 to 13 mm
Major axis: 20 kN
Minor axis: 8 kN
Open gate: 7 kN
5.25 in. x 3 in. x 1.25 in.
Weight: 3.76 oz

ITEM

90-4006



PIRATE CARABINER

by Rock Exotica®

Ideal as a belay biner at the harness, or as a connector in a system

SPECS

Gate Opening: 25 mm
MBS Major: 26 kN
MBS Minor: 12 kN
MBS Open: 7 kN
H 107 mm x W 71 mm
Weight: 3.09 oz (87 g)

ITEM

90-4003



ROCKO AUTO LOCK

Carabiner by Rock Exotica®

3-stage Autolocking carabiner with intuitive pull-down and twist action

SPECS

Gate Opening: 24 mm
MBS Major: 24 kN
MBS Minor: 11 kN
MBS Open: 6 kN
H 107 mm x W 71 mm
Weight: 2.88 oz (81.65 g)

ITEM

90-4001



ROCKO WIRE EYE

Carabiner by Rock Exotica®

Standard carabiner with the safety and capability of a lanyard pin, or captive eye

SPECS

Gate Opening: 22 mm
MBS Major: 24 kN
MBS Minor: 9 kN
MBS Open: 6 kN
H 107 mm x W 71 mm
Weight: 2.96 oz (84 g)

ITEM

90-4002



ROCKO NON LOCKING

Carabiner by Rock Exotica®

Non locking enhanced Oval biner with Snag-resistant keylock nose design

SPECS

Gate Opening: 25 mm
MBS Major: 22 kN
MBS Minor: 10 kN
MBS Open: 6 kN
H 107 mm x W 71 mm
Weight: 2.36 oz (67 g)

ITEM

90-4022



OMEGALITE 4 WIREGATE

ISO Cold Forged design

SPECS

Gate Opening: 110 in. (28 mm)
MBS Major: 26 kN
MBS Minor: 9 kN
MBS Open: 8 kN
L 4.25 in. x W 2.75 in.
Weight: 0.12 lb

ITEM

ZZ-0314



QUICK LINK SCREW LOCK

5mm Oval Stainless Steel Quick Link Carabiner

SPECS

Rod stock
Diameter: 5 mm
Gate opening: 6 mm
Major Axis: 25 kN
H 1.97 in. x W 0.98 in.
Weight: 0.7 oz

ITEM

90-4015



QUICK LINK SCREW LOCK

7mm Oval Stainless Steel Quick Link Carabiner

SPECS

Rod stock
Diameter: 7 mm
Gate opening: 16 mm
Major Axis: 25 kN
Minor Axis: 10 kN
H 3.25 in. x W 1.25 in.
Weight: 2.4 oz

ITEM

ZZ-0544



ALPHA FIX Y SLING

30/60 cm by Edelrid®

Y-Shaped self belay sling for edge restraint, DCD extender, knot passing on rappel, edge transition safety, & hasty pick-offs

SPECS

Load capacity: 23 Kn
CE 0123
Weight: 4.16 oz (118 gm)

ITEM

90-4021



*NOTE:
Not for use as a via ferrata set

ARAMID CORE BLUE SLING

120 cm by Edelrid®
High-strength sewn Aramid sling ideal as a full-spectrum anchor, harness extension, alternate prusik, & canyon quickdraw for Rig to Raise conversions

SPECS

Load capacity: 22 Kn
CE 0123
Width: 0.24 in. (6 mm)
Weight: 0.78 oz (22 gm)

ITEM

90-4011





BASIC ASCENDER

by Petzl®

Compact ascender with excellent grip for rope ascent on a fixed line or as progress capture in a hauling system



SPECS

Rope Size: 8-11 mm
Weight: 3 oz (85 g)

ITEM

ZZ-0554

GRIGRI 2 BELAY DEVICE

by Petzl®

Compact, ultra-light belay device allowing excellent control during the descent with assisted braking capability



SPECS

Rope Size: 8-11 mm
Weight: 6 oz (170 g)

ITEM

ZZ-0545

MINI MACHINED PULLEY

by Rock Exotica®

Mini pulley with recessed axle screw, machined sideplates & integrated axle for a stronger, compact, snag-free profile



SPECS

Rope Size: up to 11 mm
Weight: 2.9 oz (83 g)

ITEM

ZZ-0548

REVERSO 4

Belay/Rappel device by Petzl®

Regulates the amount of braking friction on the rope according to diameter and condition



SPECS

Weight: 2.08 oz (59 g)

ITEM

ZZ-0546

TRIRIG RIGGING PLATE

by Rock Exotica®

Lightweight, high performance 3 hole rigplate with numbered holes to assist with detailed rigging



SPECS

Break Strg: 36 kN
Weight: 1.8 oz (51 g)

ITEM

ZZ-0553

ROCKD CARABINER

by Rock Exotica®

Narrow asymmetric "D" shape machined frame with snag resistant keylock nose and large, ergonomic angled gate opening



SPECS

Gate Opening: .24 mm
Weight: 2.6 oz (73 g)

ITEM

ZZ-0543

FOOTAPE

Adjustable Foot Loop by Petzl®

Adjustable nylon foot loop for attachment to the BASIC rope clamp for rope ascents



SPECS

Weight: 2.3 oz (65 g)

ITEM

ZZ-0551

ANNEAU SLINGS

by Petzl®

Polyester sling to set up a belay or an anchor



SPECS

W 19 mm
L 60 or 120 cm
Wgt: 60 cm - 2 oz
120 cm - 3.5 oz

LENGTH	ITEM #
60 cm	ZZ-0549
120 cm	ZZ-0550

Rescue Knives

EDELRID



ROPE TOOTH KNIFE

by Edelrid®

- Ideal for all cutting requirements outdoors
- Perfectly suitable to shorten static lines or ropes
- Large swivel hole for attaching to harness or backpack with carabiner
- Sharp, locking blade (partly serrated)
- High-quality Stainless Steel

SPECIFICATIONS

L 3 in. x W 1 in. (closed blade)
Blade Length: 3 in.
Weight: 1.69 oz (48 gm)

ROPE TOOTH KNIFE

ITEM

90-4027



INTRO

M.A.R.C.H

KITS

TAC.EVAC

SURG

RESCUE

TREATMENT SUPPORT

TRAIL



BENCHMADE® KNIVES & CUTTERS

8.50 in.
8.00 in.
7.50 in.
7.00 in.
6.50 in.
6.00 in.
5.50 in.
5.00 in.
4.50 in.
4.00 in.
3.50 in.
3.00 in.
2.50 in.
2.00 in.
1.50 in.
1.00 in.
0.05 in.



1 HOOK SAFETY CUTTERS by Benchmade®

- Hook cutters designed to quickly and reliably cut through a wide array of materials
- Extensively tested and widely used by the US Military for medical and egress purposes

ITEM #	STYLE	NSN #
ZZ-0652	7 Hook	----
ZZ-0250	8 Hook	----

SPECS
 7 Hook: L 4.3 in.
 8 Hook: 6.5 in.
 Weight: 1.5 oz

2 SOCP RESCUE TOOL by Benchmade®

- Fixed 440C high-chromium stainless steel blade hook for rescue performance and less-than-lethal self-defense
- Includes carbide glass breaker, oxygen wrench, single-piece handle with finger loop & molded sheath
- Designed by Special Operations Combatives Program (SOCP)

SPECS
 Blade: 2.5 in. Overall: 6.75 in.
 Weight: 1.88 oz

ITEM #	NSN #
ZZ-0760	----



3 NAVY SAR SOML 112 SBK by Benchmade®

Approved for Navy SAR - "So Others May Live"

- N680 ultra stainless steel blade with opposing fixed bevel blade style with blunt tip, integrated rescue hook and carbide glass breaker
- Boltaron® sheath designed to fit current Navy SAR kits

SPECS
 Blade: 3.5 in. Overall: 8.02 in.
 Weight: 3.2 oz

ITEM #	NSN #
90-0078	1095-01-677-1013

4 TRIAGE 917 BK by Benchmade®

- Black drop point style CPM-S30V steel blade with a smooth edge
- 440C folding hook & carbide glass breaker

ITEM #	NSN #	SPECS
90-0092	----	Blade: 3.5 in Overall: 8.3 in. Weight: 5.24 oz

5 TRIAGE 916 SBK by Benchmade®

- Dynamic triple utility tool featuring three essential rescue tools in a single compact design
- N680 high corrosion resistant, partially serrated, opposing bevel blade, 440C folding hook & carbide glass break

SPECS
 Blade: 3.5 in Overall: 8.2 in.
 Weight: 5.1 oz

ITEM #	HANDLE	NSN #
90-0096	Orange	----
90-0060	Black	----



Leatherman® RAPTOR™ MULTI-TOOL

High-performance multi-tool equipped with a variety of mission critical tools essential to emergency medical practitioners in the field

1. 420HC Stainless Steel Folding Medical Shears
2. Strap Cutter
3. Oxygen Tank Wrench
4. Ring Cutter
5. Carbide Glass Breaker
6. 5 cm Ruler
7. Holster with Clip

RAPTOR MULTITOOL

HANDLE	ITEM #	NSN #
Black	70-0014	5110-01-626-8766
Blue	70-0264	----
Coyote	70-0265	----
Red	70-0266	----

SPECS

Shear Length: 1.9 in.
Closed Length: 5 in.
Weight: 5.8 oz

HANDLE COLORS:



See our Treatment Support Section on page 87 for the compatible Scissor Leash

SUPERTOOL 300

by Leatherman®

Heavy-duty, all purpose multi-tool featuring 19 devices in one easy to carry configuration

SPECIFICATIONS

Closed: L 4.5 in.
Primary Blade: L 3.2 in..
Weight: 9.6 oz

SUPERTOOL 300

ITEM #
ZZ-0542

1. Needlenose Pliers
2. Regular Pliers
3. Replaceable Wire Cutters
4. Hard-Wire Replaceable Cutters
5. Stranded-Wire Cutters
6. Electrical Crimper
7. Wire Stripper
8. 420HC Knife
9. 420HC Serrated Knife



10. Saw
11. Awl with Thread loop
12. 9 in. Ruler
13. Can Opener
14. Bottle Opener
15. Wood/Metal File
16. Phillips Screwdriver
17. Large Screwdriver
18. Medium Screwdriver
19. Small Screwdriver



MUT® by Leatherman®

Tactical and practical multi-tool for military, law enforcement and civilian shooters - features 16 devices

1. Needlenose Pliers
2. Regular Pliers
3. Replaceable Wire Cutters
4. Hard-Wire Replac. Cutters
5. Stranded-Wire Cutters
6. Electrical Crimper
7. 420HC Knife
8. Saw
9. Replaceable Cutting Hook
10. Hammer
11. Bolt Override Tool
12. Bronze Carbon Scraper
13. #8-32 Cleaning Rod/Brush Adapter
14. Firearm Disassembly Punch
15. Carabiner/Bottle Opener Punch
16. Large Bit Drive

SPECIFICATIONS

Closed: L 5 in.
Primary Blade: L 3 in.
Weight: 11.2 oz

MUT MULTITOOL

ITEM #
ZZ-0539



TACTIKKA HEADLAMP

ITEM #

ZZ-0552

SPECIFICATIONS 

Light Output: Min. 5 Lumens
Max: 120 Lumens
(160 Lumens on Boost mode)
Weight: (with batteries) 3 oz

TACTIKKA + RGB

Headlamp by Petzl®

- Ambient Mode: Wide beam with reduced brightness, avoids blinding others in a group situation and provides a long burn time
- Close Range Vision Mode - Mixed beam with wide, uniform flood or 3 colored LED lighting modes (red, green & blue) to remain discreet
- Boost Mode for temporary access to maximum power, 160 lumens (60 m)
- Strobe mode for signaling one's presence; accessible at any time
- Reserve Mode- Burn time management automatically switches on by signaling with blinking light when batteries are almost drained. Red lighting activates to preserve a minimal amount of usage

NAR HEADLAMP

- 100 Lumen LED Light
- Dual color (white and red)
- High and low lighting modes (White)
- Sturdy and simple to rotate single arm bracket
- Large push button for operation
- Splash and water resistant (IPX4)
- Comes with 2 AAA batteries
- 58 hour burn time



NAR HEADLAMP

ITEM #	NSN #
ZZ-0106	6515-01-542-7937

SPECIFICATIONS

Size: W 2.25 in. x D 1.5 in.
Weight: 2.25 oz

Streamlight HEADLAMP

- Ring Sealed construction – Waterproof to 1M for 1 Hr
- Up to 6 hrs, runtime on high up to 24 hrs. on low
- Impact resistant ABS housing
- Includes (2) 1.5V "AAA" batteries
- Solid State Power Regulation provides maximum light output throughout battery life
- Shatter resistant polycarbonate lens



SPECIFICATIONS

Faceshield Diameter: 1.3 in.
Maximum headlight: W 2.25 in. x D 1.6 in.
Weight: (with batteries & headband) 2.75 oz

STREAMLIGHT HEADLAMP

ITEM #

ZZ-1090

GECKO GRIP®

Multi-purpose Tape

Laminated polyethylene cloth tape featuring an aggressive natural rubber based adhesive formulation and perfect for heavy-duty applications

SPECIFICATIONS

Pack of 6:
L 12.25 in. x W 0.9 in. x D 0.7 in.
Individual Roll: W 2 in. x L 5 yd
Weight: Roll: 2.7 oz
Pack: 1 lb 8.7 oz



GECKO GRIP

UNIT	ITEM #
6 Pack	80-0131



TOMAHAWK NIGHT VISION

by First Light®

- Right angle tactical light with rotating head and 3 button control panel for one-handed thumb actuation
- Provides White, IR, and Red+Green LED mix for highlighting blood during casualty treatment
- Featuring Integrated high intensity IR signaling beacon
- Includes momentary and three intensity levels of constant-on (all output colors)
- Can be used handheld or hands-free with MOLLE TRS Blade Mount



TOMAHAWK NIGHT VISION

CASE COLOR	TASK LIGHT	ITEM #
Coyote	UV Task	70-0015
OD Green	UV Task	70-0112
Coyote	Blue Task	70-0111
OD Green	Blue Task	70-0110

SPECIFICATIONS

Output Wavelengths: White, IR, & Red+Green (Blood Highlighting), IR Signaling Beacon
 Maximum White Light Output: 425 lumens
 Max Output Runtime: 2 hours (ANSI Tested)
 Power Source: CR123 Batteries (2) | Material: Aerospace-Grade Aluminum
 With mount: H 3.4 in. x L 2.4 in. x W 1.5 in. | Weight: (with mount & batteries) 8 oz

TOMAHAWK MULTI-COLOR

by First Light®

- Right angle tactical light with rotating head and 3 button control panel
- Provides White, Red, Green, and Red+Green LED mix for highlighting blood during casualty treatment
- Includes momentary and three intensity levels of constant-on

SPECIFICATIONS

Output Wavelengths: White, Red, Green, & Red+Green (Blood Highlighting)
 Maximum White Light: 425 lumens
 Max Output Runtime: 2 hours (ANSI Tested)
 Power Source: CR123 Batteries (2)
 Material: Aerospace-Grade Aluminum
 With mount: H 3.4 in. x L 2.4 in. x W 1.5 in.
 Weight (with mount & batteries): 7.2 oz

TOMAHAWK MULTI-COLOR

CASE COLOR	TASK LIGHT	ITEM #
Coyote	None	70-0016
OD Green	None	70-0117
Coyote	UV Task	70-0116
OD Green	UV Task	70-0115
Coyote	Blue Task	70-0114
OD Green	Blue Task	70-0113



TORQ

by First Light®

- Right angle tactical light, rotating head and 3 button control
- White, Red, Green, and Red+Green LED mix for highlighting blood
- Includes momentary and three intensity levels of constant-on (all output colors)
- Features integrated high intensity white tactical strobe and red/white/green flashing safety beacon
- 10 hour Runtime on high output

TORQ

ITEM #	NSN #
70-0118	6515-01-646-2636

SPECIFICATIONS

Wavelengths: White, Red, Green, & Red+Green (Blood Highlighting)
 Max Runtime: 10 hours (ANSI Tested)
 With mount: H 3.9 in. x L 2.6 in. x W 1.3 in.
 Weight (with mount & batteries): 7.1 oz







TREATMENT

see pages 88-96

SUPPORT • PPE • TRIAGE • ACCESSORIES





INTRO
M-A-R-C-H
KITS
TAC-EVAC
MED-SURG
RESCUE
TREATMENT SUPPORT
TRAINING
INDEX

POLAR SKIN™

- The latest cooling technology for effective Human Thermo-regulation
- Polar Skin Mini-Mesh,™ a unique smart textile with a hydrofunctional polymer that allows for liquid absorptive capacity 350% greater than base weight
- ArcticAqua,™ a proprietary solution that prevents rigidity at temperatures well below 0°F while providing the requisite moisture to support evaporative, convective and conductive cooling
- Adheres to current U.S. Army and National Athletic Trainers' Association protocols for heat illness treatment
- Protracted direct cooling of high blood flow regions

COMPREHENSIVE COOLING SYSTEM

- Highly-effective portable cooling system for the treatment of heat-related illnesses
- Scientifically proven to reduce core temperature by .09°F/minute in 100°F heat when used according to protocol



COMPREHENSIVE COOLING

ITEM #

ZZ-0674

SPECIFICATIONS

L 11 in. x W 10 in. x H 4 in.

Weight: 8 lb 5.3 oz

COMPLETE CARE PACK

- Portable, highly insulated option
- Contains one Polar Skin Comprehensive Cooling System and one heavy-duty, leak & sweat proof cooler that keeps ice from melting for 24 hours in +100°F heat



COMPLETE CARE PACK

ITEM #

80-0451

SPECIFICATIONS

W 14 in. x L 12 in. x D 7 in.

Weight: 10 lb 3.2 oz

COMPLETE CARE BACKPACK

Available in either Single Patient (one comprehensive cooling system) or Double Patient (two cooling systems) configurations

COMPLETE CARE BACKPACK

SIZE	ITEM #
Single	80-0543
Double	80-0544

SPECIFICATIONS

W 19.5 in. x L 15 in. x D 7.5 in.

Single Weight: 12 lb 13 oz

Double Weight: 21 lb 13 oz



Available in Standard, Large and Backpack COOLER

Standard size accommodates one Polar Skin™ Comprehensive Cooling System and Large or Backpack Size accommodates two cooling systems



POLAR SKIN COOLERS

COLOR	SIZE	ITEM #
Red	Standard	ZZ-0682
Black	Standard	ZZ-0759
Black	Large	ZZ-0756
Black	Backpack	ZZ-0794

SPECIFICATIONS

Standard: L 14 in. x W 12 in. x H 7 in.

Weight: 1 lb 3.2 oz

Large: L 18.5 in. x W 13 in. x H 9.5 in

Weight: 2 lb 13 oz

Backpack: L 19.5 in. x W 15 in. x D 7.5 in.

Weight: 3 lb 13 oz



ICE SHEETS *Standard*

- Principle cooling component of Polar Skin™ cooling systems
- Adheres to current National Athletic Trainers' Assoc and U.S. Army protocols for heat illness treatment
- Sheet Size: W 36 in. x L 30 in.



STANDARD ICE SHEETS

UNIT	ITEM #	NSN #
Single Pack	ZZ-0668	6532-01-653-6054
Double Pack	ZZ-0669	----
Triple Pack	ZZ-0755	----

SPECIFICATIONS

Single: L 4.5 in. x W 7.75 in. x H 0.5 in. Wgt: 9.6 oz
 Double: L 9.5 in. x W 8 in. x H 0.5 in. Wgt: 1 lb 6.4 oz
 Triple: L 12.5 in. x W 7.5 in. x H 0.5 in. Wgt: 2 lb 3.2 oz

ICE SHEETS *Large*

Same characteristics as the standard size Ice Sheets, but in a larger 36 in. x 60 in. size

LARGE ICE SHEETS

ITEM #
ZZ-0758

SPECIFICATIONS

L 8.5 in. x W 5.5 in. x H 1.5 in.
 Weight: 1 lb 8 oz

C2E ICE PACK BUNDLE

- Four Comfort Gel Packs
- Facilitates Conductive, Convective, and Evaporative (C2E) cooling of high blood flow (large artery) regions of the body



C2E ICE PACK

ITEM #
ZZ-0675

SPECIFICATIONS

Indiv. Gel Pk: W 6 in. x L 10 in.
 Pkg: L 8 in. x W 9.5 in. x H 1.25 in.
 Weight: 4 lb 4 oz

COOLING COLLAR



- Two gel packs encased in integrated pouches
- Designed to align with the carotid artery

COOLING COLLAR

ITEM #
ZZ-0670

SPECIFICATIONS

Indiv. Gel Pk: W 6 in. x L 10 in.
 Pkg: L 8 in. x W 9.5 in. x H 1.3 in.
 Weight: 4 lb 4 oz



Splinting



TC SPLINT

Telescoping Collapsible Splint

- Multi-capable device for splinting a wide variety of major and minor fractures
- Lightweight and compact for easy field carry with 6 easily adjustable interconnected telescoping poles and 4 wide elastic bands with hook & loop securing
- Includes clearly marked directions on the device as well as an easy-to-use, step-by-step instructional insert

SPECIFICATIONS

Pkgd: L 8 in. x W 5 in. x D 3 in.
 Open: W 9.5 in. x L 34 in.
 Weight: 1 lb 4 oz

TC SPLINT

ITEM #
50-0039



online video available scan to view >



INTRO

M.A.R.C.H

KITS

TAC·EVAC

MED·SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX



TREATMENT SUPPORT | Splinting



TACTICAL TRACTION SPLINT (TTS) Compact Mechanical Traction Device

SPECIFICATIONS

Packaged:
L 11 in. x W 7 in. x D 3 in.
Weight: 1 lb 5 oz

CE


- Relieve pain, increase patient comfort and reduce the risk of vascular injury or nerve damage in femur fractures

- Lightweight and portable kit with clearly marked elastic straps for the proper fit to the thigh, knee & shin

TACTICAL TRACTION SPLINT

COLOR	ITEM #	NSN #
Black	50-0015	6515-01-346-9186
Coyote	50-0033	----
Digital Camo	50-0034	----
OD Green	50-0032	----

TACTICAL TRACTION SPLINT (TTS):

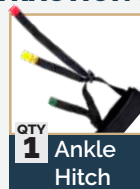

SAM® SPLINT II

- Versatile, malleable splint easily molds to fit most fractures
- Lightweight, comfortable & strong

- Waterproof, Reusable

SPECIFICATIONS

H 9.25 in. x W 4.25 in. x D 1 in.
Unfolded: L 36 in. x W 4.25 in.
Weight: 4.5 oz

SAM SPLINT

ITEM #	NSN #
50-1005	6515-01-494-1951

SAM® PELVIC SLING *for Pelvic Ring Fractures*

- Force-controlled, ideal compression every time without the guesswork
- Rapid, simple, three-step application, usually in less than one minute



SAM PELVIC SLING

ITEM #	NSN #
50-0004	6515-01-509-6866

SPECIFICATIONS

H 7.25 in. x W 7.75 in. x D 4 in.
Hip Circumference: 32 in. to 50 in.
Weight: 9.8 oz

SAM® SOFT SHELL SPLINT

SPECIFICATIONS

L 15 in. x W 4.25 in. x D 0.5 in.
Weight: 3.2 oz

SAM SOFT SHELL SPLINT

UNIT	ITEM #	NSN #
Each	50-0008	6515-01-523-5440
Box of 20	50-0009	6515-01-507-8074



TRAUMA SHEARS

SPECIFICATIONS

SMALL: L 6.25 in. x W 2.75 in. x D 0.3 in.
Weight: 15 oz
LARGE: L 7.25 in. x W 3.7 in. x D 0.4 in.
Weight: 2 oz

TRAUMA SHEARS * NATO NSN

SIZE	ITEM #	NSN #
S 6.25 in.	ZZ-0064	6510-01-622-1624*
L 7.25 in.	ZZ-0063	6515-01-538-9276



SCISSOR LEASH™

SPECIFICATIONS

H 6 in. x W 16 in. x D 0.6 in.
Cable Length - 36 in.
Weight: 16 oz

SCISSOR LEASH

ITEM #	NSN #
70-0005	6515-01-540-7226



6 oz MOLLE Retractable GEAR LEASH

- MOLLE Mount compatible / 360° rotating
- High break strength. salt water resistant self-flushing
- 36 in. Extension / Spectra Nylon Line



SPECIFICATIONS

Retracted:
H 3.5 in. x W 1.6 in. x D 1 in.
Cable Length - 36 in.
Weight: 2.0 oz

GEAR LEASH

ITEM #
70-0155



INTRO
MARCH
TAC-EVAC
MED-SURG
RESCUE
TREATMENT SUPPORT
TRAINING
INDEX

NITRILE GLOVES

Nitrile gloves available in Black or Tan colors and come either as pre-packaged pairs (kits) or in bulk (box)

BEAR CLAW®



- Form-fitting, ambidextrous design with textured fingertips
- Thick palm ensures durability
- Textured for improved wet and dry gripping
- Longer and thicker than standard nitrile gloves

BLACK TALON®



SPECIFICATIONS

Kit: L 4.75 in. x W 6 in. x D 4 in.
Weight: 1 lb 1 oz
Box: L 11 in. x W 5 in. x D 3.75 in.
Weight: 2 lb 4 oz

BEAR CLAW GLOVES

CE *NATO NSN

SIZE	PKG QTY	ITEM #	NSN #
S	Kit 25 pair	ZZ-0246	6515-01-619-5350
S	Box 50 pair	70-0009	6515-01-600-6059*
M	Kit 25 pair	ZZ-0247	6515-01-619-5689
M	Box 50 pair	70-0010	6515-01-600-6057
L	Kit 25 pair	ZZ-0248	6515-01-600-0262
L	Box 50 pair	70-0011	6515-01-600-6039
XL	Kit 25 pair	ZZ-0249	6515-01-619-5367
XL	Box 45 pair	70-0012	6515-01-600-6052*

BLACK TALON GLOVES

CE *NATO NSN

SIZE	PKG QTY	ITEM #	NSN #
S	Kit 25 pair	ZZ-0218	----
S	Box 50 pair	70-0001	6515-01-515-0195
M	Kit 25 pair	ZZ-0217	6515-01-672-7977*
M	Box 50 pair	70-0002	6515-01-515-0200
L	Kit 25 pair	ZZ-0216	6515-01-632-0059
L	Box 50 pair	70-0003	6515-01-515-0197
XL	Kit 25 pair	ZZ-0215	6515-01-672-7985*
XL	Box 45 pair	70-0004	6515-01-515-0190

NAR RESPIRATOR

- NIOSH N95 certified particulate barrier
- Adjustable nose bridge for a secure, comfortable fit

SPECIFICATIONS

Packaged:
L 5 in. x W 4.25 in. x H 0.75 in.
Mask size: Adult
Weight: 0.9 oz



NAR RESPIRATOR

ITEM #	NSN #
70-0006	6515-01-537-1531

NAR PPE KIT

1 N95 certified Face Mask & 2 pair Bear Claw Nitrile Gloves

SPECIFICATIONS

L 3.5 in. x W 4 in. x H 1.5 in.
Weight: 2 oz



NAR PPE KIT

ITEM #	NSN #
70-0007	6515-01-537-1534


BIOLOGICAL PERSONAL PROTECTION KIT

Designed to maximize personal protection when handling and managing detainees, refugees, etc.

SPECIFICATIONS

H 6 in. x W 8 in. x D 1 in.
Weight: 4.6 oz

BIOLOGICAL PERSONAL PROTECTION KIT:

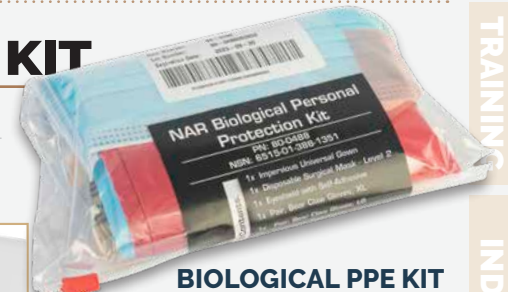

BIOLOGICAL PPE KIT

ITEM #	NSN #
80-0488	6515-01-388-1351



INTRO

M-A-R-C-H

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT SUPPORT

TRAINING

INDEX

PES Polycarbonate Eye Shield

- Meets CoTCCC treatment protocols for penetrating eye trauma
- Convex eye shield accommodating either eye
- Translucent orange shell features perforations to allow air circulation



SPECIFICATIONS

3 in. x 2.38 in. **CE**
Weight: 0.2 oz

POLYCARBONATE EYE SHIELD

UNIT	ITEM #	NSN #
Single	30-0142	6515-01-647-5456
Box of 100	30-0143	6515-01-630-5551

FOX Eye Shield

- Convex rigid eye shield made of lightweight aluminum with perforations allowing air circulation
- Meets CoTCCC treatment protocols for penetrating eye trauma
- Available with or without garter



SPECIFICATIONS

3 in. x 2.38 in.
Weight: 0.2 oz

FOX EYE SHIELD

STYLE	UNIT	ITEM #	NSN #
No Garter	Box of 12	30-0144	6515-01-449-1016
Garter	Single	30-0055	----
Garter	Box of 50	30-5550	6515-01-598-1877

NAR IPRO[®] Tactical Goggle System

High performance eye protection with Hexo-vents for maximum ventilation, fine particle filtration and fog-free viewing

SPECIFICATIONS

Packaged:
L 2.6 in. x W 7 in. x D 3 in.
Weight: 8 oz

NAR IPRO

ITEM #
ZZ-0109



AEGIS ECHO II Ballistic Glasses & Case

- Superior ballistic impact protection per US MIL-PRF-31013 standard clause 3.5.1.1 and 100% protection from UVA/UVB rays
- Includes clear & gray interchangeable lenses

SPECIFICATIONS

Glasses: W 5.5 in. x H 2 in. x D 4.75 in.
Case: L 6.5 in. x W 4.3 in. x H 3 in.
Weight: 9 oz

AEGIS ECHO II

ITEM #
90-0067

OTW GOGGLES Outside the Wire - with Case

- 100% UVA/B/C protection with maximum ballistic protection
- Tapered, neutral tint gray lenses with anti-fog and scratch resistant coating, reduced glare and correct distortion



OTW GOGGLES

ITEM #
90-0066

SPECIFICATIONS

Goggles: W 7.5 in. x H 2.75 in.
Case: L 8.75 in. x W 3.5 in. x H 3 in.
Weight: 10.1 oz

GRAY MAN EyePro

Low visibility EyePro option with classic, subtle styling

GRAY MAN EYEPRO

ITEM #
90-0068



SPECIFICATIONS

Glasses: W 5 in. x H 2 in. x D 5 in.
Case: L 6.625 in. x W 3 in. x H 2.5 in.
Weight: 9 oz



BURNTEC® Hydrogel Dressings

- Sterile non-stick hydrogel dressings for treating 1st, 2nd & 3rd-degree burns (thermal, chemical, electric)
- Also treats other skin injuries, such as sunburns, abrasions, skin irritations, bruises and swelling (insect bites)
- Applied dressings provide a cool, soothing effect for up to 24 hours
- Easy to handle without spilling or damaging the hydrogel
- Permeable to water vapor and oxygen but impermeable to bacteria
- Biocompatible and hypo-allergenic; Reduces risk of hypertrophic scars

BURNTEC DRESSINGS

CE

SIZE	ITEM #	NSN #
2 in. x 2 in.	30-0114	6510-01-631-1913
2.5 in. x 5 in.	30-0115	6510-01-631-2080
4 in. x 4 in.	30-0110	6510-01-631-1883
5 in. x 5 in.	30-0116	6510-01-631-2138
5 in. x 10 in.	30-0117	6510-01-631-1907
8 in. x 8 in.	30-0118	6510-01-631-1897
8.5 in. x 11 in.	30-0119	6510-01-631-1890
4 in. x 16 in.	30-0111	6510-01-631-2124
8 in. x 16 in.	30-0112	6510-01-626-2990
16 in. x 24 in.	30-0120	6510-01-631-1886
10 in. x 10 in. Face	30-0121	6510-01-627-1834
12 in. x 16 in. Face	30-0113	6510-01-627-1855
16 in. x 24 in. Face	30-0122	----

DRY STERILE BURN DRESSING

- Dry sterile, non adherent burn dressing & bandage
- Includes two No. 3 closed safety pins and rubber band individually packaged



DRY STERILE BURN DRESSING

SPECIFICATIONS

Bandage: 45 in. x 45 in. x 63 in.
Pkg: W 4.5 in. x H 4 in. x D 0.5 in.
Weight: 1 oz

ITEM #	NSN #
30-0162	6510-01-587-6579



BurnTec® BURN DRESSING KIT

- Convenient array of BurnTec® dressings in a handy carrying case
- Kit contains (1) each of the following dressings:

2 in. x 2 in.	5 in. x 10 in.	8 in. x 16 in.
2.5 in. x 5 in.	8 in. x 8 in.	10 in. x 10 in. (face)
4 in. x 4 in.	4 in. x 16 in.	12 in. x 16 in. (face)
5 in. x 5 in.	8.5 in. x 11 in.	

BURNTEC BURN KIT

OPTION	ITEM #	NSN #
Kit	80-0361	6515-01-630-9177
Bag Only	30-0123	6510-01-626-6079



online video available

A scan to view

SPECIFICATIONS

CE

L 12 in. x W 15.5 in. x D 3.5 in.
Weight: 6 lb 4 oz

BurnTec® MINOR BURN DRESSING KIT

- Variety of BurnTec® dressings in a handy carrying case
- Kit contains the following:
 - Two 4 in. x 4 in. dressings
 - Two 4 in. x 16 in. dressings
 - One 8.5 in. x 11 in. dressing



BURNTEC MINOR BURN KIT

OPTION	ITEM #	NSN #
Kit	80-0415	6510-01-644-5427
Bag Only	30-0150	----

SPECIFICATIONS

CE

L 14 in. x W 8.5 in. x D 2.5 in.
Weight: 2 lb 13 oz







TREATMENT SUPPORT | Triage System

T2® TAG NAR

- Triage & treatment tag utilizing the S.A.L.T. Mass Casualty Triage Standard
- To aid in documenting critical patient information throughout the entire mass casualty response process
- Liquid-resistant in all settings; made of rugged materials contained in a plastic sheath with latex-free band

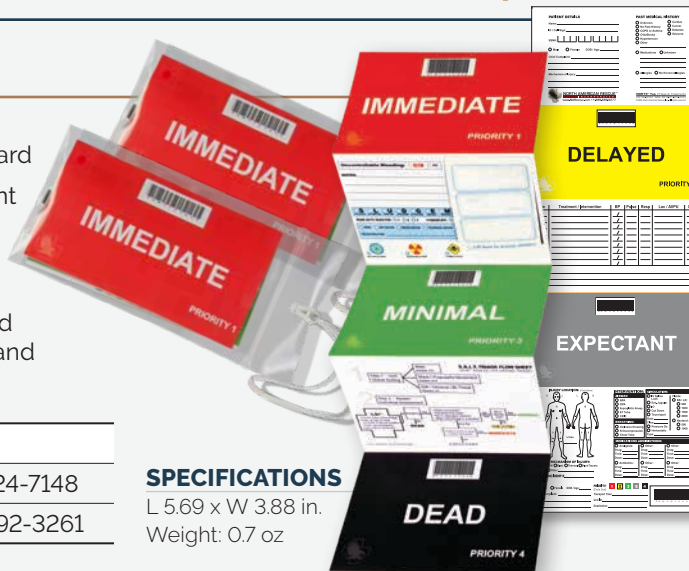

T2 TAG:

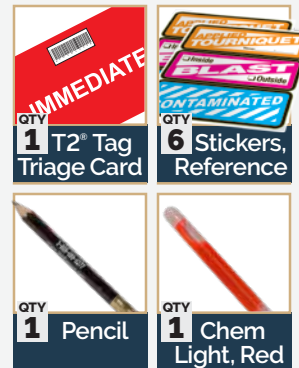

T2 TAG

UNIT	ITEM #	NSN #
20 Pack	20-0031	6515-01-624-7148
1000 Pack	20-0032	6545-01-592-3261

SPECIFICATIONS

L 5.69 x W 3.88 in.
Weight: 0.7 oz

TACTICAL TRIAGE KIT T2®

- Compact incident command system
- Streamlines medical triage
- Contained in a lightweight, rugged vacuum-sealed pouch

SPECIFICATIONS

L 5.25 x W 7 in. x D 1.25 in.
Weight: 4.9 oz

TACTICAL TRIAGE KIT

ITEM #	NSN #
20-0040	6545-01-591-2134

TACTICAL TRIAGE KIT:



COMBAT CASUALTY CARD

- Streamlined triage card that meets all TCCC standards for casualty management
- State-of-the-art ruggedized paper is tear resistant, waterproof & designed to fit into the IFAK
- Includes grommet reinforced punch with latex-free band for attachment
- Currently in use by Special Operations Forces



SPECIFICATIONS

H 6 in. x W 4.5 in. each
Weight: 0.22 oz ea | 5.5 oz 25-Pack

COMBAT CASUALTY CARD

UNIT	ITEM #	NSN #
25 Pack	ZZ-0262	6515-01-587-1700

COMBAT CASUALTY REFERENCE CARD

- Rugged reference card for emergency management featuring the current TCCC Guidelines
- Clear, concise instructions for medics and warfighters
- Information divided into four sections:
 - Care Under Fire
 - Tactical Field Care
 - TACEVAC Protocols
 - 9-line MEDEVAC request/notes
- Pocket sized for convenience
- Eighteen pages



COMBAT CASUALTY REFERENCE CARD

ITEM#	NSN #
ZZ-0057	6515-01-537-4161

SPECIFICATIONS


H 6 in. x W 5 in. x D 0.1 in.
Weight: 2.3 oz



TRIAGE OFFICER KIT T2®



T2 TRIAGE OFFICER KIT:

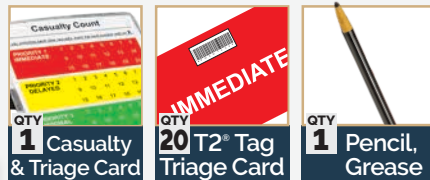

- Mid-range solution providing initial triage of up to 20 casualties
- Compact bag with mesh storage pockets equipped with T2/SALT triage cards
- Includes a laminated S.A.L.T. Triage Flowchart & Casualty Count Card

SPECIFICATIONS

Pkg: H 10 in. x W 8 in. x D 5 in.
Weight: 1 lb 14 oz

TRIAGE OFFICER KIT

COLOR	ITEM #	NSN #
OD Green	20-0023	6545-01-537-4189
Digital Camo	20-0075	----
Red	20-0076	----

COMMAND KIT T2®

- Organized, rapid response system for Incident Commanders during a mass casualty/tactical triage incident
- 3 laminated, detachable Command Boards
- Pad of 50 Patient-Tracking Boards
- 3-ring, clear storage sheets for bar codes & mesh pockets to store extra accessories

SPECIFICATIONS

Case Closed: L 21 in. x W 17 in. x D 6 in.
Weight: 4 lb 14 oz



COMMAND KIT

COLOR	ITEM #	NSN #
OD Green	20-0025	6545-01-537-4197
Digital Camo	20-0074	----
Red	20-0026	----

T2 COMMAND KIT:



Mass Casualty Incident Management MCIM KIT T2®



- Comprehensive triage system ideal in disaster response, field hospitals & military troop concentrations.
- Lightweight, modular, man-portable triage kit.
- includes user-friendly tools for every aspect of the uniformed S.A.L.T. model for disaster triage

SPECIFICATIONS

H 21 in. x W 17 in. x D 9 in. Weight: 15 lb 2 oz

MCIM KIT

COLOR	ITEM #	NSN #
OD Green	20-0029	6545-01-537-3138
Digital Camo	20-0073	----
Red	20-0077	----

T2 MCIM:





OPERATORS CAP NAR

- 100% cotton
- Adjustable strap
- 3.5 in. by 2.5 in. cut loop on the front designed for attaching patches (sold separately)
- 1 in. square of cut loop on the top for attaching an IR square

SPECIFICATIONS

One size fits most
Weight: 3.2 oz

OPERATORS CAP

ITEM#

ZZ-0092



EMBROIDERED LOGO CAP NAR

- All American Cap
- Made in the USA
- Comfortable, attractive Fit-4-You™ design fits most wearers
- American flag on the front, NAR logo on back

SPECIFICATIONS

Fit-4-You™ S/M
Weight: 3.2 oz

EMBROIDERED LOGO CAP

ITEM#

ZZ-0176



MEDIC CROSS PATCHES PVC

- Soft flexible PVC rubber patch
- Backed with hook material for easy attachment and removal
- Ideal for use on medical kits such as the M-FAK and T.O.R.K.
- 1 in. x 1 in. available at NARescue.com

MEDIC CROSS PATCHES

COLORS	ITEM #
Coyote & Black	ZZ-0797
Black & Gray	ZZ-0798
OD Green & Black	ZZ-0800
Black & Orange	ZZ-0799



SPECIFICATIONS

L 2 in. x W 2 in. x D 0.125 in.
Weight: 0.6 oz

IR FLAGS NAR

- Clearly visible under IR light
- Hook-and-loop attachment
- Sold in Packs of 2

SPECIFICATIONS

L 2 in. x W 2 in. x D 0.125 in.
Weight: 0.6 oz

IR FLAGS

COLOR	STARS	ITEM #
Full Color	L	ZZ-0093
Coyote	L	ZZ-0094
Coyote	R	ZZ-0235
OD Green	L	ZZ-0107



MED ID PATCH LUMINOUS

Glow in the Dark

- Soft flexible custom rubber PVC
- Ultimate glow-in-the-dark identification technology
- 1 in. block letters
- 2 medical patches per package

SPECIFICATIONS

Each: L 3 in. x W 2 in. x D 0.125 in.
Weight: 0.4 oz

LUMINOUS MED ID PATCH

ITEM#

ZZ-0121



IR SQUARES NAR

- Infrared patch printed on honeycomb reflective material
- Black hook fasteners on the back for easy attachment
- Clearly visible under IR light illumination
- Six squares per pack

SPECIFICATIONS

Each: 1 in. x 1 in.
Weight: 0.2 oz.

IR SQUARES

ITEM#

ZZ-0108





NORTH AMERICAN RESCUE®

Education and Training




North American Rescue Education and Training (NAR Training)

has roots that run deep and wide in the world of pre-hospital medicine and technical rescue. Our team was built from the ground up on the vision and experience of seasoned military medics seeking to fill the training gaps commonly encountered in the operational environment.

Since then, we have continued to broaden our vision to include other invaluable first responder perspectives such as pre-hospital healthcare providers in civilian, law enforcement, emergency medical services and community public safety initiatives.

CORE Competencies

- Continuing Medical Education (CMEs)
- Prolonged Casualty Care
- Austere Medicine
- Technical Rescue
- Emergency Medical Services
- Critical Care Flight Medicine

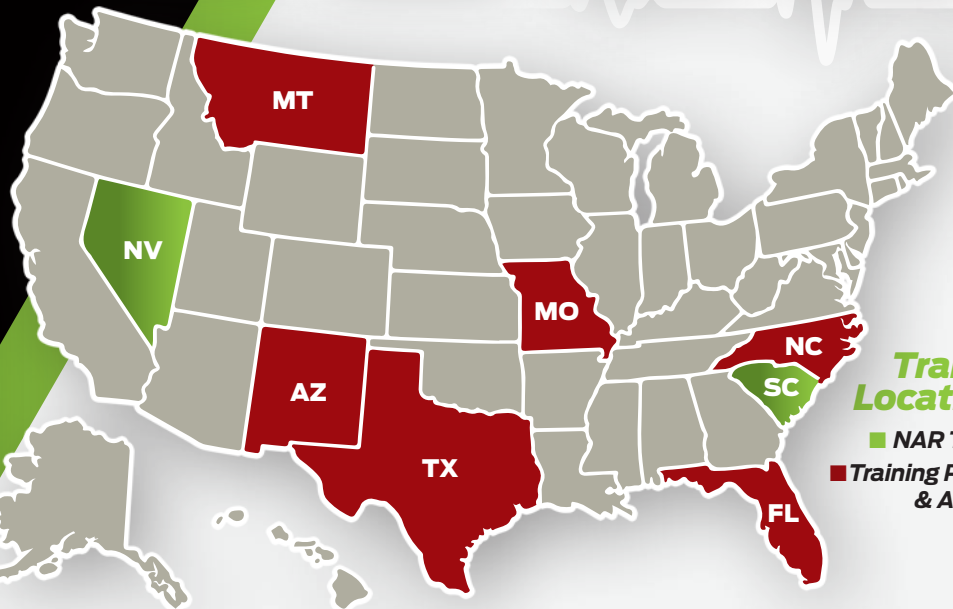

Training Locations:

- NAR Training
- Training Partners & Affiliates

ACCREDITATIONS:



American Heart Association®

Now SAME DAY Cards Issued At End Of Course!



Explore training opportunities online

NARTraining.com

➤ Phone: 888.689.6277 ➤ Email: NARTraining@NARescue.com

INTRO

MAR-CCH

KITS

TAC-EVAC

MED-SURG

RESCUE

TREATMENT-SUPPORT

TRAINING

INDEX

COURSES

PUBLIC ACCESS



➤ **NAR BCON**

6-hour course oriented towards civilians; students are exposed to realistic trauma scenarios and taught how to treat major hemorrhaging injuries, receiving critical knowledge on how respond to preventable bleeding deaths prior to arrival of EMS. Stop the Bleed certification issued upon completion.

➤ **First Responder - Work Place Safety**

8-hour course designed for civilians. Students will learn CPR, bleeding control, and basic first aid to the most common workplace injuries. AHA BLS and NAR First Responder certification upon successful test completion.

➤ **Civilian Medical**

3-day medical course designed for civilians. Students will learn basic emergency medicine, med bag packing, and med supply procurement. AHA BLS and NAEMT TECC certification upon successful test completion.

➤ **First Responder - Trauma Care**

8-hour course for law enforcement and first responders that will provide students the ability and confidence to apply trauma care in the field. Upon completion of testing, students receive NAEMT TECC certification.

IBSC

➤ **NAR C3**

42-hour in-house course for IBSC critical care course certification or recertification (FP-C or CCP-C) and national registry recertification. Students will learn or refresh all skills and knowledge needed to certify/re-certify, and will earn 42 hours of IBSC and National Registry Continuing Education.

➤ **Vent Hero**

Supporting literature that provides an introduction into the theory of successful ventilation management geared towards the critical care provider and clinician responsible for managing complex patient ventilation support.

➤ **Essentials In Critical Care 4.0**

Online refresher course which exceeds the minimum requirements (above 16 hrs) for approved IBSC certifications.

NATIONAL REGISTRY

➤ **EMT Hybrid**

Self-paced (over period of one year), online learning designed to prepare entry-level student for written and practical application testing for National EMT testing and certification.

➤ **AEMT Hybrid**

Self-paced (over period of one year), online learning designed to prepare EMT student for written and practical application testing for National AEMT testing and certification.

LAW ENFORCEMENT

➤ **TECC-LE**

8-hour course specifically designed for law enforcement. Covers hemorrhage control, airway control, splinting, extraction techniques, and stop the bleed scenarios. Upon passing test, students receive NAEMT TECC-LE certification.

➤ **FEDMED**

5-day course designed specifically for federal law enforcement officers. Students learn how to diagnose and treat various illnesses and injuries utilizing role players and simulations. Upon successful completion, AHA BLS and NAEMT TECC certifications are achieved; for EMTs, this course also serves as a full EMT re-certification.



NORTH AMERICAN RESCUE®

Education and Training



➤ **Phone:** 888.689.6277 ➤ **Email:** NARTraining@NARescue.com



MILITARY

➤ TCCC-AC

This 40-hour course is designed for non-medical military personnel that provides a 4-year NAEMT TCCC certification upon completion. Students will learn emergency battlefield medicine focusing on preventable causes of death.

➤ TCCC-MP

This 16-hour course is designed for non-medical military personnel that provides a 4-year NAEMT TCCC certification upon completion. Includes standard training with addition of scenarios utilizing role players and moulage.

➤ EGT-AC

This 5-day course is designed for small military units that provides a 4-year NAEMT TCCC-AC certification. Provides students with TCCC medicine and drills, field maneuvers, and trauma lanes with role players and CASEVAC trainers.

➤ EGT-MP

This 5-day course is designed for military medical provider units that awards a 4-year NAEMT TCCC-AC certification upon completion. Provides students with TCCC medicine and drills, field maneuvers, and trauma lanes with role players and CASEVAC trainers.

➤ SOFMED

This 8-day course is specifically designed for Special Operations Medics that provides the following certifications: BLS, ACLS, PALS, TCCC-MP, and a full NCCP 60-hour CEU/NRP refresher. Students focus on medicine while enhancing tactics via use of simulations, role players, and trainers.

➤ SOTM

This 8-day course is specifically designed for Special Operations Medics that provides the following certifications: BLS, ACLS, PALS, TCCC-MP, and a full NCCP 60-hour CEU/NRP refresher. Students focus on medicine while enhancing tactics via use of perfused cadavers, simulations, role players, and trainers.

➤ Toxic MED 101

2-day course for anyone who may be involved in a CBRN incident. Students will receive critical knowledge and skills on how to respond to incidents involving chemical/biological/radiological/nuclear contaminants. NAEMT TECC certification awarded upon completion.

➤ Toxic MED 201

5-day course designed for medical providers/first responders. Provides background, knowledge and methodologies associated with responding to or supporting CBRN operations with casualties. Students will have the opportunity to receive the AHLS certification if desired. All students will earn NAEMT TECC certification upon successful completion.

ADVANCED

➤ Austere Rescue Course

5-day course aimed at military medical providers and paramedics. Students will receive periods of instruction covering austere medicine, remote rescue, and patient transport in rugged terrain with limited resources. Course will involve treatment and movement in remote environments coupled with challenging problems involving emergency medicine, transport, and elements of PCC.

➤ Austere Surgical Course

5-day course exclusively designed for military surgical teams. Students will utilize perfused cadavers to hone surgical skills in a variety of environments; battlefield conditions and equipment will be incorporated into training to increase proficiency and team skills. SOPs will be tested in preparation for full mission profile event.

➤ Prolonged Critical Care

5-day course that gives military medical providers the knowledge and skills necessary to treat patients in an austere environment for longer periods of time with limited resources when patient evacuation is not tenable.

Students will be instructed on how to keep patients stable with whatever resources are available until transport arrives.



Explore training opportunities online at [NARTraining.com](https://www.NARTraining.com)

Tactical Operation
Medical Manikin SERIES

TOMManikin™

BASIC | BLAST | GSW | WATER | CONFINED SPACE | BURN | CBRNE | FEMALE

Rugged, state-of-the-art training simulators, NARS TOMManikin Series provides the comprehensive solution for enhanced training in Tactical Combat Casualty Care (TCCC) from Point of Injury to transfer of higher care!

- **All-in-one, hyper-realistic trauma simulation solution** enabling instructors to build **multiple realistic traumatic casualty scenarios** with **reactive patients**
Ideal for use in classroom settings, trauma lanes and field training scenarios
- **Multiple, interchangeable extremities** present a variety of wounding patterns to provide **full mission profiles** in **combat scenarios** and **disaster/mass casualty re-creations**

- **Solid metal bones** and **joints** constructed of **welded steel** for **unparalleled durability** and **quick exchange leg joints** for **flexibility** in the most **austere environments**

- **Proprietary silicone skins** available in **light, medium** or **dark skin tones** that **will not break down** during **prolonged UV exposure** or **extreme temperatures**
- **180 lb** manikin provides extremities capable of **withstanding up to 6,000 lb** of pressure -
Ideal for **building collapse/entrapment training scenarios** beyond conventional exercises

- **Ruggedized remote controls** provide **simultaneous long-range communications** with **multiple NARS devices** all from a **single device!**

SPECIAL FEATURES

Interactive Animatronics:

- Lifelike Breathing Mechanisms
- Functional Anatomical Airway
- Variable Radial Pulse Rates & Pressure
- High Fidelity Bleeding from Multiple Wounds
- Customizable, Multi-language Voice Tracks
- Rugged Softech® Skin in a variety of skin tones

Wireless Controls:

- MILSPEC Wireless Remote Instructor Tablet & Software
- Stand alone point-to-point communications with .5 mile line of site range

Flexible Wound Configurations:

Basic | Blast | GSW | Water | Confined Space | Burn | CBRNE
Custom configurations including Female Anatomy on request



Clinical Interventions available on ALL TOMManikins:

MASSIVE HEMORRHAGE:

- **Tourniquet Application** (proximal, junctional & abdominal)
- **Active Bleeding Control** by **Direct Pressure** on Pressure Points
- **Wound Packing**

ADVANCED AIRWAY:

- **Nasopharyngeal/Oral/Supraglottic Airway** Insertion
- **Endotracheal Intubation**
- **Surgical Airway** and **Tissue Dissection** with Realistic Bleeding

RESPIRATIONS:

- **Chest Seal Application**
- **Bilateral Needle Decompression**
- **Chest Tube Capable**

CIRCULATION:

- **IO** (sternal & humeral)
- **Antecubital IV Access** on uninjured upper limbs



NORTH AMERICAN RESCUE®

Simulation



Every TOMManikin™ includes:

- **ACCESSORY KIT**
 - 2 x NiMH rechargeable batteries (5-amp & 10-amp) with independent chargers
 - 1 x T-handle wrench
 - Replacement connectors
 - 5-pack of surgical airway neck skins (35 op sites)
 - 8 oz concentrated sim blood
 - NARS mil spec-rated tablet

- 2 x 3-wheeled transport and storage carriers designed for TOMM
- 1 x NAR IFAK
- **VELCRO UNIFORM**
 - Digital pattern
 - Sand color scheme
 - Blouse features Velcro sleeves and lateral seams
 - Trousers feature Velcro inseam and lateral seams
 - Button fly



MILSPEC INSTRUCTOR TABLET with Scenario Preset Builder Software

Features:

- Features proprietary communication link that requires NO Wi-Fi, Bluetooth, or LAN connections
- Simultaneous long-range communications with numerous NARS devices
- NO interference risk with telemetry, avionics, or weapon systems



MILSPEC Rated:

- **Operating:** -20 to 60°C (-4 to 140°F) MIL-STD-810G, Method 501.5 Procedure II
- **Storage:** -40 to 70°C (-40 to 158°F) MIL-STD-810G, Method 501.5/502.5 Procedure I
- **Water:** IP65, IEC 60529
- **Drop/Shock:** 26 drops from 1.2 m (4 ft) MIL-STD-810G, Method 516.6 Procedure IV
- **Vibration:** MIL-STD-810G, Method 514.6 Procedures I, II
- **Sand/Dust:** IP65, IEC 60529
- **Humidity:** 0 ~ 95% (non-condensing) MIL-STD-810G, Method 507.5
- **Altitude:** 4572 m (15,000 ft) MIL-STD-810G, Method 500.5 Procedure I



- **Tourniquet Application** (proximal, junctional, and abdominal tourniquets)
- **Active bleeding** controllable by **direct pressure** on pressure points
- **Wound packing** with realistic bleeding control
- **Nasopharyngeal and oral airways, supraglottic airway** insertion and **endotracheal intubation**
- **Surgical airway** with realistic bleeding and tissue dissection 35 operative sites included
- Remote controlled **breathing rates:** Slow, Normal, and Fast
- **Bilateral Needle Decompression** at either 2nd/3rd intercostal (mid-clavicular), or 5th/6th intercostal (mid-axillary) locations
- Realistic **chest seal application** & adhesion
- **Chest Tube** capable (bilateral)
- **Radial pulses** in uninjured upper extremities & bilateral carotid pulses
- **Variable pulse rate & pressure** (Normal, bounding, thready)
- **IO** (sternal & humeral) and **antecubital IV access** on uninjured upper limbs




All Manikins include a MILSPEC-rated instructor tablet

BASIC TOMManikin™

BASIC Wounding Patterns:

- 1. Partial Arm Amputation below the elbow** with massive bleeding and dramatic soft tissue damage, exposed hard bone and discernible tissue layers Responds to tourniquet application or brachial artery direct pressure application
- 2. Full Amputation Leg** with massive bleeding, extensive soft tissue damage, exposed hard bone, visible tissue layers. Responds to direct junctional pressure and tourniquet application

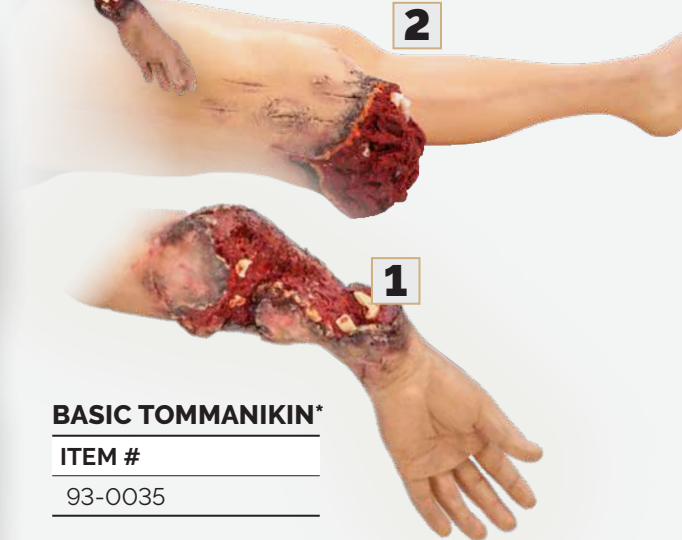

BASIC TOMMANIKIN*

ITEM #

93-0035

*Other configurations are available upon request

BLAST TOMManikin™

BLAST Wounding Patterns:

- 1. Open Jaw Fracture Head Skin** with massive maxillofacial bone and soft tissue injuries, actively bleeds. Supports OP/NP/ET/Surgical airways
- 2. Packable Abdominal Evisceration** with intricate handcrafted attention to anatomical detail with visible tissue layers, and 10 in. of eviscerated bowel with peripheral ecchymosis surrounding large open wound. Massive bleeding with "clampable" vessel
- 3. Shrapnel Wound Arm** with active bleeding, multiple soft tissue wounds, responds to tourniquet application
- 4. Full Amputation Leg** with massive bleeding, extensive soft tissue damage, exposed hard bone, visible tissue layers. Responds to direct junctional pressure and tourniquet application



BLAST TOMMANIKIN*

ITEM #

93-0036

*Other configurations are available upon request

GSW (Gun Shot Wound) TOMManikin™

GSW Wounding Patterns:

1. Sucking Chest Wound with anterior GSW entry & posterior exit. Actively bleeds from both locations. Anterior wound mixes with air from onboard compressor to create bubbling bleeding effect. Supports chest seal application and pleural needle decompression

2. GSW to Groin with GSW anterior entry/posterior exit with soft tissue trauma to penis and avulsed scrotum with prolapsed testicle. Wounds present with peripheral wound site ecchymosis, exposed sub-q fat layer and muscle striation. Entry wound is packable and responds to both peripheral and junctional tourniquet application and femoral pressure point application

3. GSW Leg with through & through GSW to leg above the knee. Entry wound presents with small entrance hole with circular powder burn. Exit wound has exposed fat and wound site ecchymosis. Variable bleeding rate. Responds to both peripheral and junctional tourniquet application and femoral pressure point application. Can be used as an impaled object using steel rebar, or wooden stake

4. Full Amputation Leg (Option)** with massive bleeding, extensive soft tissue damage, exposed hard bone, visible tissue layers. Responds to direct junctional pressure and tourniquet application

*Other configurations are available upon request

Addressing injury locations based on US Army study suggesting females who suffer penetrating injuries to chest, abdomen, pelvis have a higher mortality rates than male colleagues

TAMikin GSW

Female version of TOMManikin GSW, with anatomically correct female features and voice

- Injuries located in sensitive areas such as under breast & pubic area as well as the leg
- Parallels GSW TOMManikin's capabilities & wounding patterns

**Leg Amputation Option
93-0369/93-0947

GSW TOMMANIKIN*

OPTION	ITEM #
GSW TOMM (Basic)	93-0040
GSW TOMM (w/Leg Amputation)	93-0369
GSW TAMikin*** (Basic)	93-0039
GSW TAMikin** (w/Leg Amputation)	93-0947

***Female version



*Other configurations are available upon request

WATER TOMManikin™

Wounding Patterns:

Partial Leg Amputation below the knee with massive soft tissue trauma and exposed hard bone, muscle and tissue layers.

Options:

Sucking Chest Wound with anterior GSW entry & posterior exit. Actively bleeds from both locations. Anterior wound mixes with air from onboard compressor to create bubbling bleeding effect. Supports chest seal application and pleural needle decompression

-or-

Packable Abdominal Evisceration with intricate handcrafted attention to anatomical detail with visible tissue layers, and 10 in. of eviscerated bowel with peripheral ecchymosis surrounding large open wound. Massive bleeding with "clampable" vessel

Closed Fracture Arm with ecchymosis and palpable crepitus.

Designed to be submerged up to 33 feet in either fresh or salt water

WATER TOMMANIKIN*

OPTION	ITEM #
Basic (Sucking Chest Wound)	93-0041
ABD Evisceration Closed Arm Fracture	93-0108



**Other configurations are available upon request*

CONFINED SPACE TOMManikin™

Wounding Patterns:

Sucking Chest Wound with anterior GSW entry & posterior exit. Actively bleeds from both locations. Anterior wound mixes with air from onboard compressor to create bubbling bleeding effect. Supports chest seal application and pleural needle decompression

Partial Leg Amputation below the knee with massive soft tissue trauma and exposed hard bone, muscle and tissue layers.

Shrapnel Wound Arm with active bleeding, multiple soft tissue wounds, responds to tourniquet application



CONFINED SPACE TOMMANIKIN*

ITEM #
93-0110

**Other configurations are available upon request*



BURN TOMManikin™

Wounding Patterns:

Burn Face* with 1st, 2nd, 3rd-degree burns to face that supports all airway adjuncts. *Head skin is interchangeable with other TOMManikins*

Burn Arms with 1st, 2nd, 3rd-degree circumferential burns to both arms

BURN TOMMANIKIN*

ITEM #

93-0037

**Other configurations are available upon request*

****Head skin is interchangeable with other TOMManikins**



CBRNE TOMMANIKIN*

ITEM #

93-0038

**Other configurations are available upon request*



CBRNE TOMManikin™

Wounding Patterns:

Burn Face* with 1st, 2nd, 3rd-degree burns to face that supports all airway adjuncts. *Head skin is interchangeable with other TOMManikins*

Sucking Chest Wound with anterior GSW entry & posterior exit. Actively bleeds from both locations. Anterior wound mixes with air from onboard compressor to create bubbling bleeding effect. Supports chest seal application and pleural needle decompression

Abdominal Evisceration (Wearable Moulage) with intricate handcrafted attention to anatomical detail with visible tissue layers, and 10 in. of eviscerated bowel with peripheral ecchymosis surrounding large open wound. Massive bleeding with 'clampable' vessel.

GSW to Groin with GSW anterior entry/posterior exit with soft tissue trauma to penis and avulsed scrotum with prolapsed testicle. Wounds present with peripheral wound site ecchymosis, exposed sub-q fat layer and muscle striation. Entry wound is packable and responds to both peripheral and junctional tourniquet application and femoral pressure point application

K9 Medical Manikin

K9 Features:

- Redesigned all-metal frame
- Expanded joint range of motion
- CPR chest
- 5 interchangeable limbs with various injuries: Complete Amputation, Open Fracture, Soft Tissue Injuries, Road Rash/Burns, & Gunshot Wound
- Improved airway
- Barks, Bleeds & Breathes
- 1L onboard blood capacity
- Wireless operation via NARS tablet



K9 MANIKIN

ITEM #

93-0027

INTRO
M-A-R-C-H
KITS
TAC-EVIC
MED-SURG
RESCUE
TREATMENT SUPPORT
TRAINING
INDEX

Deluxe ALS KIT

Designed to simulate various vital signs for advanced life support monitoring and simulation

Features:

- Simulated NIBP cuff, SpO₂ finger probe, and ETCO₂ inline detector with waveform and digital value
- Simulated cardiac monitor screen with Lead II tracing rhythm
- ETCO₂ & SpO₂ waveforms, NIBP. Fully functional monitor screen for defibrillation/AED mode, alarms, and rhythm synchronization
- All student-visible values controlled via NARS Instructor tablet (included)
- NARS software includes exportable student evaluation performance documentation



Disconnect your student's dependency on instructor "over-the-shoulder" input

PATIENT MONITORING SUITE with NARS Instructor Tablet

NARS Patient Monitoring Suite of simulators can be used in all phases of patient care for any human (role-player), or ANY patient simulator


Features include:

- **NARS Capnograph Simulator**
Designed to be used as a mainstream carbon dioxide monitor that can be connected in-line with a mask or tube. Simulator displays include carbon dioxide partial pressure, respiratory rate and CO₂ waveform.
- **NARS Pulse Oximeter Simulator**
Utilizes the widely used finger-style probe for monitoring blood oxygen saturation. Simulator display includes Heart Rate, SpO₂ and SpO₂ waveforms.
- **NARS NIBP Monitor Simulator**
Designed to replicate a portable NIBP monitor. This simulator allows medical responders to gather vital signs without instructor intervention. Onboard controls allow various modes of operation. Display includes systolic & diastolic pressures (mmHg) and heart rate.

DELUXE ALS KIT

ITEM #

93-0071

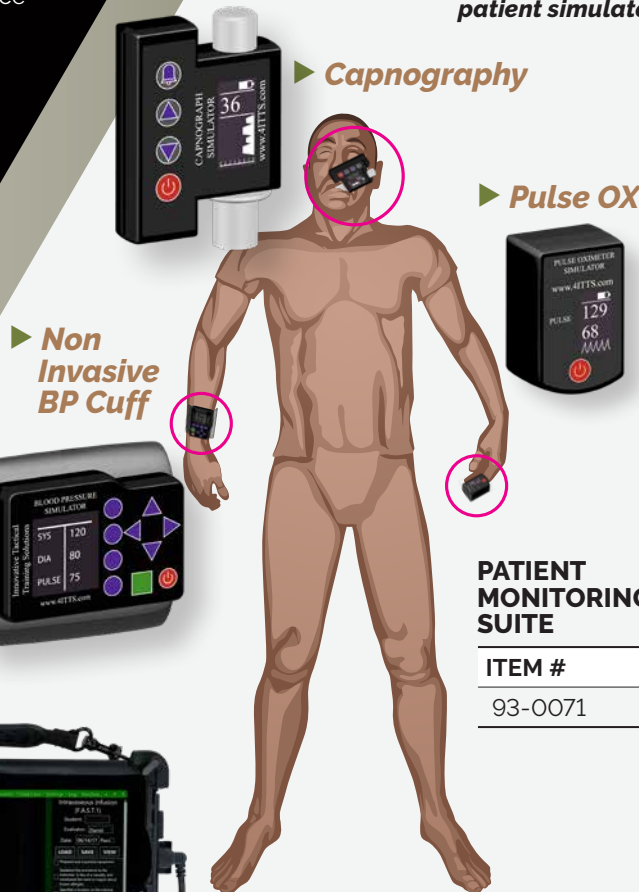

▶ Capnography

▶ Pulse OX

▶ Non Invasive BP Cuff

PATIENT MONITORING SUITE

ITEM #

93-0071

◀ Includes NARS Instructor Tablet





MOBILE FEEDBACK SYSTEM

**requires AC power*

Standalone mobile platform to capture simulation events*

Features include:

- Pelican 1650 case with custom-cut foam to protect all system components
- 2 x Cameras
- 2 x 72 in. Telescopic Tripods with camera mounts
- 3 x Terabyte stand-alone recorder
- 1 x NARS mil-spec rated Instructor Tablet for wireless communications
- Proprietary software with instructor bookmark/replay feature, simulator status capture, trainee biometric capture, trainee evaluation software (exportable reports)
- 2 x 25 ft. CAT-6 retractable cord reels
- Trainee-worn stress monitor
- Features streaming participant heart and respiratory rate



Wearable Trainee Stress Monitor



CASEVAC Simulators **HELOS**

Features:

- Available in **H60, CV22 & MH6** configurations
- Reusable, all-weather systems
- Includes loudspeakers for inflight noise, smoke generator to replicate inflight emergencies, wind machine to simulate rotor wash
- DOT-rated for over-the-road towing
- Platform-specific options include HD cameras (with playback system), FRIES bar, litter stanchions, blackout window drapes, ICS Kit (with headsets and cords)



NEW H-60 ▲
Cockpit with Realistic Pilot/Co-Pilot Seats and Graphics Panel ▶



MOBILE FEEDBACK SYSTEM

ITEM #

93-0945



TCCC SKILLS SYSTEM

All Combatant

- All the necessary tactical medical and casualty extraction/evacuation equipment needed when teaching the TCCC – All Combatant Skills required tasks for completing the Tactical Combat Casualty Care (TCCC) course
- A turnkey solution for training centers and educational institutions designed to equip up to twenty-five (25) students
- Complete system consists of three self-contained Hard Storage Case modules: Skill Stations | IFAKS | Evacuation

SPECIFICATIONS

Hard Storage Case:

External – L 32.5 in. x W 21 in. x D 13 in.

Internal – L 29 in. x W 18 in. x D 10 in.

Weights: TCCC Skills Set: 55 lb

TCCC IFAK Skills Set: 70 lb

TCCC Evacuation Skills Set: 50 lb

TCCC SKILLS SYSTEM

UNIT	ITEM #
3 Modules	80-0970

TCCC SKILLS SET

All Combatant

- Train up to twenty-five (25) students at individual skill stations in the treatment of: *extremity hemorrhage, wound packing, chest seal application, needle chest decompression, basic airway adjuncts, splinting, hypothermia prevention & management, and eye injuries*
- Self-contained, wheeled, rigid case with internal dividers to keep the equipment organized and easy to inventory

SPECIFICATIONS

Hard Storage Case:

External – L 32.5 in. x W 21 in. x D 13 in.

Internal – L 29 in. x W 18 in. x D 10 in.

Weight: 55 lb

TCCC SKILLS SET

MODULE	ITEM #
TCCC Skills Set	80-0968



TCCC SKILLS SET MODULE:



All Combatant TCCC IFAK SKILLS SET

TCCC IFAK SKILL SET MODULE:

- Twenty-five (25) pre-packed Training IFAKs for use during trauma lanes or training scenarios
- IFAKs come packaged in a self contained wheeled hard, case that allows for easy transport and storage



QTY 25 TCCC IFAK Skills Kit

Kit contents below ▼

SPECIFICATIONS

External Case – L 32.5 in. x W 21 in. x D 13 in.
Internal Case – L 29 in. x W 18 in. x D 10 in.
Weight: 70 lb

TCCC IFAK SKILLS SET MODULE

MODULE	ITEM #
TCCC IFAK Skills Set	80-0969



TCCC IFAK SKILLS KIT All Combatant

SPECIFICATIONS

H 7 in. x W 4.5 in. x D 2.75 in.
Weight: 2 lb 8 oz

TCCC IFAK SKILLS KIT

OPTION	ITEM #
TCCC IFAK Skills Kit	85-1958
TCCC IFAK ReSupply*	80-1028

TCCC IFAK SKILLS KIT:

 QTY 1 C-A-T® Tourniquet, Trainer	 QTY 1 NARETD™ 6 in. Flat	 QTY 1 Gauze S-Rolled	 QTY 1 Combat Gauze® Trainer
 QTY 1 28F Naso Airway & Lube	 QTY 1 HyFin® Trainer, Twin Pk	 QTY 1 ARS® Needle, 14g	 QTY 1 Surgical Tape, 2 in.
 PAIR 1 Nitrile Gloves, Lg	 QTY 1 Permanent Marker, Large	 QTY 1 Combat Casualty Card	

* TCCC IFAK ReSupply Kit contains all components of the IFAK SKILLS Kit, excluding the Nylon Bag



TCCC EVACUATION SKILLS SET

- Various high threat extraction devices and casualty evacuation systems to train students in different options to move casualties during the all three TCCC Phases of Care
- Packaged in a self-contained, wheeled, rigid case that allows for easy transport and storage

TCCC EVACUATION SKILL SET MODULE:

 QTY 1 HPMK®	 QTY 1 Talon® 90C Litter	 QTY 1 R.E.E.L. Litter Black
-----------------	-----------------------------	---------------------------------

SPECIFICATIONS

External Case – L 32.5 in. x W 21 in. x D 13 in.
Internal Case – L 29 in. x W 18 in. x D 10 in.
Weight: 50 lb

 QTY 1 Hasty Harness	 QTY 1 Rat™ Strap	 QTY 1 Dragon Handle System	 QTY 1 Carabiner Non-locking
-------------------------	----------------------	--------------------------------	---------------------------------

TCCC EVACUATION SKILLS SET MODULE

MODULE	ITEM #
TCCC Evacuation Skills Set	80-0971





NAVAL FIRST AID BOX RESPONSE KIT *All Service Member Trainer*

- Designed to meet the TCCC requirements of the U.S. Naval Ships, Military Sealift Command ships and afloat commercial vessels
- Provides the All Service Members and first responders with rapid access and organization of critical lifesaving equipment for rendering on scene emergency care at sea in a training environment
- Lightweight, sand/water resistant, slimline pouch constructed of Berry Amendment 500D Cordura solution dyed nylon, Ghille Tex hardware, self repairing zippers & silent 550 cord pulls
- Pouch features securable anti-pilferage seals and a clear window to readily view contents
- Easily fits into the standard Naval First Aid Box

Naval First Aid Box (not included)



NAVAL FIRST AID BOX RESPONSE KIT TRAINER

ITEM #	NSN #
85-3204	6910-01-695-9363

SPECIFICATIONS

L 15 in. x H 6 in. x D 4 in.
Weight: 2 lb

NAVAL FIRST AID BOX RESPONSE KIT TRAINER:

 QTY 1 C-A-T® Tourniquet, Trainer	 QTY 1 Gauze, Compressed	 QTY 1 Gauze Wound Packing	 QTY 2 ETD™6 in. Dressing	 QTY 1 HyFin® Trainer, Twin Pk
 QTY 1 28F Naso Airway	 QTY 1 CPR Mask, Res-Cue Key	 QTY 1 Survival Blanket	 QTY 1 BurnTec 4"x16" Dressing	 QTY 1 SAM® Splint II
 QTY 1 Triangular Bandage	 QTY 2 Elastic Bandage, 3in.	 QTY 1 Fox Eye Shield w/garter	 QTY 1 Surgical Tape, 2 in.	 QTY 1 Trauma Shears, Lg
 PAIR 1 Nitrile Gloves, Lg	 QTY 1 Surgical Lubricant	 QTY 1 TCCC Triage Card	 QTY 1 Permanent Marker, Large	 QTY 1 Green ChemLight™





NAVAL GUN RESPONSE AID KIT *Trainer*

- Designed to meet the TCCC requirements of the U.S. Navy for the emergency care equipment onboard U.S. Naval Ships, Military Sealift Command Ships and afloat commercial vessels
- Ergo dynamic, low profile, casualty response system featuring 2 identical removable equipment pouches that attach via hook and loop to a lightweight, durable outer bag
- Meets the requirements for the U.S. Navy's Naval Gun Response Aid Kit - Trainer
- Lightweight, sand/water resistant, slimline pouches constructed of Berry Amendment 500D Cordura solution dyed nylon, Ghille Tex hardware, self repairing zippers & silent 550 cord pulls
- Both outer bag & 2 inner pouches are securable using anti-pilferage seals
- Slim-line, lightweight kit stores easily on a wall, can be carried over the shoulder or bandolier-style

NAVAL GUN RESPONSE AID KIT TRAINER

ITEM #	NSN #
85-3205	6910-01-695-9287

SPECIFICATIONS



Outer Bag:
H 16 in. x W 11 in. x D 5.5 in.
Inner Pouch:
H 15 in. x W 6 in. x D 4 in.
Weight: 1 lb 15 oz

NAVAL GUN RESPONSE AID KIT TRAINER:



QTY **2** Naval First Aid Box Response Bag

QTY **1** Naval Gun Resp.Outer Bag

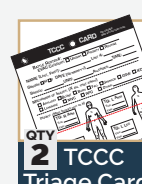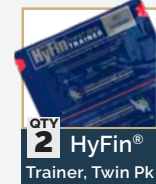



TRAINING | Individual Components

C·A·T® TRAINER

Combat Application Tourniquet

Fully functioning C-A-T® tourniquet colored blue to distinguish it as a training device as opposed to a personal/operational life-saving device



See page 6 for Tactical & Rescue versions of the C-A-T®

CAT TRAINER

ITEM #	NSN #
30-0033	6910-01-560-2972

SPECIFICATIONS

H 6.5 in. x W 2.4 in. x D 1.5 in.
Weight: 2.7 oz

HYFIN® TRAINER Twin Pack

Unique adhesive gel formulation does not damage mannequins or leave behind excessive residue

SPECIFICATIONS

Packaged: L 7.5 in. x W 4.5 in. x D 0.25 in.
Chest Seal: H 6 in. x W 6 in. x D 0.1 in.
Weight: 2.5 oz

HYFIN TRAINER

ITEM #	NSN #
10-0032	6910-01-908-8679

See page 12 for other HyFin® Chest Seal products



COMBAT GAUZE® Moulage Trainer

Completely inert (no hemostatic properties) - training purposes only in a re-sealable package for multiple uses

SPECIFICATIONS

Packaged: H 6 in. x W 4.5 in. x D 1.25 in.
Gauze: L 4 yd x W 3 in.
Weight: 0.75 oz

COMBAT GAUZE TRAINER

ITEM #
30-0063

Responder TRAINING ETD

Non-sterile, elastic pressure bandage for an economic, low cost training solution

RESPONDER TRAINING ETD

ITEM #
30-0145

SPECIFICATIONS

Bandage: L 60 in. x W 4 in.
Pad: H 6 in. x W 4 in.
Packaged: L 4.5 in. x W 2.5 in. x D 1.75 in.
Weight: 2.0 oz



Intraosseous Trainers | TRAINING



SAM IO TRAINING KIT Single Bone

SAM IO SINGLE BONE

ITEM #
30-0256

- All of the equipment needed to familiarize trainees with the SAM IO Access hardware & usage
- Featuring single training humerus bone with tissue block & rugged nylon carrying case

SPECIFICATIONS

H 7 in. x W 11.25 in. x D 5.5 in. | Weight: 2.32 lb

SAM IO TRAINING KIT - SINGLE BONE:

 QTY 1 SAM IO Driver	 QTY 1 SAM IO Adaptor	 QTY 1 SAM IO 15mm Needle	 QTY 1 SAM IO 25mm Needle
 QTY 1 SAM IO 45mm Needle	 QTY 1 SAM IO Stabilizer	 QTY 1 Humerus Bone Trainer	 QTY 1 Nylon Carry Bag

SAM IO FIELD STORAGE CASE only



SAM IO FIELD STORAGE CASE

ITEM #
30-0254

SPECIFICATIONS

H 7 in. x W 11.25 in. x D 5.5 in.
Weight: 1 lb

INTRO

M·A·R·C·H

KITS

TAC·EVAC

MED·URG

RESCUE

SUPPORT

TREAT

TRAINING

INDEX



G3 Molded Hard Side Carrying Case #30-0077 sold separately

G3 POWER DRIVER TRAINER

by EZ-IO®

- Battery-powered revolutionary technology
- Provides access within 10 seconds
- Multiple insertion sites with three needle sets
- Trigger guard prevents accidental activation of the device when stored in kits
- Standard Luer-Lock catheter for fast, easy removal
- Lightweight, compact design
- Red coloring to immediately indicate it as a training only device

G3 POWER DRIVER TRAINER

ITEM #
30-0076

SPECIFICATIONS

H 4.5 in. x W 6.5 in. x D 2.5 in.
Weight: 15 oz

G3 EDUCATOR KIT by EZ-IO®

- Includes all of the Arrow® EZ-IO® System training equipment
- Allows clinicians to train in the correct use of the EZ-IO® System on adult and pediatric patients

G3 EDUCATOR KIT

ITEM

30-0074

KIT CONTENTS:

- G3 Training Power Driver
- EZ-IO® Training Needle Sets: 15mm Pediatric, 25mm Adult, 45mm Large Adult
- 4 Wristbands
- 4 EZ-Connects
- G3 Vascular Access Pack
- 4 EZ-Stabilizers
- Training Bones: Infant, Tibia, Adult Proximal, Humerus, Adult Tibia
- Softside Nylon Case
- Training CD

EZ-IO Products are restricted to Federal Purchases only. Please contact Customer Service for details



EZ-IO® G3 Power Driver TRAINING NEEDLES

Designed for use with any EZ-IO® Power Driver




G3 TRAINING NEEDLES

SIZE	ITEM #
15 mm (Pediatric)	30-0082
25 mm (Adult)	30-0084
45 mm (Adult, large)	30-0086

T.A.L.O.N. MANUAL TRAINING NEEDLE

Train in manual I/O insertion at all seven optional sites, including the sternum

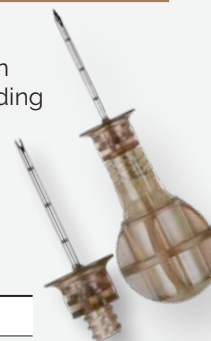

T.A.L.O.N. MANUAL TRAINING NEEDLE

ITEM #
30-0133

MANUAL TRAINING NEEDLE by EZ-IO®

For training in six different FDA cleared sites for easily accessible insertion



EZ-IO MANUAL TRAINING NEEDLE

ITEM #
30-0079

TRAINING | Part Task Trainers



WOUND MANAGEMENT SIMULATOR

WOUND MGT. SIMULATOR

ITEM #

93-0069

- Rugged part-task trainer for learning essential lifesaving skills of tourniquet application, wound packing and applying pressure bandages
- Multiple wounding patterns to enhance training objectives
- Simulates severe laceration, penetrating trauma (entrance and exit wounds), along with the ability to apply tourniquets
- Proprietary silicon that is rugged and durable and will not break down due to environmental exposure

SPECIFICATIONS

L 15 in. x Diameter 6 in.
Weight: 9.3 lb

TRUECLOT Large Caliber GSW TRAINER:



SPECIFICATIONS

Task Trainer:
L 7.25 in. x H 2.5 in.
Storage Case: L 13.5 in.
x W 10 in. x D 5.5 in.
Liquid Capacity: 1 liter
Weight:
5 lb

TRUECLOT® Hemorrhage Training Kit™ Large Caliber GSW with Bone Insert

- Large caliber gunshot wound with hidden bleed points and a simulated fracture/defect hard plastic bone insert to aid in teaching deep wound packing techniques
- Sturdy silicone construction with no moving parts, reinforced tubing mount, and a more realistic wound pattern and outer anatomy
- Ideal for the repeated use of a training environment



TRUECLOT LARGE CALIBER GSW TRAINER

ITEM #

93-0020

TRUECLOT GSW & Laceration TRAINER:



TRUECLOT® Hemorrhage Training Kit™ with 1 Gun Shot Wound & 1 Knife/Laceration

- Task Trainer for teaching realistic hemorrhage control in wounds including a large caliber gunshot & knife/laceration
- Training capabilities include general wound packing techniques & advanced hemostatic gauze training

TRUECLOT GSW & LACERATION TRAINER

ITEM #

93-0932



SPECIFICATIONS

H 17 in. x W 20 in. x D 7.5 in.
Weight: 14 lb



HEMORRHAGE CONTROL TRAINING KIT

with Bio Feedback

- Training tool designed to teach effective wound packing skills with immediate results for both trainer and student
- Verifies correct or incorrect technique

BIO FEEDBACK

ITEM

93-0019

SPECIFICATIONS

Case: 23 in. x 8 in. x 15 in. | Weight: 17 lb 3.2 oz

HEMORRHAGE CONTROL TRAINING KIT

Combat Gauze®

- Comes with 10 Combat Gauze Trainers & 5 Combat Gauze LE Demo Samples
- Includes a "training leg" that features an entry/exit bullet wound and a deep laceration

COMBAT GAUZE

ITEM

93-0017



SPECIFICATIONS

Case: H 15 in. x W 23 in. x D 8 in.
Weight: 15 lb 8 oz

TRUECLOT® Tourniquet Simulator

- Wearable device that allows full application of a tourniquet for training purposes
- Includes Trainer Sling Assembly (Left Arm Fit Only), Liter Squeeze Bottle, Quick Disconnect Fitting & Plastic Case

TRUECLOT TQ SIMULATOR

ITEM

93-0931

SPECIFICATIONS

H 12 in. x W 16 in.
x D 7 in.
Weight: 5 lb



CUSTOMER SUPPORT

- **COMMITMENT:** We pride ourselves in providing an exceptional customer experience beginning with your initial inquiry up to order placement and beyond. Our customer service staff members are trained by certified instructors on all NAR products and are standing by to help with any questions you may have
- **SUPPORT:** Our in-house technical support professionals consist of an elite group of Nationally Registered Paramedics and Special Operations Medics from military, law enforcement, tactical law enforcement, fire-rescue and EMS backgrounds

ORDERING

- **ONLINE:** www.NARescue.com
SPECIAL PRICING: automatically available ONLINE for all military customers with a .mil e-mail address! Register your account in just minutes by creating an online account using your .mil e-mail address

ESTABLISHED CUSTOMERS:

Want to order online with special pricing? Contact Customer Service at 888.689.6277 to have your online account set up

NEW CUSTOMERS: Simply create an online account and you will have access to customer specific pricing – anytime, anywhere you are logged in to NARescue.com

QUOTE REQUESTS: Online accounts have no waiting to receive your custom pricing to place an order or, we will gladly accept your submittal for a pricing quote request

- **CONTRACT:**
Tax ID: 27-1024029
Cage Code: 06ST7
DUNS: 83-24-26782
SAM Registered: Yes
ECAT: dmmonline.dscp.dla.mil
Prime Vendor: dmmonline.dscp.dla.mil/primevendor/primehome
- **EMAIL:** info@narecue.com
- **PHONE:** 888.689.6277
- **ADDRESS:**
35 Tedwall Ct | Greer, SC 29650
- **FAX:** 864.675.9880

SHIPPING

- In-stock merchandise ships within two business days
- Need it ASAP? Ask Customer Service about overnight or second day delivery



PRODUCTS
WITH A
MISSION



North American Rescue® • 35 Tedwall Court • Greer, SC 29650 • 888.689.6277 • NARescue.com

MILITARY CATALOG | Volume 18 | Number 1 | REV081021