

The Quantum is a Game-Changing innovation that has succeeded in closing a critical capability gap previously left exposed by other commercial off-the-shelf fluid warmers currently fielded. Size, weight, performance, and ease of use are no longer collective barriers to the delivery of whole blood in far forward and near POI casualty management.

RUCATION System

REV050622 The Quantum exceeds all healthcare and logistical requirements for the delivery of warmed blood or fluids across the defined spectrum of care. The system consists of unique thermal administration sets that incorporate sensing elements with an integrated heating system that does not disrupt the fluid path. With advanced intelligence throughout the system, the lightweight Controller modulates and moderates the energy flow to the integrated heating elements to ensure consistent normothermic fluid delivery to the casualty. The 44V battery is the smallest in category,

able to rapidly adjust power distribution based on the temperature and flow rate of the input fluid. This creates a highly efficient system capable of previously unrealized performance optimizations. The Quantum system has an operational weight of just 22 total ounces yet can deliver two units of 39°F (4°C) Whole Blood (WB) at 100mL/min with an output temperature of 100.4°F (38°+/- 2°C) with only a 24 second warm-up period. The Quantum Blood and Fluid Warmer provides the Tactical Health Care Professional with significant technological overmatch in the fight to eliminate preventable combat death.

#### CONTUNER™ UCALIFEWARMER™ IV Fluid & Blood Warming System

- **Lightweight** The Quantum<sup>™</sup> has the lightest operational weight of any portable blood & fluid warming system (Total Operational Weight: 22 oz)
- **Powerful** Warms two units of 39°F (4°C) refrigerated whole blood at 100mL/min+ on a single charge
- **Compact** Fewest components to carry (three)
- **Fast** Warm-up time of only 24 seconds
- **Intuitive** The thermal transfusion/ infusion set looks, handles, and sets up/primes like any conventional IV set

**IV & BLOOD WARMING SYSTEM CONTENTS:** 

Quantum

Controller

 Safe Tested safe for blood, blood products & IV solutions

1

**COMPLETE QUANTUM** 

Quantum

Battery

- **Effective** The system is plug-and-play for ease of setup under stressful conditions such as tactical field care (only two connection points, with no plug orientation or set-up order requirements – system status constantly visible at a glance)
- **Complete** Does not require additional standard transfusion or infusion set all other warmers do
- **Tactical** Only system FDA cleared allowing easy disable of audible alerts to protect operators in combat

SPECIFICATIONS

(Quantum<sup>™</sup>Controller,

2

Thermal

**Fransfusion Set** 

Total Operational Weight: 22 oz

Battery & TTS-B)

Battery: 16 oz Controller: 3 oz Battery Charger: 22 oz Controller Cover: 0.23 oz Thermal Transfusion Set (TTS-B): 2.7 oz

#### **BLOOD** or **BLOOD PRODUCTS**

From a *single battery charge*, capable of warming *1000 mL* of *39°F (4°C) Whole Blood* to an output temperature of *100.4°F (38°+/- 2°C)*. Warm up time: *24 seconds*. Delivery Flow rate: *100mL/min*.

#### QUANTUM COMPLETE SYSTEM

ITEM #	NSN #	
35-0006	6515-01-693-7810	

#### IV FLUIDS

From a *single battery charge*, capable of warming **1700 mL** of **68°F (20°C) IV fluids** to an output temperature of **100.4°F (38°+/- 2°C)**. Warm up time: **24 seconds**. Delivery Flow rate: **200mL/min**.

## Crĕdo **PROMED**<sup>™</sup> SERIES FOUR CONTAINER

1

Quantum

Charger

- Temperature-controlled portable medical transport bags that thermally protect the integrity of valuable medical supply payloads for 48 - 96 hours
- Outer Container: Ballistic nylon fabric with durable web strapping and secure buckles
- Featuring the patented Thermal Insulation Chamber (TIC<sup>™</sup>) panels with phase change material and VIP components contained within to consistently protect medical materials such as blood and platelets within 2-8°C.
- Available in 3 capacities: 2L, 4L & 8L

#### CREDO PROMED™ SERIES FOUR CONTAINER

CAPACITY	ITEM #
2 units	35-0011
4 units	35-0012
8 units	35-0013

SPECIFICATIONS

**2 Units:** Volume: 2L Duration: 48 hours L 8 in. x W 9 in. x D 10 in. Weight: 7.64 lb

**4 Units:** Volume: 4L Duration: 96 hours L 10 in. x W 10.5 in. x D 10.5 in. Weight: 11.45 lb

**8 Units:** Volume: 8L Duration: 72 hours L 11 in. x W 15.5 in. x D 11 in. Weight: 19.6 lb

PROMED

Ideal for use with the Quantum and LifeFlow devices to provide a complete blood delivery package

# TTS-B Thermal Transfusion Set

ITEM #

35-0005

- Compatible with blood, blood products, and IV solutions
- Similar size and weight to traditional IV tubing
- Intelligent sensing & heating elements integrated into the walls of the tubing set
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller
- 20 gtts/mL drip rate

\*Retain a ready charge by keeping Battery plugged into Charger. Always re-charge after each use, even if not fully discharged.

#### QUANTUM BATTERY

ITEM # NSN #

35-0001

Quantum CONTROLLER

- Moderates the energy flow from the Quantum Battery to the Thermal Tubing
- Designed for ease of connecting and durability
- Colored LED lights constantly show during administration of blood/fluid to show status
- Lightweight & durable
- IP53 rated for dust & water resistance with jacket

#### **QUANTUM CONTROLLER**

CART OF CONTROLLER			
OPTION	ITEM #	NSN #	
Controller	35-0003		
Cover (Replacement)	35-0010		

**SPECIFICATIONS** 

TTS-B package (each): 6.5 in. x 11 in. Tubing Length: 80 in. TTS-B Wgt (each): 2.6 oz

Packaged Wgt (2 TTS-B & box): 11 oz

TTS-B (2 Pack)

6515-01-695-4444

NSN #

Quantum **BATTERY** 

environment- water resistant,

subdued color, and the option

to disable the audible alarm

44V Lithium-Polymer Battery

• Intelligent, state-of-the-art

charge demand and will

**SPECIFICATIONS** 

overcharge\*

state at any time

Weight: 16 oz

Replacement cover

is also available

battery technology senses

charge as needed but never

Status button to check charge

H 4.5 in. x W 3.25 in. x D 1.25 in.

Minimum service life -

W 1.25 in. x H 3 in. x D 0.75 in.

Controller Total Length: 51.25 in.

Weight (Controller only): 2.5 oz

**SPECIFICATIONS** 

Packaged Weight: 5 oz

Controller Body:

1000 insertions/removals

Designed for the tactical

**TIS** Thermal Infusion Set

- Similar size and weight to traditional IV tubing
- Intelligent sensing and heating elements integrated into the walls of the tubing set **SPECIFICATIONS**
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller
- Compatible with standard IV fluids such as colloids and isotonics
- 20 gtts/mL drip rate

#### TIS (2 Pack)

ITEM # NSN # 35-0004 ----

#### Quantum HARDCASE

- Rugged hard case designed to protect your Quantum from dust, dirt, and moisture
- Features pressure relief valve for automatic pressure equalization
- Custom laser cut high density foam holds components securely in place
- Flexible storage capacity for multiple batteries, controllers and thermal tubing sets

#### QUANTUM HARDCASE


NSN # ITEM # 35-0008 ----

- Fully charge the Quantum Battery in < 90 minutes
- Smart, safe charging
- Includes 6 ft AC power cord
- Colored LED lights to inform the user of charging and battery status
- IP 22 rated for dust and water resistance
- 100V to 264V AC input; 50Hz. 60Hz
- Minimum service life-2500 insertions/removals

#### QUANTUM CHARGER

NSN # ITEM #

35-0002 \_\_\_

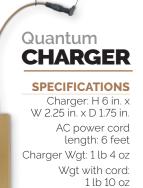



TIS packaged (each): 6.5 in. x 11 in. Tubing Length: 80 in. Weight: TIS (each): 2.1 oz Packaged (2 TIS & box): 10 oz







# NUCATUM Field Blood Transfusion Kit

 Comprehensive kit designed for use with the Quantum Blood & Fluid Warming System allowing rapid, consistent normothermic field blood transfusion in far forward and near point of injury (POI) prehospital casualty management.

#### QUANTUM FIELD BLOOD TRANSFUSION KIT CONTENTS:

- 1 2
- Includes the Quantum Thermal Transfusion Set -(Blood) that incorporates an integrated heating system that does not disrupt the fluid path
- Equipped with Donor/Recipient Modules to verify blood types of the donor and recipient
- Also allows collecting of blood samples from the donor for additional surveillance testing

Each Donor Module Kit contains





1 450ml Blood



### **QUANTUM FIELD BLOOD TRANSFUSION KIT**

RECIPIENT

**ITEM #** NSN # 83-0028 \_\_\_

All the needed supplies to allow the user to rapidly perform a field blood transfusion as outlined in the **Joint Theater Trauma System Clinical Practice** Guidelines

COMPANY OF THE OWNER.

DONOR

# Q.U.I.C.K. Quantum Ultimate Intensive Resuscitation Carry

- Combines the blood and fluid warming system with individually packaged NAR Saline Locks
- MOLLE/PALS configured nylon pouch, equips providers with an optimized solution for carriage to the point of need.

#### **PRODUCT WEIGHTS**

Batterv: 16 oz Controller: 3 oz Thermal Transfusion Set (TTS-B): 2.7 oz Charger: 22 oz Total Operational Weight: 22 oz (Controller, Battery & TTS-B)



www.NARescue.com • 888.689.6277





Q.U.I.C.K.	Q.U.I.C.K. CONTENTS:				
ory 1 Quntum	ary 1 Quantum	atry 1 Quantum	ary Thermal	aty 2 NAR	
Battery	Controller	Charger	Transfusion Set	Saline Lock Kit	
	Q.U.I.C.K. Advanced also contains:				



#### Q.U.I.C.K.

ITEM #	OPTION	NSN #
35-0026	Basic	
35-0025 Advanced		

